


Wear & Care Instructions

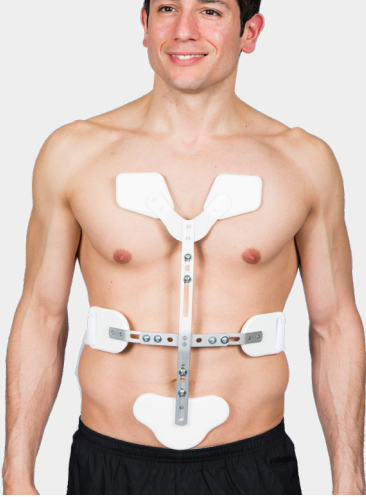

CASH Hyperextension Brace

Steeper manufactures a quality range of orthotic products and to help you care for the product, this leaflet provides information on how and when to use your prescribed orthosis, how to care for it and how to contact us if you have a problem or require a replacement.

Care & Cleaning

Wash with a sponge dampened in warm water and mild soap. Dry with soft cloth. Check the screws periodically and tighten any loose ones. Contact an orthopaedic expert for replacement of worn parts. Do not dispose the orthosis or any parts carelessly.

Applying the Orthosis

1. Open the waist belt on one side.
2. Position the brace against the chest, pass the belt behind the back and fasten it through the opening.
3. Position the brace so that the sternal plate is centred on the chest, and the triangular pelvic pad is not uncomfortable when sitting.

IMPORTANT

The advice and information contained in this document is not intended as a substitute for medical counselling. If you have any questions about treatment for a specific condition, please check with a health care professional.



Unit 3 Stourton Link,
Intermezzo Drive,
Leeds LS10 1DF
Tel: +44 (0) 113 207 0435
Email: orthoticenquiries@steepergroup.com

www.steepergroup.com

 **steeper**
Creating life's turning points, together