do more™

taskaprosthetics.com

10 Nelson Street, Riccarton, Christchurch 8011, New Zealand.

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 **TASKA™**

Introduction

This guide is designed to help you get the most from your TASKA Hand. There are over 20 pre-programmed grips, and you can create your own using the MyTASKA app, which you can find on your preferred app store.

Scan the QR Code for video guides and other resources.

For more information on your TASKA Hand talk with your prosthetist, or visit taskaprosthetics.com

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Core grips

The grips you will use for most of your daily activities based on three fundamental grasps.

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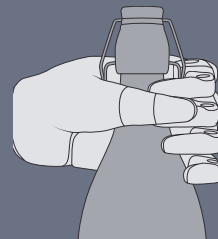
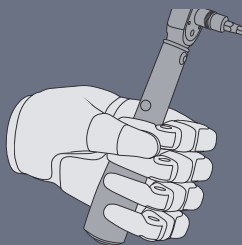
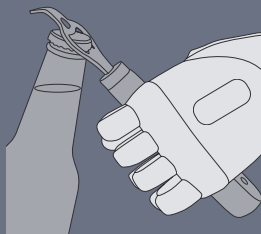
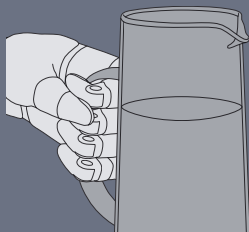
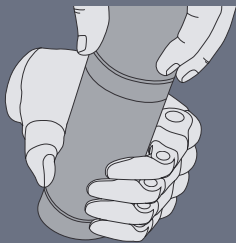
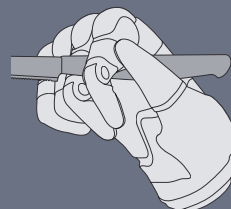
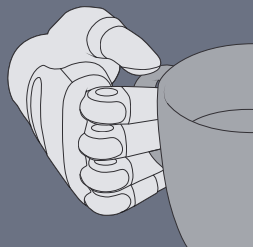
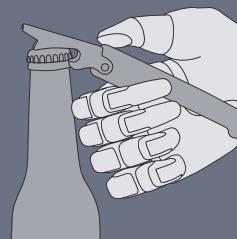
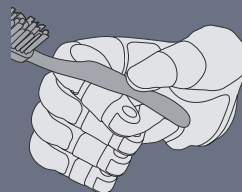
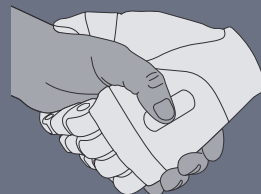
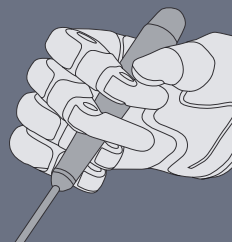
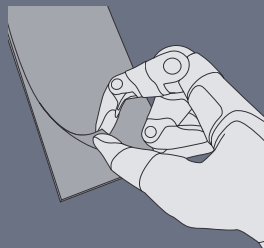
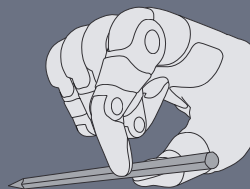
Melanie Rivera Waldman

“

I've been using my TASKA to **find my way back** to the things I love the most.”

Melanie Rivera Waldman
Actor, yoga teacher, differently
limbed advocate, TASKA Champion.

At her core, Melanie is someone who finds a way. She is on a journey of adaptation and connection with her authentic self.

Read Melanie's full story and
meet our other Champions at
taskaprosthetics.com

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Dynamic grips

Grips intended to focus on single movements of digits to keep grasping simple and predictable.

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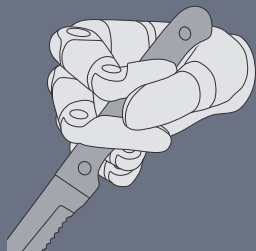
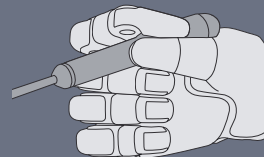
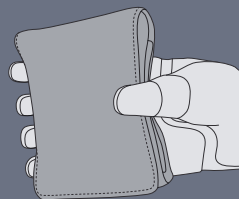
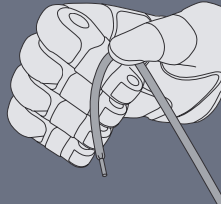
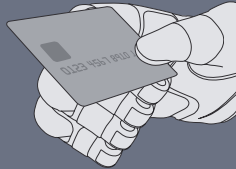
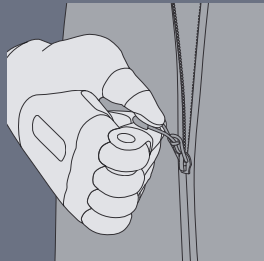
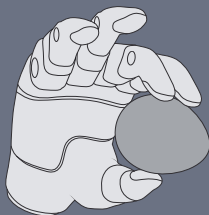
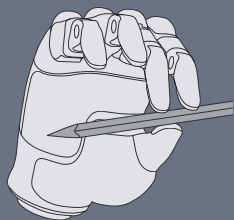
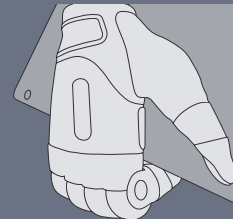
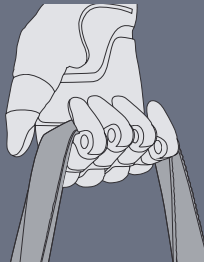
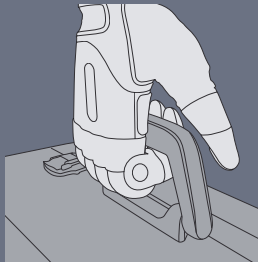
Bryden Zaloum

“

There's simply no other myoelectric hand on the market **that can do what it does.**”

Bryden Zaloum
Passionate car restorer,
TASKA Champion.

Achieving the right balance between comfort and functionality eluded Bryden for more than three decades - until he discovered the TASKA Hand.

Read Bryden's full story and meet our other Champions at taskaprosthetics.com

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Static grips

Natural and resting positions used for supporting objects where dynamic movement is not required.

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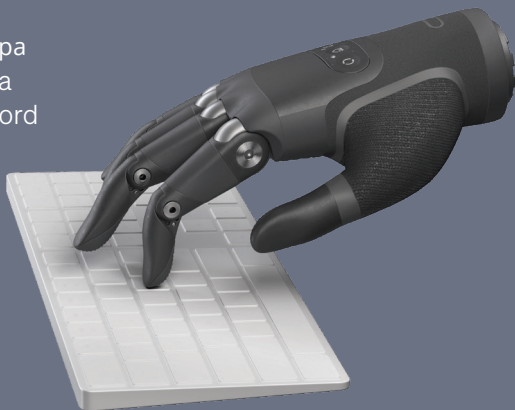
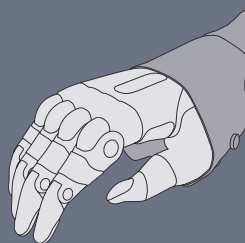
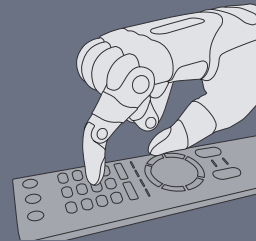
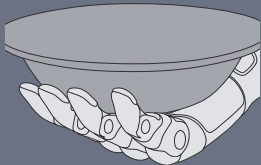
Nic Hite

“

The TASKA Hand can **hold on to a lot more weight** than I thought!”

Nic Hite
Engineer, gym enthusiast,
TASKA Champion.

Nic has always been a confident person, and believes his TASKA helps him transfer that confidence to others.

Read Nic's full story and
meet our other Champions at
taskaprosthetics.com

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