

Aequi Ankle Brace



Features and Benefits

- The Push Ortho Aequi is an effective ankle orthosis which strikes a balance between efficient mechanical support of a disabled ankle and good wearer comfort
- The rigid support is situated exclusively on the medial side of the ankle and has a recess for the malleolus
- The lateral side is made from preformed foam, which easily adapts to the individual anatomy

Indications for Use

- Serious (residual) instability (under high stress)
- Treatment of acute lateral ankle ligament lesions
- After-care of conservative or operated ankle fractures
- Arthrosis with instability without accumulation of fluid
- Secondary prevention of ankle ligament lesions

Size	Circ. of Short Heel (cm)	Product Code Left	Product Code Right
1	27 - 31	320111	320121
2	31 - 34	320112	320122
3	34 - 40	320113	320123

Prim Dorsi Assist




Features and Benefits

- Consists of two separated pieces, a supramalleolar strap over the ankle and a non-slip support union to fix between the tongue and shoelaces by means of a central anchorage hook type. Both parts are joined by an elastic strap, which allows variable amounts of dorsiflexion assist
- On the back of the ankle band, there is a gel pad which is cushioned for increased comfort on the achilles tendon
- Discrete, flexible, it takes no space on the shoe, being easy to fit
- To be fitted when lacing the shoe
- Has the function to produce dorsiflexion on the take-off phase of gait, indicated in flaccid palsy

Indications for Use

- Drop foot
- Peroneal neuropathies
- Dysfunctional gait

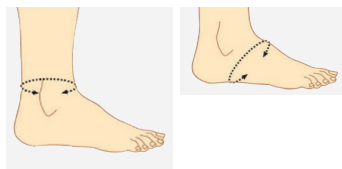

Measurement

Size	Measurement (cm)	Product Code
X-Small	13 - 17	ATXO1/XS
Small	17 - 21	ATXO1/S
Medium	21 - 25	ATXO1/M
Large	25 - 29	ATXO1/L

Plantar Band		
Size	Measurement (cm)	Product Code*
Small	17 - 21	AXT11/S
Medium	21 - 25	AXT11/M
Large	25 - 29	AXT11/L

*Please specify right or left.