



Tred-lite
Modular Footwear Collection

 **steeper**
by eqwal.

Diabetic specifications

Softee leather, moleskin lining with minimum seams, rim toe puff and padded collars can be added to all styles to create a complete diabetic specification.

Lightweight, practical and customisable footwear

Available in 17 leather (Calf Softee) colour options, Nubuck and Neoprene, the lightweight Tred-lite range is practical, hard wearing and customisable, offering a mixture of styles, fastenings and colour combinations.

Adaptations on Tred-lite

Footwear adaptations and trial fitting service available, along with repairs and adapts.

Specialist sole alternatives

Standard sole and heel provided in lightweight micro EVA, with specialist Vibram® hard-wearing commando also available.

Bespoke insoles ensure best clinical outcome

With a range of insole densities, complementary insoles manufactured at the same time as the footwear ensure the best possible clinical outcome, particularly for total-contact insoles to distribute load and reduce peak pressures on vulnerable areas.

Shock absorbing, lightweight insole support

Using the highest quality EVA, our bespoke insoles retain their shape.

The semi-bespoke Tred-lite collection offers a variety of fashionable styles to provide comfortable and lightweight footwear for men and ladies. As an introduction to the range, this guide includes the current collection and presents the wide range of colours, options and combinations to assist you with your order.

Colour and Material Options

The Tred-lite footwear range is available in Calf Softee Leather, nubuck and neoprene material. Certain styles may be limited to one type of material depending on the design. The shoes in this brochure are shown in standard colours and additional colour options (outside of the standard range) to illustrate a mixture of colour and style combinations. A small surcharge applies to colours (e.g. patent) that are outside of the standard collection.

- **Calf Softee Leather** - Chianti, Merlot, Easter Rose, Ruby, Jaffa, Daffodil, Shamrock, Bottle Green, Waterfall, Baby Blue, Life Boat, Madeira Blue, Purple, Lilac, Barbie Pink, Candy Floss, Baby Pink, White, Milkwood, Taupe, Caramel, Acorn, Chestnut, Willow, Havana, School Grey and Black.
- **Nubuck** - Blush Nubuck, Cerise Nubuck, Mustard Nubuck, Wet Sand Nubuck, Milk Chocolate Nubuck, Black Coffee Nubuck, French Navy Nubuck, Midnight Blue Nubuck, School Grey Nubuck and Jet Black Nubuck.
- **Neoprene** - Beige Neoprene, Brown Neoprene, Navy Neoprene and Black Neoprene.

For accurate depiction of colours and shades, please contact our Customer Services department for a swatch sample. A digital swatch palette is also available at www.steepergroup.com, however please note that accuracy may vary between computer screens, therefore a physical swatch should always be referenced.

Strap Options

Round laces are supplied as standard, and most styles are available in hook and loop fastening, however certain styles are limited to lace only or hook and loop only depending on the design. Low opening styles are available as standard and some styles can be made into a low opening option on request - by 15mm or by 25mm, please specify this on your order form. See icons for details.

Sole Options

We provide the standard sole and heel in micro EVA, however we can also provide Vibram® hard-wearing commando as an alternative sole option (please note there is an additional cost for this). If you require a different sole from the standard, please specify this on the order form.

Default Sole Colours

The following table outlines the default sole colour to match the upper:
Please note: all trainer soles are white by default.

Upper	Default Sole Colour	Upper	Default Sole Colour
French Navy	Black	Cherry Red	Black
Navy	Black	Pillar Box Red	Black
Black	Black	Purple	Black
White	White	Burgundy	Black
Light Beige	Brown	Tan	Brown
Dark Beige	Brown	Dark Brown	Brown
Taupe	Brown	Mid Brown	Brown
Mid Grey	Black	Bottle Green	Black
Fuschia	Black		

Insoles

Tred-lite footwear is provided with 2 x 3mm Poron® material insoles and 1 x 3mm top cover.

Sizing and Lasts

Different last fittings are available. Please see the towards the back of these pages for sizing charts.

Diabetic Footwear

The majority of styles are made to a diabetic specification, which incorporate softee leather uppers, antelope quarters with moleskin lining, minimal seams, padded collars, rim toe puff, and leather quarter lining. We are able to manufacture your preferred style in a diabetic specification, please clearly state your requirements on the order form.

Neoprene Footwear

Moleskin vamp linings and Rim toe puff are available on request. Please state if you would like these adaptations when placing an order.

Footwear Adaptations

We carry out a full range of footwear adaptations and also offer a repair and adaptation service for patients' own footwear. Please get in touch with our Customer Services department for further details.

Tred-lite Insoles

Steeper can provide custom made insoles using either foam impressions, casts or scans of the foot that are manufactured to the highest standards. Frequently, clinicians use total-contact insoles alongside modular footwear to distribute load and reduce peak pressures on vulnerable areas. As the interface between footwear and the foot, the correct insole can ensure the best possible clinical outcome. By manufacturing footwear and insoles together, we ensure that each device complements the other and functions in conjunction.

Base Materials for CAD-CAM

Our bespoke manufactured insoles use the highest quality EVA (Ethylene Vinyl Acetate). This lightweight material offers shock absorption and retains its shape through its hard wearing properties.

The densities available are: low, medium, high, soft dual and firm dual density. Other options are also available such as Nora, Plastazote and Polyurethane. Please contact Customer Services for a complete list of materials.

Top Covers Available:

- Grey poron and medical grade porons
- Multi-layered materials including Nora™
- Memory foams, sponge rubber
- Anti-bacterial covers
- Leather, blue oyster/pearl, gas meter, black cape

Please contact our Customer Services department on orthotics@steepergroup.com for a complete list of top cover materials available to order.



Ordering & Delivery Information

How to Order

Creating footwear exactly to our customers' specifications provides an increased likelihood of a perfect fit for the patient. Therefore, please note that providing fully correct information is essential at the initial order placing stage, as re-working such as relasting or crafting new uppers after the rough fitting will be charged for.

Please also note that all footwear made to order or with any customisation that is incorporated during manufacture will render the footwear non-returnable.

When ordering a product from the Tred-lite footwear range, please include all the information below to enable us to produce your order accurately to the required specifications, as well as deliver the item on time.

- Style (name of shoe)
- Colour (if applicable, the colour combinations and numbers)
- Fastening option
- Material
- Size and width measurements
- Diabetic specifications
- Any other adaptations or amendments

Please send your completed order form to orthotics@steepergroup.com, including all the required details. You may wish to refer to the sizing charts for specific Tred-lite shoe, trainer and boot sizes. For any order queries or to request an order form, please contact us on 0113 207 0435.

Trial Fitting Service

As part of our customer service, we offer a free trial fitting service on all our modular footwear. This allows for minor alterations before the footwear is finished. Charges may apply for additional trial fittings.

It is essential that you specify on your order form if a trial fitting is required, otherwise the footwear will be finished. Major alterations after a trial fit e.g. Uppers re-cut, will incur a charge.

Returns & Refunds

Shoe Styles

The shoes listed below are returnable for a full refund providing that they are returned within 28 days of receipt, in all of their original packaging and in full re-saleable condition.

- Ladies' styles on standard lasts (sizes 3-6) without adaptations/modifications - in black, brown and navy
- Men's shoe styles on standard lasts (sizes 8-11) without adaptations/modifications - in black and brown

Boot Styles

Boots provided in all of their original packaging and in full resalable condition may be returned but will incur a 30% restocking charge, providing they are returned within 28 days of receipt.

- Ladies' styles on standard lasts (sizes 3-6) without adaptations/modifications - in black, brown and navy
- Men's styles on standard lasts (sizes 8-11) without adaptations/modifications - in black and brown

Non-Returnable Items

All patient-specific footwear (shoes and boots) including made-to-measure, non-standard last fittings, non-standard size footwear (see sizes above), split sizes, fittings, split fittings, special colours, last or upper altered or adapted in any other way are non-returnable.

Single shoes or boots, or worn footwear are non-returnable.

Any footwear requested with any other sole variant to that specified in the sole options table, will not be returnable unless the sole unit is standard for that type of footwear.

Note: If you've received a non-returnable product in a damaged/defective condition, please contact customer services on receipt of your item.

Please contact customer services if you require order forms, any further information or assistance with measurements, adaptations and ordering.


Ladies' Collection

These designs are available in a range of colours, materials, soles and fastening options, offering a stylish and comfortable solution.


The Tred-lite range covers shoe, trainer and boot style options.

Please note: an * indicates a surcharge which is applied to colours that are not part of the standard colour collection. Please refer to the colour swatch for confirmation. Please also note that photos are to demonstrate the style of shoe, and that upper and sole combinations may not reflect the default.


Key




Real leather




Nubuck leather




Neoprene




Lace up round




Lace up flat




Hook & Loop fastening




Low opening 15mm




Low opening 25mm



Diabetic design










Available as bespoke made to measure




Corella

Standard sole: Half bridge micro EVA

*Shown in fuchsia pink leather











Eleania

Standard sole: Half bridge micro EVA

Please specify contrasting colour combinations












Mandarin

Standard sole: Half bridge micro EVA

Please specify contrasting colour combinations










Oriole

Standard sole: Half bridge micro EVA









*Shown in purple nubuck






Starling







Standard sole: Half bridge micro EVA






Serin

Standard sole: Half bridge micro EVA











Ibis

Standard sole: Half bridge micro EVA









Low opening style as standard






Kea








Standard sole: Half bridge micro EVA






Wheatear









Standard sole: Half bridge micro EVA





Robin

Standard sole: Half bridge micro EVA



www.steepergroup.com

To order call 0113 207 0435

Ladies' Collection



Minivet

Standard sole: Half bridge micro EVA





Nightingale


Standard sole: Half bridge micro EVA





Jacamar

Standard sole: Half bridge micro EVA
A low opening style as standard.
Please specify contrasting colour combination.






Hummingbird


Standard sole: Half bridge micro EVA





Nunlet

Standard sole: Half bridge micro EVA





Piculet


Standard sole: Half bridge micro EVA





Murrelet

Standard sole: Trainer micro EVA
Please specify contrasting colour combinations





Myna

Standard sole: Trainer micro EVA
Please specify contrasting colour combinations





Aracari

Standard sole: Half bridge micro EVA






Gila


Standard sole: Half bridge micro EVA






Jay


Standard sole: Trainer micro EVA
Please specify contrasting colour combinations





Lora

Standard sole: Half bridge micro EVA



Ladies' Collection



Phoebe
Standard sole: Square lightweight commando
Please specify contrasting colour combinations.

Real leather

Nubuck leather

Lace up flat

Hook & Loop fastening

Low opening 15mm

XLO

D

B



Rosefinch
Standard sole: Half bridge micro EVA
Contrasting colour on side panel optional.
Please specify contrasting colour combinations.

Real leather

Nubuck leather

Lace up flat

Hook & Loop fastening

B

Men's Collection

The Tred-lite men's collection covers a range of styles from casual wear to smart, offering comfort and options for an active lifestyle. This collection is available in a wide range of colours, soles, materials and fastening options.

Please note: an * indicates a surcharge which is applied to colours that are not part of the standard colour collection. Please refer to the colour swatch for confirmation. Please also note that photos are to demonstrate the style of shoe, and that upper and sole combinations may not reflect the default.

Key

Real leather

Nubuck leather

Neoprene

Lace up round

Lace up flat


Hook & Loop fastening

Low opening 15mm

XLO

D

B



Falcon
Standard sole: Square micro EVA

Real leather

Nubuck leather


Lace up flat

Hook & Loop fastening

Low opening 15mm

D

B



Jacana
Standard sole: Square micro EVA

Real leather

Nubuck leather

Lace up flat

Hook & Loop fastening

Low opening 15mm

XLO

D

B



Fulmar
Standard sole: Square micro EVA

Real leather

Nubuck leather

Lace up flat

Hook & Loop fastening

Low opening 15mm

D

B



Kestrel
Standard sole: Square micro EVA

Real leather

Nubuck leather

Neoprene

Lace up flat


Hook & Loop fastening

Low opening 15mm

D









B

Men’s Collection



Kite

Standard sole: Square micro EVA





Shearwater

Standard sole: Square micro EVA






Hobby






Standard sole: Trainer micro EVA
Please specify contrasting colour combinations.
*Shown in grey nubuck





Swift




Standard sole: Trainer micro EVA
Please specify contrasting colour combinations.






Manx







Standard sole: Square micro EVA
Low opening style as standard.






Whistler









Standard sole: Square micro EVA






Curlew









Standard sole: Half bridge micro EVA






Harrier









Standard sole: Square lightweight commando





Prion

Standard sole: Square micro EVA





Peregrine

Standard sole: Trainer micro EVA
Low opening as standard.
Please specify contrasting colour combinations.



Sizing Chart - Ladies' Tred-*lite*

	Size	Heel Width (mm)	Joint Width (mm)	Joint Circumference (mm)	Instep Circumference (mm)	Toe Depth (mm)	Last Length (mm)
Standard	2	60.25	89	216 - 234	206 - 224	16.25 - 25.25	233
	3	61.50	91	222 - 240	212 - 230	17.50 - 26.50	241
	4	62.75	93	228 - 246	218 - 236	18.75 - 27.75	249
	5	64.00	95	234 - 252	224 - 242	20.00 - 29.00	257
	6	65.25	97	240 - 258	230 - 248	21.25 - 30.25	265
	7	66.50	99	246 - 264	236 - 254	22.50 - 31.50	273
	8	67.75	101	253 - 271	243 - 261	23.75 - 32.75	281
	9	69.00	103	259 - 277	249 - 267	25.00 - 34.00	289
Narrow	10	70.25	105	265 - 283	255 - 273	26.25 - 35.25	297
	2	56.25	81	204 - 222	196 - 214	16.25 - 25.25	233
	3	57.50	83	210 - 228	202 - 220	17.50 - 26.50	241
	4	58.75	85	216 - 234	208 - 226	18.75 - 27.75	249
	5	60.00	87	222 - 240	214 - 232	20.00 - 29.00	257
	6	61.25	89	228 - 246	220 - 238	21.25 - 30.25	265
	7	62.50	91	234 - 252	226 - 244	22.50 - 31.50	273
	8	63.75	93	241 - 259	233 - 251	23.75 - 32.75	281
Wide	9	65.00	95	247 - 265	239 - 257	25.00 - 34.00	289
	10	66.25	97	253 - 271	245 - 263	26.25 - 35.25	297
	2	62.25	96	229 - 247	216 - 234	16.25 - 25.25	233
	3	63.50	98	235 - 253	222 - 240	17.50 - 26.50	241
	4	64.75	100	241 - 259	228 - 246	18.75 - 27.75	249
	5	66.00	102	247 - 265	234 - 252	20.00 - 29.00	257
	6	67.25	104	253 - 271	240 - 258	21.25 - 30.25	265
	7	68.50	106	259 - 277	246 - 264	22.50 - 31.50	273
Deep	8	69.75	108	266 - 284	253 - 271	23.75 - 32.75	281
	9	71.00	110	272 - 290	259 - 277	25.00 - 34.00	289
	10	72.25	112	278 - 295	265 - 283	26.25 - 35.25	297
	2	60.25	89	226 - 244	211 - 229	22.25 - 31.25	233
	3	61.50	91	232 - 250	217 - 235	23.50 - 32.50	241
	4	62.75	93	238 - 256	223 - 241	24.75 - 33.75	249
	5	64.00	95	244 - 262	229 - 247	26.00 - 35.00	257
	6	65.25	97	250 - 268	235 - 253	27.25 - 36.25	265
Wide & Deep	7	66.50	99	256 - 274	241 - 259	28.50 - 37.50	273
	8	67.75	101	263 - 281	248 - 266	29.75 - 38.75	281
	9	69.00	103	269 - 287	254 - 272	31.00 - 40.00	289
	10	70.25	105	275 - 293	260 - 278	32.25 - 41.25	297
	2	62.25	96	238 - 256	218 - 236	22.25 - 31.25	233
	3	63.50	98	244 - 262	224 - 242	23.50 - 32.50	241
	4	64.75	100	250 - 268	230 - 248	24.75 - 33.75	249
	5	66.00	102	256 - 274	236 - 254	26.00 - 35.00	257


Please note: Half sizes are available on request, please specify on your order form.

Sizing Chart - Men's Tred-*lite*

	Size	Heel Width (mm)	Joint Width (mm)	Joint Circumference (mm)	Instep Circumference (mm)	Toe Depth (mm)	Last Length (mm)
Standard	5	63.25	94	235 - 253	226 - 244	22.25 - 31.25	257
	6	64.50	96	241 - 259	232 - 250	23.50 - 32.50	265
	7	65.75	98	247 - 265	238 - 256	24.75 - 33.75	273
	8	67.00	100	253 - 271	244 - 262	26.00 - 35.00	281
	9	68.25	102	259 - 277	250 - 268	27.25 - 36.25	289
	10	69.50	104	265 - 283	256 - 274	28.50 - 37.50	297
	11	70.75	106	271 - 289	263 - 281	29.75 - 38.75	305
	12	72.00	108	277 - 295	269 - 287	31.00 - 40.00	313
Narrow	13	73.25	110	283 - 301	275 - 293	32.25 - 41.25	321
	14	74.50	112	289 - 307	281 - 299	33.50 - 42.50	329
	5	61.25	88	223 - 241	218 - 236	22.25 - 31.25	257
	6	62.50	90	229 - 247	224 - 242	23.50 - 32.50	265
	7	63.75	92	235 - 253	230 - 248	24.75 - 33.75	273
	8	65.00	94	241 - 259	236 - 254	26.00 - 35.00	281
	9	66.25	96	247 - 265	242 - 260	27.25 - 36.25	289
	10	67.50	98	253 - 271	248 - 266	28.50 - 37.50	297
Wide	11	68.75	100	259 - 277	255 - 273	29.75 - 38.75	305
	12	70.00	102	265 - 283	261 - 279	31.00 - 40.00	313
	13	71.25	104	271 - 289	267 - 285	32.25 - 41.25	321
	14	72.50	106	277 - 295	273 - 291	33.50 - 42.50	329
	5	66.25	102	248 - 266	236 - 254	22.25 - 31.25	257
	6	67.50	104	254 - 272	242 - 260	23.50 - 32.50	265
	7	68.75	106	260 - 278	248 - 266	24.75 - 33.75	273
	8	70.00	108	266 - 284	254 - 272	26.00 - 35.00	281
Deep	9	71.25	110	272 - 290	260 - 278	27.25 - 36.25	289
	10	72.50	112	278 - 296	266 - 284	28.50 - 37.50	297
	11	73.75	114	284 - 302	273 - 291	29.75 - 38.75	305
	12	75.00	116	290 - 308	279 - 297	31.00 - 40.00	313
	13	76.25	118	296 - 314	285 - 303	32.25 - 41.25	321
	14	77.50	120	302 - 320	291 - 309	33.50 - 42.50	329
	5	63.25	94	242 - 260	229 - 247	28.25 - 37.25	257
	6	64.50	96	248 - 266	235 - 253	29.50 - 38.50	265
Wide & Deep	7	65.75	98	254 - 272	241 - 259	30.75 - 39.75	273
	8	67.00	100	260 - 278	247 - 265	32.00 - 41.00	281
	9	68.25	102	266 - 284	253 - 271	33.25 - 42.25	289
	10	69.50	104	272 - 290	259 - 277	34.50 - 43.50	297
	11	70.75	106	278 - 296	266 - 284	35.75 - 44.75	305
	12	72.00	108	284 - 302	272 - 290	37.00 - 46.00	313
	13	73.25	110	290 - 308	278 - 296	38.25 - 47.25	321
	14	74.50	112	296 - 314	284 - 302	39.50 - 48.50	329
	5	66.25	102	254 - 272	237 - 255	28.25 - 37.25	257
	6	67.50	104	260 - 278	243 - 261	29.50 - 38.50	265
	7	68.75	106	266 - 284	249 - 267	30.75 - 39.75	273
	8	70.00	108	272 - 290	255 - 273	32.00 - 41.00	281
	9	71.25	110	278 - 296	261 - 279	33.25 - 42.25	289
	10	72.50	112	284 - 302	267 - 285	34.50 - 43.50	297
	11	73.75	114	291 - 309	274 - 292	35.75 - 44.75	305
	12	75.00	116	297 - 315	280 - 298	37.00 - 46.00	313
	13	76.25	118	303 - 321	286 - 304	38.25 - 47.25	321
	14	77.50	120	309 - 327	292 - 310	39.50 - 48.50	329

Last Features & Modular Amendments

Standard

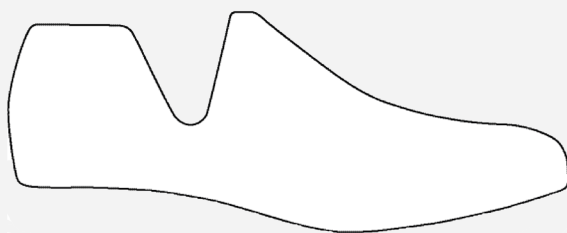


Front view



Heel

Narrow

Narrow last features compared to Standard: slimmer in heel and joint width, shallower in joint depth, toe depth identical to Standard.

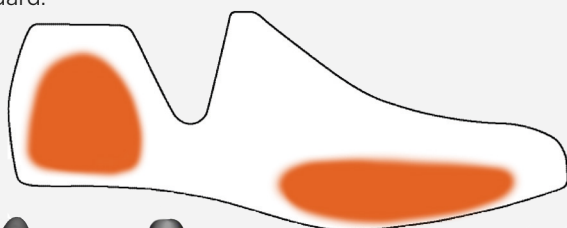


Front view




Heel

Wide

Wide last features compared to Standard: wider at heel and joint width, deeper at instep and circumference, toe depth identical to Standard.

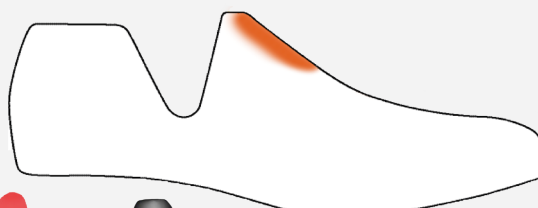


Front view

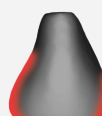

Heel

Deep

Deep last features compared to Standard: identical heel and joint widths to Standard, deeper at the joint and toe depth (defining feature), instep marginally deeper than standard (3mm).

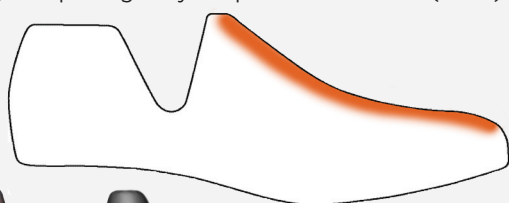

Front view



Heel

Wide & Deep

Wide and Deep - key landmarks wider and deeper than that of the Standard last.



Front view



Heel

Amend Options

Shoes and boots are available without the sole unit for a rough fitting at no additional cost. The following adaptations are also available:

- Raises
- Rocker soles
- Heel flares
- Heel wedges
- Caliper sockets
- T-straps
- Y-straps
- Stirrup straps
- Heel elongations
- Trainer sole (through wedge)

Increase length

Up to 1 or ½ size

Extra depth over toes

Raised by 3mm, 6mm, 9mm or 12mm

Increase lateral joint width

By 3mm, 6mm, or 9mm

Increase depth over throat

By 3mm, 6mm, 9mm or 12mm

Increase instep depth

By 3mm, 6mm, 9mm or 12mm

Straighten back of last

Please note: not available on styles with patterned vamps.

Increase medial joint width

By 3mm, 6mm or 9mm

Bunion pockets

3mm, 6mm, 9mm or 12mm, medial or lateral Moving facings forward 15mm (low) or 25mm (extra low)

