



LIMB-art
Prosthetic Leg Covers



LIMB-art Prosthetic Leg Covers

Allowing patients to stand out and stand proud



#standoutstandproud



CORE

WAVE

VENT

ULTRALIGHT

LIMB-art leg covers combine high-grade Nylon material with cutting-edge design and technology to provide a robust, affordable cosmeses for all users of above and below knee prostheses.

These covers restore balance to the silhouette of the leg and play an instrumental role in raising patient confidence and self-esteem with an optional bright or patterned design, they achieve the LIMB-art mission of 'stand out stand proud'.

Features & Benefits

Lightweight

The cosmetic accessory weighs less than 250g which is balanced across the limb.

Durable & Robust

Made with maximum durability in mind, these covers are constructed from robust high-grade Nylon, known for its strength and toughness whilst retaining its flexibility for ease of application.

Multiple Designs and Custom Options

Offering a range of over 20 solid colour and stock designs, as well as a bespoke design service. LIMB-art covers can be customised to suit any style. If your patient has a specific pattern, photograph or finish in mind, please contact our Customer Service Team.

Quick & Easy to Fit

With its flexible body, the cover can be easily opened and securely strapped to prostheses.



LIMB-art Models

With four different models and over 20 different designs to choose from, there's an option for everyone. See the different models below:

CORE

This elegant smooth design creates a great leg/calf shape under trousers, leggings, or long boots, providing maximum protection for prosthetic limbs and an ideal base for any custom vinyl wrap. Solid colours, stock designs and completely custom options are available for this model.



WAVE

The Wave leg cover has curved ridges that give it a modern, sleek appearance. To protect the prosthetic limb, perforations are kept to a minimum, and the design gives the leg a more natural shape. Choose from solid colours, stock designs or fully bespoke options.



VENT

Designed with the 'super hero' in mind with vented and raised sections. The design creates a great leg shape under leggings, trousers or long boots. Available in a variety of solid colours, stock designs as well as custom finishes.



ULTRALIGHT

The lightest LIMB-art cover features a diamond cut-out pattern that looks great with a contrasting inner colour and provides a natural leg shape under clothing. Finish with a solid colour or a multicoloured stock design for added drama.



Personalisation

Solid Colours

Choose from 9 solid colour finishes, see rear page for more details.

Stock Designs


Select from a range of 'Ready to Go' stock designs, see latest price list for full range of designs.

Custom & Bespoke Designs

Does your patient have a design in mind? With our bespoke design service, we can bring this idea to life.

Custom Vinyl Wrap & Hydrodipped Finish

Using these applications allows covers to be produced with the highest quality finish.

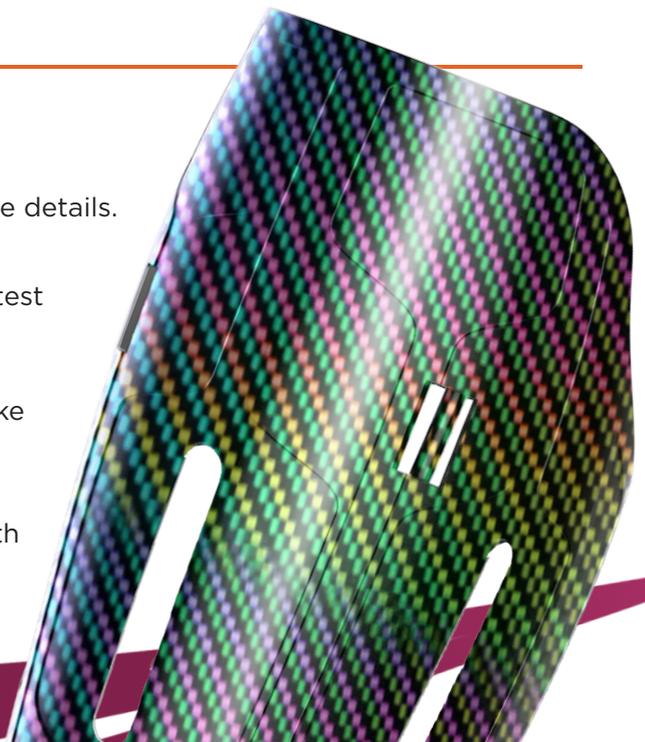

Ordering Information

Example		NC330-370BLUS									
Model ID		Fit		Length		Design				Finish	
Core	NC	Slim (330mm)	330	Short (330mm)	330	Black	BLK	Pink	PINK	Gloss	G
Wave	NW	Standard (360mm)	360	Medium (350mm)	350	Blue	BLU	Yellow	YELL	Satin	S
Vent	NV			Long (370mm)	370	Green	GRN	Maroon	MAR		
Ultralight	NU			†X Long (407mm)	407	Red	RED	White	WHT		
								Teal	TEAL		

†Only Available in Standard Fit

CHOOSE FROM

SOLID COLOURS



STOCK DESIGNS*



BESPOKE DESIGNS

Does your patient have a design idea in mind? With our fully bespoke design service, we can bring this idea to life. Contact us with your design and speak to our friendly Customer Service Team.



*For a complete list of part numbers, designs and custom finishes, please see the latest LIMB-art price list from Steeper

SCAN FOR
PRICELIST

Alternatively, visit our website for a full product overview.

www.steepergroup.com





After receiving his LIMB-art prosthetic leg cover, Tony confidence grew dramatically. Choosing the vent model, thanks to its distinct look, Tony added a bespoke 'stormtrooper' design to reflect his love of mechanics, looking particularly sleek when riding his motorbike.

Manufacturing

LIMB-art leg covers are produced to the highest quality using the latest technology in Computer Aided Design and Manufacture (CAD & CAM) processes.

The durable Nylon material is highly resistant to substances, this means the cover can be submerged within water, detergents, acids or bases without the worry of deteriorating. It's also resilient to cracking and stress, which is ideal for when flexing is required when applying the cover.

These covers are made from 100% recyclable Nylon, which means all covers can be reformed to be made into new products.



Technical Information

How To Measure for Ordering:

LIMB-art covers are available in four different lengths, short, medium, long and X Long (*X Long only available in standard fit*). As well as two different fits, standard and slim. Please follow the steps below to determine which length is required:

Weight of Unit (g)	Material	Build Height (mm)	Warranty
180 - 290	High Grade Nylon	330 - 407	1 Year

	Above the Knee	Below the Knee
Length	Measure from the top of the foot shell to the top of the knee, remove 2cm to leave enough clearance between the top of the foot shell and bottom of the cover.	Measure from the top of the foot shell up to Mid Patella Tendent, remove 2cm to leave enough clearance between the top of the foot shell and the bottom of the cover.
Circumference	Measure the circumference around the widest part of the sound limb/leg socket.	

Steeper Group
Unit 3 Stourton Link
Intermezzo Drive
Leeds
LS10 1DF

Tel: +44 (0) 870 240 4133
Email: customerservices@steepergroup.com
www.steepergroup.com

 **steeper**
Creating life's turning points, together

STPPR280 Issue 1 Jan 2022