


RevoFit Lanyard

Using the BOA M3 dial with 2:1 gearing to provide the mechanical advantage required, ensures that the liner and residual limb can be easily pulled into the socket in correct alignment, negating the need for patients to stand and force their limb into the device.

The fitting of a simple expulsion valve and optional suction sleeve can enable the system to achieve an extremely high level of suspension, allowing the user to engage in high level activities.

Within the Click Medical range, there are a number of products suitable for orthotic application. For further details, please visit our website.

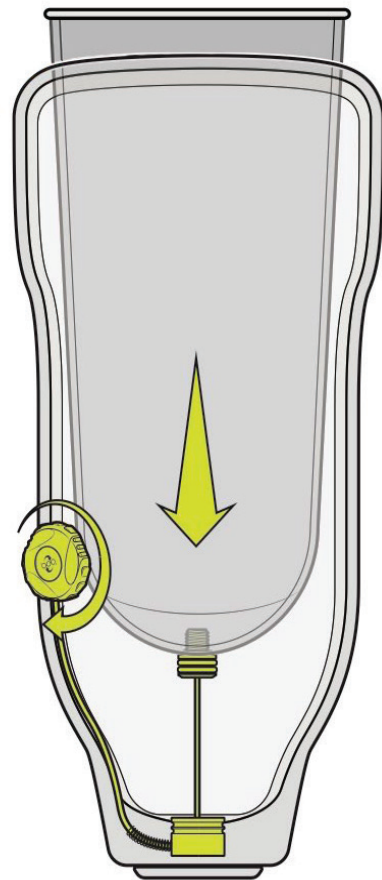


Features and Benefits


- Easy and rapid donning, offering patients the ability to don whilst seated
- The lanyard automatically positions the liner in the correct alignment, eliminating pin lock frustrations
- Dial can be mounted anywhere to accommodate the patient's functional needs and is easily adjustable through clothes and during activity
- High lanyard tension through 2:1 gearing
- Positive lanyard lock
- 'Air lock' created as once drawn into the socket, the lanyard connector seals into the distal adapter
- Easy lock release
- Waterproof system, making it ideal for use in a swimming or bathing prosthesis
- Spectra cord lanyard rated to 100kg loading



Waterproof

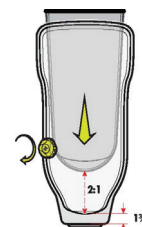

2:1 mechanical advantage

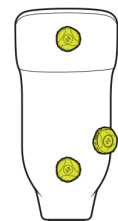

Numerous dial mounting options

Ordering Information

Description	Part Number
RevoFit Lanyard Kit (includes spare lanyard)	NCMPK3000-220-05

Service Parts Ordering Information

Description	Part Number
Lanyard Re-lacing Kit	NCMRP3001-000-10
Lanyard Dial Replacement Set	NCMRP3000-220-05
Lanyard - Lower Limb Ex-Threaded Insert	NCMRP3002-000-00
Lanyard - Upper Ex-Threaded Insert	NCMRP3003-000-05

RevoFit Versa Kit

Particularly beneficial in cases where a patient has limited motion or mobility to reach the dial, the RevoFit Versa allows the reel to be placed anywhere on the socket, with the cord ports exiting the H3 BOA reel parallel to each other.



Ordering Information


Description	Part Number
RevoFit Versa Kit	NCMPK0150-320-05

Service Parts Ordering Information

Description	Part Number
Lamination Dummy	NCMPK0101-320-10
RevoFit Lamination Tubing with Sleeve	NCMPK0080-000-10
Restringing Kit for Direct and Versa	NCMPK0060-000-10

RevoFit Direct

RevoFit Direct applies to cases where it is necessary for the cord ports to be placed directly opposite each other. As with the Versa, the Direct uses the H3 BOA reel.

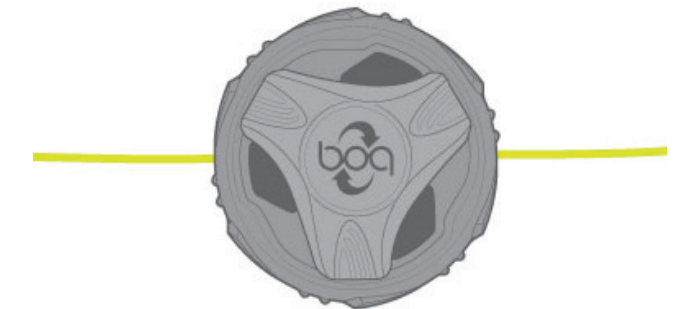

Ordering Information

Description	Part Number
RevoFit Direct Kit	NCMPK0100-320-05

Service Parts Ordering Information

Description	Part Number
Lamination Dummy	NCMPK0151-320-10
RevoFit Lamination Tubing with Sleeve	NCMPK0080-000-10
Restringing Kit for Direct and Versa	NCMPK0060-000-10