


RevoFit2™

The RevoFit2 enables built-in user adjustability anywhere on a socket. An adjustable multi-panel socket design provides global and targeted compression, with every 'click' of the dial providing 1mm of movement.


For use on every lamination type, the RevoFit2 uses a high-power dial with a 4:1 mechanical advantage. It has been ideally formulated to withstand forces and repetitive stress found in an adjustable socket, incorporating fatigue-resistant Spectra lace and low-friction Teflon tubing.

RevoFit2 incorporates both configurations as shown:

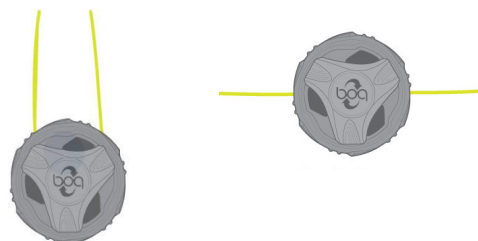

Features and Benefits

- Ease of releasing tension with the system allows for rapid donning and doffing for bulbous residual limb ends or large volume fluctuations
- Micro/macro adjustments provide 1mm of movement for compression, and can adjust for up to 10ply of sock changes.
- Dial can be mounted anywhere to accommodate the patient's functional needs and is easily adjustable through clothes or during activity
- Offers patients the ability to don whilst seated
- Uses strong and fatigue-resistant Spectra lace
- Low-friction Teflon tubing supplied
- Waterproof system, making it ideal for use in a swimming or bathing prosthesis
- Rated to 100kg loading. For a higher load, a second system can be employed to share excess weight across multiple systems

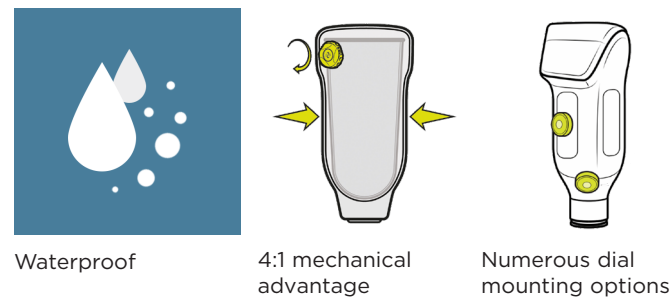

Ordering Information

Description	Part Number
RevoFit2 Lamination Kit	NCMPK2000-320-05

Service Parts Ordering Information

Description	Part Number
Re-lacing Pack	NCMPK0060-000-10
Replacement Dial & Lace Pack	NCMRP0100-320-05
RevoFit Bulk Spectra Lace	NCML0K014
Bulk Tubing and Sheathing	NCMPK0080-000-10
Dummy with Silicone Overmold	NCMPK2000-000-05
Bulk Tubing with Sheathing	NCMPK0080-000-10

Structural and Socket Components - Click Medical - RevoFit2™

RevoFit2™ Diagnostic Kit

The RevoFit2 diagnostic kit allows for socket adjustment early on in development stages. This kit can be easily applied to diagnostic sockets after their successful initial fitting.



Features and Benefits

- User can influence during fitting the location of the dial, ensuring it is within easy reach or hidden beneath clothing
- Designed to be added after a diagnostic socket has been produced and the initial socket fitting has taken place, for easy and simple fabrication
- Introducing the kit at socket stage allows flexibility in socket design and configuration of compression windows
- Multiple kits can be applied to the same socket based on the user's needs

Ordering Information

Description	Part Number
RevoFit2 Diagnostic Kit	NCMPK2150-320-05

Service Parts Ordering Information

Description	Part Number
Re-lacing Pack	NCMPK0060-000-10
Replacement Dial & Lace Pack	NCMRP0100-320-05
RevoFit Bulk Spectra Lace	NCML0K014

Service Parts Information

Parts	Details
RevoFit Re-lacing Pack	Includes: 2.1m of RevoFit Spectra Lace (capable of supporting 226kg), Lace Feeder Tool, and T6 Tool. Replacing lace every 6 months or at sign of abrasion is recommended.
RevoFit Replacement Dial & Lace Pack	Includes: RevoFit2 Dial, 2.1m of RevoFit Lace (capable of supporting 226kg), Lace Feeder Tool, and T6 Tool. Replacing lace at time of dial replacement is recommended.
RevoFit Bulk Spectra Lace	10m of RevoFit Spectra Lace.
Bulk Tubing with Sheathing	1.8m of low-friction tubing and Teflon sheathing.