



Lindhe Xtend
Xtend Connect

 **steeper**
Creating life's turning points, together



Xtend Connect by Lindhe Xtend



The Xtend Connect gives prosthetic users more flexibility by allowing them to easily switch between a prosthesis or socket. The component part can be fitted anywhere with a variety of industry standard connections. A simple wheel lock mechanism and safety button give true flexibility.

The Xtend Connect can be used to connect different prosthetic components, depending on the lifestyle and usage needs of the user. Users can now wear the same prosthetic socket across multiple limbs. Dependent upon activity, users can simply swap their prosthetic foot, whilst keeping the same socket on.

Suitable for both transtibial and transfemoral applications, the system is ideal for quickly disconnecting the prosthetic foot to make dressing easier, or allowing the same limb to be worn with feet set to varied heel heights - without the need to change footwear.

Features and Benefits:

- Accommodates for weights up to 150kg*
- Quickly exchange feet whilst keeping the same socket, and maintaining optimal alignment of prosthesis
- Change sockets whilst keeping limb builds intact
- Designed to be used with industry standard 30mm systems
- Greater flexibility when dressing - such as removing feet to pull on trousers
- Waterproof in fresh water, sea water and chlorinated water with appropriate care**
- Less baggage - users can pack their activity specific prosthesis without having to pack the entire limb build
- Safety button ensures secure connection from the first step

* For activities reaching K4 levels, patient weight is limited at 100kg

** Rinse immediately with fresh water and dry following exposure to salt water or chlorinated water.



Xtend Connect Components

The Xtend Connect can be rotated in any direction and freely combined with compatible prosthetic components. It is fitted in-line and can be fitted above or below the prosthetic knee, with a variety of connection options to choose from. Available in two wheel colour options: Black or Orange.

An Xtend Connect set consists of two components:

Quick Fit-Connect


The Quick-Fit Connect unit features an easy-to-use locking wheel and safety button for connecting and disconnecting.

Base Unit

The Base Units easily fit into the Quick-Fit Connect and are available as a Pyramid or Pyramid Receiver.



Tips for assembling Xtend Connect

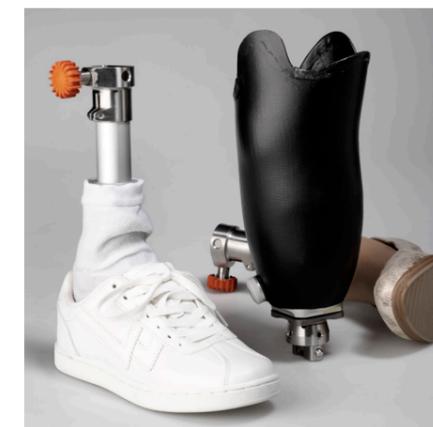


Changing feet

Place the Xtend Connect as high up as possible on the connecting tube, to allow for prosthetic feet of varying heights.



Changing trousers

The lower the Xtend Connect is on the prosthetic limb, the easier it is for the user to pull on tight trousers.

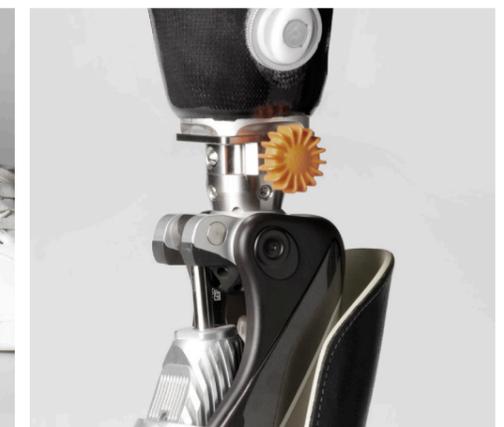

Changing the knee joint

To ensure correct alignment every time, use the female pyramid as a connecting component to the knee joint.

Xtend Connect Technical Information

Product Specification


Material	Stainless Steel	Activity Level	K3 - K4*
Weight Limit	150kg	Warranty	2 Years



* For activities reaching K4 levels, patient weight is limited at 100kg

Xtend Connect Component Ordering Information

Xtend Connect Components				
Part Number	Description	Wheel Colour	Weight Limit	
1	NA100-Q30-150-O	Quick Fit Unit, Tube Clamp 30 mm	Orange	150kg
	NA100-Q30-150	Quick Fit Unit, Tube Clamp 30 mm	Black	150kg
2	NA100-Q4-150-O	Quick Fit Unit, 4-Hole Adapter	Orange	150kg
	NA100-Q4-150	Quick Fit Unit, 4-Hole Adapter	Black	150kg
3	NA100-QF-150-O	Quick Fit Unit, Female Pyramid Adapter	Orange	150kg
	NA100-QF-150	Quick Fit Unit, Female Pyramid Adapter	Black	150kg
4	NA100-QP-150-O	Quick Fit Unit, M36 Threaded Device	Orange	150kg
	NA100-QP-150	Quick Fit Unit, M36 Threaded Device	Black	150kg
5	NA100-BF-150	Base Unit with Female Pyramid Receiver	-	150kg
6	NA100-BM-150	Base Unit with Male Pyramid Adapter	-	150kg

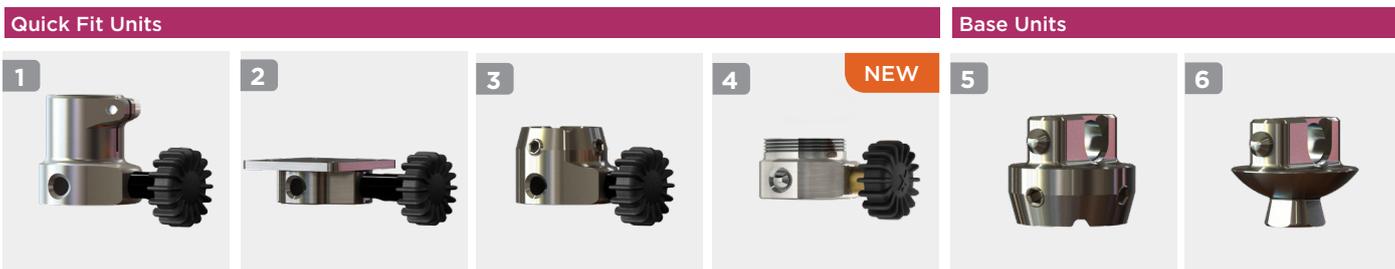


All parts are supplied with a two year warranty.

For more information on the Xtend Connect range contact your local Product Manager, or to place an order please contact our Customer Services team via customerservices@steepergroup.com or 0870 240 4133.

INTERACTIVE

Lindhe Xtend Connect Configurator

Use the Xtend Connect component selection tool to ensure the most effective combination is selected.

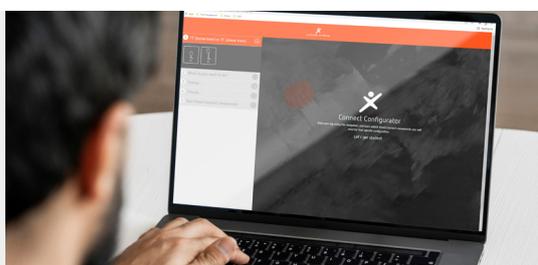

SCAN TO TRY

Steeper Group
Unit 3 Stourton Link
Intermezzo Drive
Leeds, LS10 1DF

Tel: +44 (0) 870 240 4133
Email: customerservices@steepergroup.com



www.steepergroup.com
STPPR097 May 2022