

LIMB-art leg covers combine high-grade Nylon material with cutting-edge design and technology to provide a robust, affordable cosmeses for all users of above and below knee prostheses.

These covers restore balance to the silhouette of the leg and play an instrumental role in raising patient confidence and self-esteem. With optional bright or patterned designs, they achieve the LIMB-art mission of 'stand out and stand proud'.



**\*Waterproof note:** Following exposure to chlorinated water, rinse immediately with fresh water and dry.

## Features and Benefits

- This prosthetic accessory weighs less than 250g which is balanced across the limb
- Made with maximum durability in mind, these covers are constructed from robust high-grade Nylon, known for its strength and toughness whilst retaining flexibility for ease of application
- Offering a range of over 20 solid colour and stock designs, as well as a bespoke design service. LIMB-art covers can be customised to suit any style. If your patient has a specific pattern, photograph or finish in mind, please contact our Customer Service Team
- With its flexible body, the cover can be easily opened and securely strapped to prostheses

## Models



**CORE**  
An elegant smooth design which gives an excellent leg/calf shape under trousers, leggings or long boots.

**WAVE**  
The Wave leg cover features curved ridges for a modern, sleek appearance.



**VENT**  
Designed with the 'super hero' in mind with vented and raised sections.

**ULTRALIGHT**  
The most lightweight LIMB-art cover features a diamond cutout pattern that is particularly popular with a contrasting inner colour.

## Ordering Information


Ordering Example - NC330-370BLUS											
Part Number	Fit	Length	Design	Finish							
NC	330	370	BLU	S							
Model ID	Fit	Length	Design	Finish							
Core	NC	Slim (330mm)	330	Short (330mm)	330	Black	BLK	Pink	PINK	Gloss	G
Wave	NW	Standard (360mm)	360	Medium (350mm)	350	Blue	BLU	Yellow	YELL	Satin	S
Vent	NV			Long (370mm)	370	Green	GRN	Maroon	MAR		
Ultralight	NU			*X Long (407mm)	407	Red	RED	White	WHT		
								Teal	TEAL		

\*Only Available in Standard Fit

### CHOOSE FROM



**BESPOKE DESIGNS** Does your patient have a design idea in mind? With our fully bespoke design service, we can bring this idea to life. Speak to our Customer Services Team and reference **NLIMB-ART-CUSTOM** along with your design.

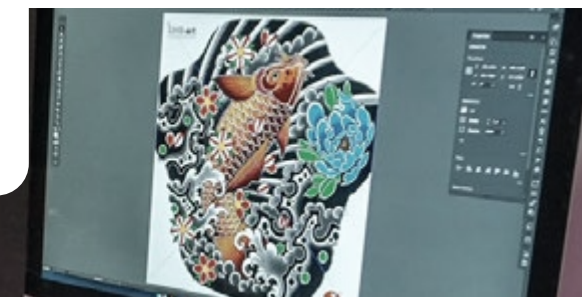



Technical Information

Weight of Unit (g)	Material	Build Height (mm)	Warranty
180 - 290	High Grade Nylon	330 - 407	1 Year

How To Measure for Ordering:

LIMB-art covers are available in four different lengths, short, medium, long and X Long (X Long only available in standard fit) as well as two different fits, standard and slim. Please follow the steps below to determine which length is required:

Amputation Level	Length	Circumference
Above Knee	Measure from the top of the foot shell to the top of the knee, remove 2cm to leave enough clearance between the top of the foot shell and bottom of the cover.	Measure the circumference around the widest part of the sound limb/leg socket.
Below Knee	Measure from the top of the foot shell up to Mid Patella Tendent, remove 2cm to leave enough clearance between the top of the foot shell and the bottom of the cover.	

For more information or to place an order, please contact our Customer Services Team.