

Accent® IP

Providing 5cm of heel height adjustment through a simple push of a button for users to control to fit their daily needs, Accent IP offers a maintenance free option for users.

A cosmetically appealing foot, the Accent IP is provided with a contoured, soft shin fairing that locks into the shell and prevents cosmetic buckling when the foot is adjusted at various heights.

This version comes with an integrated pyramid for easy connection and has dynamic response and compliance to meet with a vast array of low to moderate activities.

**CES CPI F04:** Use this code to read the Accent IP clinical evaluation summary on the Steeper website.

Features and Benefits

- User adjustable heel height with 5cm range and 13 increments of adjustment
- Easy to operate height adjustment button located in the ankle
- Ankle fairing provides superior cosmetic finishing
- Maintenance free



Height adjustment button

Cosmetics

The Accent IP's bonded sandal toe foot shell is reinforced with Impact Dispersion Technology™. The Ankle Fairing interlocks with the foot shell and conforms to both Endo and dynamic pylons. The foot shell and ankle fairing are anatomically designed for use with a skin to produce a natural-like finish.



Sandal toe



Ankle fairing

Ordering Information

Trial Period - This foot has a 30 day trial period.

Impact Level Descriptions

**Low:** Daily activities include mostly level ground walking, moving around in the home and the community.


**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

Accent IP foot includes Ankle Fairing and Foot Shell. Wider versions from sizes 23-28 are available - see Order Forms - page 37 for additional ordering information.

Product Specification

Mounting	Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty	Foot Shell
Endo (IP)	100kg		532g	21-24	8.5	2 years	 Caucasian or Brown Narrow or wide width
Endo (IP)	114kg		603g	25-28	9.9	2 years	 Caucasian or Brown Narrow or wide width

Technical Information

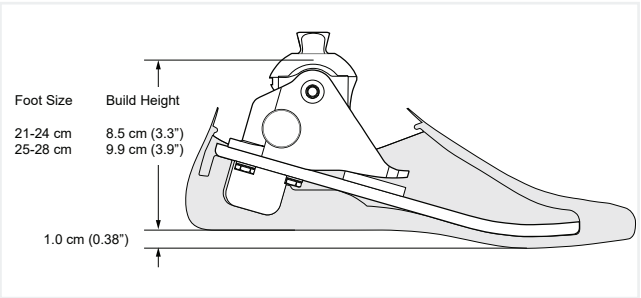

Feet - College Park - Accent IP

Shelltread - Ordering Information

Shelltread	Size
NACT  N	21-28cm
NACT  N  W	23-28cm

Shelltread Ordering Guide  
 = "L" or "R" for Left or Right  
 = Foot size in cm (21-28)  
W = Wide