

The Tribute foot offers a cost-effective option that provides stable footing and dependable performance on a variety of terrain, making it ideal for low to moderate impact individuals. Precisely gait matched, the multi-axial design and full-length toe lever provide the user with improved control and stability.

The Tribute provides anatomical plantar flexion, dorsiflexion, transverse rotation, and inversion/eversion through true multi-axial movement in all three planes, allowing for ease of movement, excellent ground compliance, and comfort.

Stride Control™ allows optional customisation of the plantar flexion/dorsiflexion characteristics with a simple adjustment of the posterior bumper.


CES CPI F02: Use this code to read the Tribute clinical evaluation summary on the Steeper website.

Features and Benefits

- Tri-axial technology to adapt to uneven terrain
- Full range of motion whilst maintaining ground contact
- Adjustable Stride Control™
- Removable foot shell
- Easy assembly and disassembly with one pin



Stride Control™ adjuster

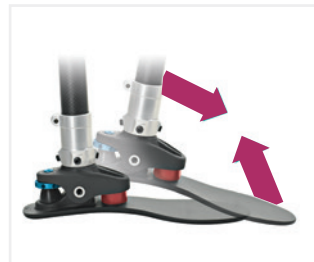

Dynamic response

Cosmetics



CAP



Removable foot shell with sandal toe

Note: Please see the following page for Tribute spare parts ordering information.



Ordering Information

Trial Period - This foot has a 30 day trial period.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around the home and the community.

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

Ordering Example: NTB-LCS-27W-G2

Part ID	Side	Shell Colour	Mounting	Toe	Size	Width**	Firmness Category
NTB	L	<input type="text"/>	EN	<input type="text"/>	<input type="text"/>	<input type="text"/>	G <input type="text"/>
NTB	R	<input type="text"/>	EN	<input type="text"/>	<input type="text"/>	<input type="text"/>	G <input type="text"/>

Includes foot shell, CPI Spectra Sock and cosmetic attachment plate.

Caucasian	C	Regular Toe	Leave blank
Tan	T	Sandal Toe	S
Brown	B		
Jet Black*	J		

*Regular Toe only

22-30cm

Narrow

Wide 24-30cm only

** Leave Blank When Ordering a Sandal Toe

N

W

1-4


Firmness Category

Weight (kg)	0-45	46-54	55-63	64-72	73-81	82-91	92-100
Size (cm)	21-30						
Low Impact	1	1	2	2	3	3	4
Moderate Impact	1	2	2	2	3	4	4

Product Specification

Mount	Weight Limit	Sizes	Impact Level	Build Height	Weight of Unit	Warranty	Foot Shell
Endo (IP)	100kg	21-30cm		6.1cm	502g Size 26 inc. foot shell	2 years Foot shell - 3 months	 Caucasian, Tan, Brown or Jet Black Narrow or Wide width

Technical Information

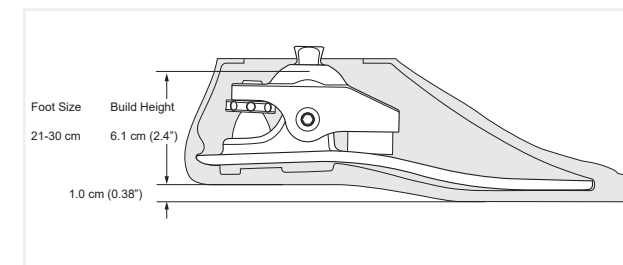

Spare Parts Ordering Information

Foot Shell	Size	Toe
NV-S ■ ● -▲N	21-30	Regular, Narrow
NV-S ■ ● -▲W	24-30	Regular, Wide
NCS-S ■ ● -▲	21-30	Sandal



- = "L" or "R" for Left or Right
- = Colour (C = Caucasian, T = Tan, B = Brown, J = Jet Black*)
- ▲ = Foot size in cm (21-30)
- *Regular toe only

Shellread	Size (cm)
NV-T ■ N-▲N	21-30 (Narrow)
NV-T ■ N-▲W	24-30 (Wide)



- = "L" or "R" for Left or Right
- ▲ = Foot size in cm (21-30)
- *Regular toe only

CAP	Size
NV-C ■ ● -2124	21-24
NV-C ■ ● -2526	25-26
NV-C ■ ● -2730	27-30



- = "L" or "R" for Left or Right
- = Colour (C = Caucasian, T = Tan, B = Brown)

CPI Spectra Sock	Size
NV-SS-2124	21-24
NV-SS-2526	25-26
NV-SS-2730	27-30



Endo Ankle	Description
NTB-ENA-2124	Endo Ankle, 21-24cm
NTB-ENA-2526	Endo Ankle, 25-26cm
NTC-ENA-2730	Endo Ankle, 27-30cm

Gait Match					
Firmness Category	Front Bumper	Ankle Bushings	Rear Bumper	*Extra Front Bumpers	*Extra Rear Bumpers
1	S	M	S	M	M
2	M	M	M	S F	S F
3	M	M	F	S F	M XF
4	F	F	XF	M	F

*Foot orders include extra front and rear bumpers (one category firmer, one category softer).

Bumpers & Bushings	Description	Bumper Colour
NV-FBS-●	Front Bumper, Soft	S
NV-FBM-●	Front Bumper, Medium	M
NV-FBF-●	Front Bumper, Firm	F
NV-RBS-▲	Rear Bumper, Soft	S
NV-RBM-▲	Rear Bumper, Medium	M
NV-RBF-▲	Rear Bumper, Firm	F
NV-RBX-▲	Rear Bumper, Extra Firm	XF
NTS-ABS	Ankle Bushing, Soft, 21-30cm	S
NTS-ABM	Ankle Bushing, Medium, 21-30cm	M
NTS-ABF	Ankle Bushing, Firm, 21-30cm	F
NV-BK-2130	Bumper Kit	

- S Soft M Medium F Firm XF Extra Firm
- = Front Bumper (FB) ▲ = Rear Bumper (RB)

Pocket Size	Foot Size		
	21-24cm	25-26cm	27-30cm
● (FB)	C	B	A
▲ (RB)	D	C	B

Assemblies	Description
NTB-AK-2124	Assembly Kit, 21-24cm**
NTB-AK-2526	Assembly Kit, 25-26cm**
NTB-AK-2730	Assembly Kit, 27-30cm**
NTB-FH ■ -▲	Foreheel
NTB-FHS-▲	Foreheel for sandal toe foot shell

- = "L" or "R" for Left or Right
- ▲ = Foot size in cm (21-30)

- **Assembly Kit Contents:
- 1 x NTL-LUB (Tru lube)
 - 1 x NV-FBS-(size) - Front Bumper, Soft
 - 1 x NV-FBM-(size) - Front Bumper, Medium
 - 1 x NV-FBF-(size) - Front Bumper, Firm
 - 1 x NCP-AP-L (Axial Pin, long)
 - 1 x NV-RBS-(size) - Rear Bumper, Soft
 - 1 x NV-RBM-(size) - Rear Bumper, Medium
 - 1 x NV-RBF-(size) - Rear Bumper, Firm
 - 1 x NCP-APS (Axial Pin screw)

