

# Tribute®

The Tribute offers a cost-effective option, providing stable footing and dependable performance on varied terrain and is ideal for low to moderate impact individuals. Precisely gait matched, the multi-axial design and full-length toe lever provide the user with improved control and stability.

With true multi-axial movement occurring on all three planes, the Tribute provides anatomical plantar flexion, dorsiflexion, transverse rotation and inversion/eversion - allowing for ease of movement, excellent ground compliance and comfort.

Stride Control™ allows optional customisation of the plantarflexion/dorsiflexion characteristics through a simple adjustment of the posterior bumper.

**CES CPI F02:** Use this code to read the Tribute clinical evaluation summary on the Steeper website.

## Features and Benefits

- Tri-axial technology to adapt to uneven terrain
- Full range of motion whilst maintaining ground contact
- Adjustable Stride Control™
- Removable foot shell
- Easy assembly and disassembly with one pin



Stride Control™ adjuster



Dynamic response



Transverse rotation inversion-eversion

## Cosmetics



CAP



Removable foot shell with sandal toe



## Ordering Information

Trial Period - This foot has a 30 day trial period.

### Impact Level Descriptions

**Low:** Daily activities include mostly level ground walking, moving around the home and the community.

**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

### Part Number

Left	NTB-L	EN	◆	●	●	-G	▲
Right	NTB-R	EN	◆	●	●	-G	▲

Includes foot shell, CPI Spectra Sock and cosmetic attachment plate.

- - Colour of Cosmesis (C = Caucasian, T = Tan, B = Brown, J = Jet Black - regular toe only)
- ◆ - Toe style (S = Sandal Toe, Regular leave blank)
- - Foot Size in cm (21-30)
- - Width (N = Narrow\*, W = Wide\*\*, Standard leave blank)
- ▲ - Firmness Category (1-4)

\*Not available in Sandal Toe

\*\*For sizes 24-30cm only. Not available in Sandal Toe. Wide option suggested for male patients.

### Firmness Category


Weight (kg)	Sizes 21-30cm	
	Low Impact	Moderate Impact
0-45	1	1
46-54	1	2
55-63	2	2
64-72	2	2
73-81	3	3
82-91	3	4
92-100	4	4

## Product Specification

Mount	Weight Limit	Sizes	Impact Level	Build Height	Weight of Unit	Warranty	Foot Shell
Endo (IP)	100kg	21-30cm		6.1cm	502g (Size 26 inc. foot shell)	2 years Foot shell - 3 months	Caucasian, Tan, Brown or Jet Black Standard, Narrow or Wide width

## Feet - College Park - Tribute

## Technical Information

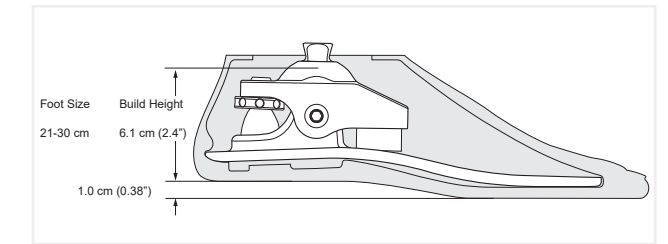

### Gait Match

Firmness Category	Front Bumper	Ankle Bushings	Rear Bumper	*Extra Front Bumpers	*Extra Rear Bumpers
1	S	M	S	M	M
2	M	M	M	S F	S F
3	M	M	F	S F	M XF
4	F	F	XF	M	F

\*Foot orders include extra front and rear bumpers (one category firmer, one category softer).

Note: Please see the following page for Tribute spare parts ordering information.

Spare Parts Ordering Information



Foot Shell	Size	Toe
NV-S ■●-▲N	21-30	Regular, Narrow
NV-S ■●-▲W	24-30	Regular, Wide
NCS-S ■●-▲	21-30	Sandal



Shelltread	Size (cm)
NV-T ■N-▲N	21-30 (Narrow)
NV-T ■N-▲W	24-30 (Wide)



Neutral

- = "L" or "R" for Left or Right
- = Colour (C = Caucasian, T = Tan, B = Brown, J = Jet Black\*)
- ▲ = Foot size in cm (21-30)
- \*Regular toe only



CAP	Size
NV-C ■●-2124	21-24
NV-C ■●-2526	25-26
NV-C ■●-2730	27-30

- = "L" or "R" for Left or Right
- = Colour (C = Caucasian, T = Tan, B = Brown)



CPI Spectra Sock	Size
NV-SS-2124	21-24
NV-SS-2526	25-26
NV-SS-2730	27-30

Endo Ankle	Description
NTB-ENA-2124	Endo Ankle, 21-24cm
NTB-ENA-2526	Endo Ankle, 25-26cm
NTC-ENA-2730	Endo Ankle, 27-30cm

Bumpers & Bushings	Description
NV-FBS-◆†	Front Bumper, Soft
NV-FBM-◆†	Front Bumper, Medium
NV-FBF-◆†	Front Bumper, Firm
NV-RBS-◆*	Rear Bumper, Soft
NV-RBM-◆*	Rear Bumper, Medium
NV-RBF-◆*	Rear Bumper, Firm
NV-RBX-◆*	Rear Bumper, Extra Firm
NTS-ABS	Ankle Bushing, Soft, 21-30cm
NTS-ABM	Ankle Bushing, Medium, 21-30cm
NTS-ABF	Ankle Bushing, Firm, 21-30cm
NV-BK-2130	Bumper Kit

◆ = Pocket Size:

	Foot Size		
	21-24cm	25-26cm	27-30cm
† Code	C	B	A
* Code	D	C	B

Assemblies	Description
NTB-AK-2124	Assembly Kit, 21-24cm**
NTB-AK-2526	Assembly Kit, 25-26cm**
NTB-AK-2730	Assembly Kit, 27-30cm**
NTB-FH ■-▲	Foreheel
NTB-FHS ▲	Foreheel for sandal toe foot shell

- = "L" or "R" for Left or Right
- ▲ = Foot size in cm (21-30)

- \*\*Assembly Kit Contents:
- 1 x NTL-LUB (Tru lube)
  - 1 x NV-FBS-(size) - Front Bumper, Soft
  - 1 x NV-FBM-(size) - Front Bumper, Medium
  - 1 x NV-FBF-(size) - Front Bumper, Firm
  - 1 x NCP-AP-L (Axial Pin, long)
  - 1 x NV-RBS-(size) - Rear Bumper, Soft
  - 1 x NV-RBM-(size) - Rear Bumper, Medium
  - 1 x NV-RBF-(size) - Rear Bumper, Firm
  - 1 x NCP-APS (Axial Pin screw)

