

The durable Truper foot was designed to combine stability and dynamic response to take on the toughest of childhood challenges. In use, the foot stores and releases energy for a smooth and controlled transition from standing still to running. The inversion/eversion and superior transverse rotation provides exceptional ground compliance; improving performance and stability, whether walking or running.

The foot is available in two size ranges, which allows the foot size to be increased simply by purchasing the next size foot shell and sole plate, providing a cost-effective and flexible solution for a growing child.

Features and Benefits


- Versatile finishing capability
- Sandal toe design allows for a choice of shoes
- Reduces fatigue and protects the residual limb
- Adjustable Stride Control™ to customise child's gait
- Sealing boot interlocks with the foot shell to prevent any debris from entering the foot



Dynamic response



Inversion/eversion

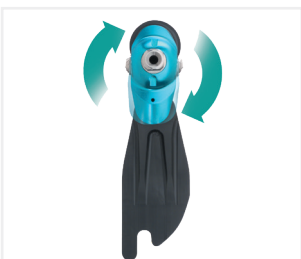

Transverse rotation



Stride Control™ Adjuster

Cosmetics



Sandal toe foot shell



Sealing boot



Exo Block

Ordering Information

There is a 30-day trial period on this foot.

Impact Level Descriptions

High: Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

Part Number			
Left	NTP-L	Colour	Firmness
Right	NTP-R	Colour	Firmness

- - Colour of Cosmesis (C = Caucasian, T = Tan, B = Brown)
- - Mounting (EN = Endo, ALX = Exo)
- - Foot Size in cm (16-21)
- ▲ - Firmness Category (1-3)

Includes foot shell, CPI Spectra Sock and sealing boot.

Firmness Category		
Weight (kg)	Sizes 16-21cm	Sizes 19-21cm
	High Impact	High Impact
0-22	1	-
23-33	2	-
34-45	3	-
46-60	3	3

Product Specification

Mounting	Sizes	Weight Limit	Impact Level	Weight	Build Height	Warranty	Foot Shell
Endo (IP)	16-21cm	45kg (16-18cm)		215g (size 16)	5.3-6.1cm (Endo IP)	1 year Foot shell - 6 months	C T B Caucasian, Tan or Brown Standard width
Exo (ALX)		60kg (19-21cm)		314g (size 21)	11-11.8cm (Endo ALX)		

Foot Shell Ordering Information


Foot Shell	Size (cm)	Shell
NTP-S ■ ● -▲	16-21	S
NTP-XS ■ ● -▲	16-21	XS*

*Exo only, add 1 day for processing.

Shelltread	Size (cm)
NTP-T ■ N-▲	16-21
	Neutral

- = "L" or "R" for Left or Right
- = Colour (C = Caucasian, B = Brown, T = Tan)
- ▲ = Foot size in cm (16-21)

Technical Information

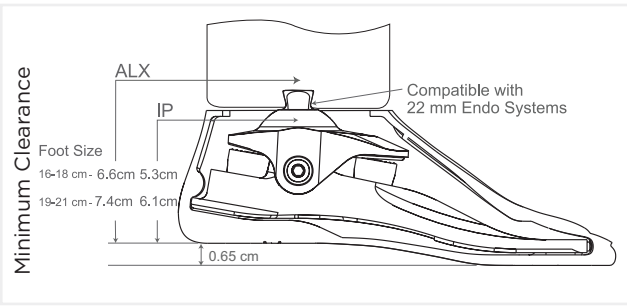


Gait Match			
Firmness Category	Front Bumper	Ankle Bushings	Rear Bumper
1	S	One Firmness	S
2	M		M
3	F		F

Note: Please see the following page for Truper spare parts ordering information.




Spare Parts Ordering Information

Bumpers & Bushings	Description
NTP-BS-1618	Bumper, Soft, 16-18cm
NTP-BS-1921	Bumper, Soft, 19-21cm
NTP-BM-1618	Bumper, Medium, 16-18cm
NTP-BM-1921	Bumper, Medium, 19-21cm
NTP-BF-1618	Bumper, Firm, 16-18cm
NTP-BF-1921	Bumper, Firm, 19-21cm
NTP-ABM-1618	Ankle Bushing, 16-18cm
NTP-ABM-1921	Ankle Bushing, 19-21cm

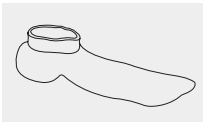


CPI Spectra Sock	Size (cm)
NTP-SS-1621	16-21
NTP-SSX-1618	16-18 (Exo)

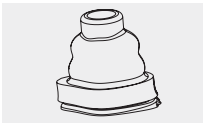


Sealing Boot	Size (cm)
NTP-SB ■●-1618	16-18
NTP-SB ■●-1921	19-21

■ = "L" or "R" for Left or Right
● = Colour (C = Caucasian, B = Brown, T = Tan)

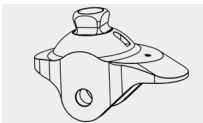


Endo Ankle	Size
NTP-ENA-1618	Endo Ankle, 16-18cm
NTP-ENA-1921	Endo Ankle, 19-21cm

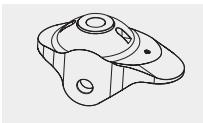


Exo Ankle	Size
NTP-ALXA-1618	Exo Ankle, 16-18cm
NTP-ALXA-1921	Exo Ankle, 19-21cm



Foreheel	Description
NTP-FHA ■-▲	Foreheel
■ = "L" or "R" for Left or Right ▲ = Foot size in cm (16-21)	

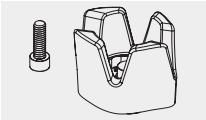

Includes Mounting Bolt

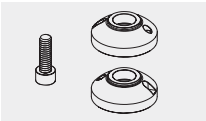



Includes 2 Growth Plates, Mounting Bolt and Finishing Ring

Exo Block Kit	Description
NTP-XBK ■-1618	Alignable Exo Block Kit, 16-18cm
NTP-XBK ■-1921	Alignable Exo Block Kit, 19-21cm

Part Number	Description
NTP-GK ■ C-1618	Growth Plate Kit, Caucasian, 16-18cm
NTP-GK ■ B-1618	Growth Plate Kit, Brown, 16-18cm
NTP-GK ■ C-1921	Growth Plate Kit, Caucasian, 19-21cm
NTP-GK ■ B-1921	Growth Plate Kit, Brown, 19-21cm
NTP-GP-1618	Growth Plate, 16-18cm
NTP-GP-1921	Growth Plate, 19-21cm

■ = "L" or "R" for Left or Right

Growth Range

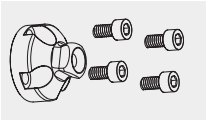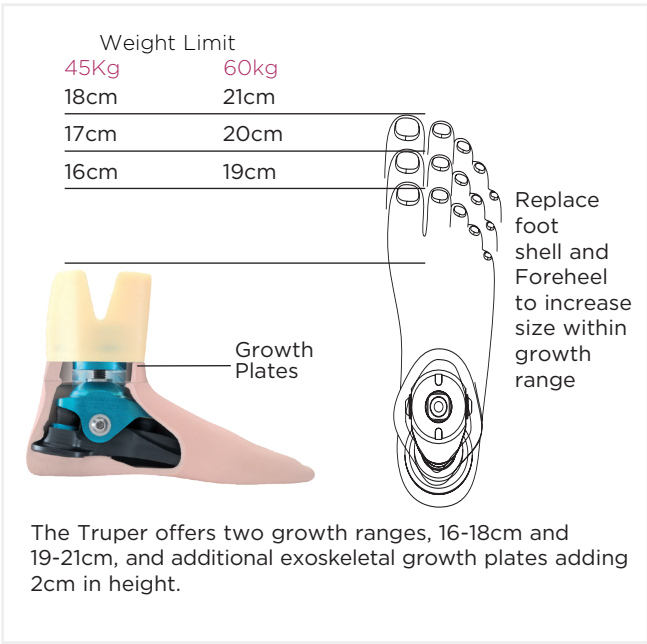

Part Number	Description
NTP-EXPK-22	Exo Pyramid Tool Kit, 22mm* *For temporary endoskeletal alignment during test fittings

Exo Mounting Blocks	Description	Application
NTP-SH-5X10	Exo Mounting Bolt, M5 x 10mm	Standard exo mounting bolt to attach an exo foot to the exo block
NTP-SH-8X16	Exo Mounting Bolt, M8 x 16mm	Used if one growth plate is being used with the exo block
NTP-SH-8X30	Exo Mounting Bolt, M8 x 30mm	Used if there are two growth plates being used with the exo block

Finishing Rings		
Left	Right	Description
NTP-FRLC-1618	NTP-FRRC-1618	Caucasian, 16-18cm
NTP-FRLC-1921	NTP-FRRC-1921	Caucasian, 19-21cm
NTP-FRLB-1618	NTP-FRRB-1618	Brown, 16-18cm
NTP-FRLB-1921	NTP-FRRB-1921	Brown, 19-21cm

Hardware	Description
NTP-AP-1618	Axial Pin, 16-18cm
NTP-AP-1921	Axial Pin, 19-21cm
NCP-APS	Axial Pin Screw
NTP-RBK	Pyramid Tool Retaining Bolt Kit (4pc) M5 x 10mm
NTP-DP-4X12	Dowel Pin, M4 x 12mm
NTP-DP-4X8	Dowel Pin, M4 x 8mm

