


The Venture foot has a full-length toe lever for a more dynamic response for active users.

Its two-stage response enables the dorsiflexion bumper to compress providing comfort during walking; this compression is controlled to engage the composite keel which improves the foot's response to high-impact activities. The heel bumper can be fine-tuned to achieve the ideal heel strike compression.

Designed for a demanding lifestyle, whether being used in the gym or the office, the precise engineering of the Venture generates a dynamic performance without sacrificing comfort.

Features and Benefits

- Tri-axial technology easily adapts to any terrain
- Energy storing to reduce fatigue
- Adjustable Stride Control™ - using 4mm hex key to customise the user's gait timing, even with the foot shell on
- Smooth transition from high to low activities

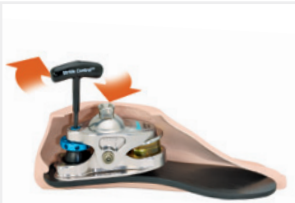

Stride Control™ adjuster



Energy storing



Transverse rotation inversion-eversion

Cosmetics



CAP



Solid toe foot shell

Note: Please see the following page for Venture spare parts ordering information.

Ordering Information

Trial Period - This foot has a 30 day trial period.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around the home and the community.

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

High: Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.


| Ordering Example: NV-LCEN-24W-G4 | | | | | | | | | |
|--|---|-----------|--------------|----------|-----|---------|---------------|-------------------|-----|
| Part ID | | Side | Shell Colour | Mounting | | Size | Width* | Firmness Category | |
| NV | - | L | | | - | | | G | |
| NV | - | R | | | - | | | G | |
| Includes foot shell, CPI Spectra Sock and cosmetic attachment plate. | | Caucasian | C | Endo | EN | 21-30cm | Narrow | N | 1-5 |
| | | Tan | T | Exo | ALX | | Wide | W | |
| | | Brown | B | | | | *24-30cm only | | |
| | | Jet Black | J | | | | | | |

| Firmness Category | | | | | | | | | |
|-------------------|-------|-------|-------|-------|-------|-------|--------|---------|---------|
| Weight (kg) | 0-45 | 46-54 | 55-63 | 64-72 | 73-81 | 82-91 | 92-100 | 101-113 | 114-125 |
| Size (cm) | 21-24 | | | | | | | | |
| Low Impact | 1 | 1 | 2 | 3 | 3 | 4 | 4 | 5 | - |
| Moderate Impact | 1 | 2 | 3 | 3 | 4 | 4 | 5 | 5 | - |
| High Impact | 2 | 3 | 3 | 3 | 4 | 4 | 5 | 5 | - |
| Size (cm) | 25-30 | | | | | | | | |
| Low Impact | 1 | 1 | 2 | 3 | 3 | 4 | 4 | 5 | 5 |
| Moderate Impact | 1 | 2 | 3 | 3 | 4 | 4 | 5 | 5 | - |
| High Impact | 2 | 3 | 3 | 4 | 4 | 4 | 5 | 5 | - |

Product Specification

| Mounting | Build Height | Weight Limit | Impact Level |
|---------------------------------|--------------|----------------------------------|---|
| Endo (IP) | 6.4-6.9cm | 113kg (Sizes 21-24cm) | LMH |
| Exo (ALX) | 12.7-13.2cm | 125kg (Sizes 25-30cm) | |
| Weight of Unit | Sizes (cm) | Warranty | Foot Shell |
| 585g Size 26 inc. foot shell | 21-30cm | 3 years Foot shell - 6 months | C T B J Caucasian, Tan, Brown or Jet Black Narrow or Wide width |

Technical Information

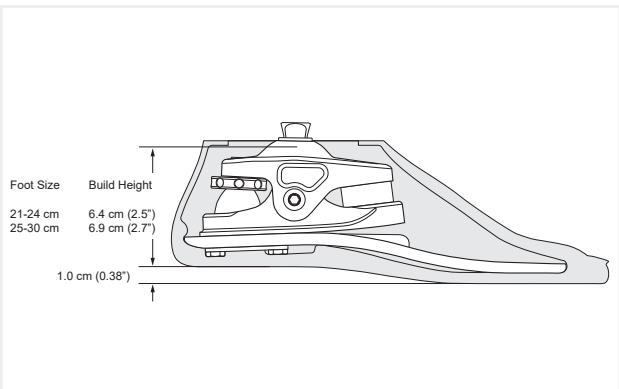

Spare Parts Ordering Information

| Foot Shell | Size (cm) |
|---------------|-----------|
| NV-S ■ ● -▲ N | 21-30 |
| NV-S ■ ● -▲ W | 24-30 |



| CAP | Size (cm) |
|----------------|-----------|
| NV-C ■ ● -2124 | 21-24 |
| NV-C ■ ● -2526 | 25-26 |
| NV-C ■ ● -2730 | 27-30 |



■ = "L" or "R" for Left or Right
● = Colour (C = Caucasian, T = Tan, B = Brown, J = Jet Black)
▲ = Foot size in cm (21-30)
N = Narrow, W = Wide

| Shelltread | Size (cm) |
|---------------|-----------|
| NV-T ■ ● -▲ N | 21-30 |
| NV-T ■ ● -▲ W | 24-30 |



■ = "L" or "R" for Left or Right
● = Colour (N = Neutral, J = Jet Black)
▲ = Foot size in cm (21-30)
N = Narrow, W = Wide




| Part Number | Size (cm) |
|-------------|----------------|
| NV-SS-2124 | 21-24 |
| NV-SS-2526 | 25-26 |
| NV-SS-2730 | 27-30 |
| NV-SSK-2124 | 21-24 (6 pack) |
| NV-SSK-2526 | 25-26 (6 pack) |
| NV-SSK-2730 | 27-30 (6 pack) |



| Foreheel | Description |
|----------------|-------------|
| NV-FHA ■ -▲ -◆ | Foreheel |

■ = "L" or "R" for Left or Right
▲ = Foot size in cm (21-30)
◆ = Firmness (M = Medium, X = Extra Firm)

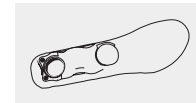


| Bumpers & Bushings | Description | Bumper Colour |
|--------------------|-----------------------------------|---------------|
| NV-FBS- ● | Front Bumper, Soft | S |
| NV-FBM- ● | Front Bumper, Medium | M |
| NV-FBF- ● | Front Bumper, Firm | F |
| NV-FBX-A | Front Bumper, Extra Firm, 27-30cm | XF |
| NV-RBS- ▲ | Rear Bumper, Soft | S |
| NV-RBM-▲ | Rear Bumper, Medium | M |
| NV-RBF- ▲ | Rear Bumper, Firm | F |
| NV-RBX- ▲ | Rear Bumper, Extra Firm | XF |
| NTS-ABS | Ankle Bushing, Soft, 21-30cm | S |
| NTS-ABM | Ankle Bushing, Medium, 21-30cm | M |
| NTS-ABF | Ankle Bushing, Firm, 21-30cm | F |
| NV-ABS-2530 | Ankle Bushing, Soft, 25-30cm | S |
| NV-ABM-2530 | Ankle Bushing, Medium, 25-30cm | M |
| NV-ABF-2530 | Ankle Bushing, Firm, 25-30cm | F |
| NV-BK-2130 | Bumper Kit | |


S Soft M Medium F Firm XF Extra Firm
● = Front Bumper (FB) ▲ = Rear Bumper (RB)

| Pocket Size | Foot Size | | |
|-------------|-----------|---------|---------|
| | 21-24cm | 25-26cm | 27-30cm |
| ● (FB) | C | B | A |
| ▲ (RB) | D | C | B |

| Exo Ankle | Description |
|--------------|--------------------|
| NV-ALXA-2124 | Exo Ankle, 21-24cm |
| NV-ALXA-2530 | Exo Ankle, 25-30cm |

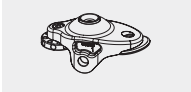

| Endo Ankle | Description |
|-------------|---------------------|
| NV-ENA-2124 | Endo Ankle, 21-24cm |
| NV-ENA-2530 | Endo Ankle, 25-30cm |

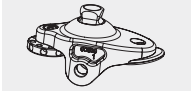


| Growth Plate Kit | Description |
|------------------|------------------------|
| NV-GPK-7 | Growth Plate Kit, 7mm |
| NV-GPK-14 | Growth Plate Kit, 14mm |
| NV-GPK-21 | Growth Plate Kit, 21mm |
| NV-GPK-28 | Growth Plate Kit, 28mm |

| Exo Pyramid Tool Kit | Description |
|----------------------|----------------------------|
| NTS-EXPK-30 | Exo Pyramid Tool Kit, 30mm |



Feet - College Park - Venture

| Exo Block Kit | Description |
|---------------|----------------------------------|
| NV-XBK-2124 | Alignable Exo Block Kit, 21-24cm |
| NV-XBK-2526 | Alignable Exo Block Kit, 25-26cm |
| NV-XBK-2730 | Alignable Exo Block Kit, 27-30cm |

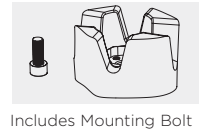

| Hardware | Description |
|------------|--|
| NCP-AP-SH | Axial Pin, Short, 22-24cm |
| NCP-AP-L | Axial Pin, Long, 25-30cm |
| NCP-APS | Axial Pin Screw |
| NTS-S10-25 | Exo Mounting Bolt, M10 x 25mm |
| NV-B10-35 | Exo Mounting Bolt, M10 x 35mm |
| NV-B10-45 | Exo Mounting Bolt, M10 x 45mm |
| NTS-RBK | Pyramid Tool Retaining Bolt Kit (4pc), M6 x 14mm |

