

Freedom Agilix™

The Agilix™ is a multi-axial, shock-absorbing flexible foot system designed to manage loading impacts, reduce socket shear forces and improve comfort while walking on nearly any terrain. The ultra-lightweight design offers K3 ambulators comfortable performance.


Features and Benefits

- Vertical shock absorption and carbon fibre flexibility optimised for low to moderate impact K3 ambulators; reduces socket shear forces and improves overall comfort.
- Split keel design for increased durability
- Multi-axial function increases ground compliance and stability.
- Full-length dynamic heel attached at the toe ensures a seamless roll-over, for superior comfort and a more symmetrical gait.
- The low profile and lightweight design of the Agilix allows users to do more and feel less tired at the end of the day.



Water-resistant - safe for contact with water and submersion*

Cosmetics

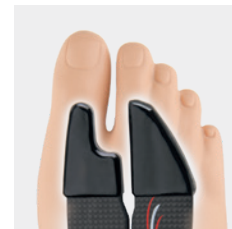

Sandal toe: sizes 22-28cm only

Ordering Information

Trial Period - This foot has a 30 day trial period - only on sizes 23-30 and categories 2-8.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around in the home and the community.

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

Ordering Example - NF15122L

Part Number	Category	Size	Left / Right
NF15	5	22	L

Sandal Toe Ordering Example - NF15122SL

Part Number	Category	Size	Sandal Toe	Left / Right
NF15	5	22	S	L

Sandal toe is only available in sizes 22-28.


Includes foot shell and Spectra® Sock

Category Selection

Weight (kg)	Impact Level	
	LOW	MOD
44-52	1	1
53-59	1	2
60-68	2	3
69-77	3	4
78-88	4	5
89-100	5	6
101-116	6	7
117-130	7	8
131-147	8	9
148-166	9	

*Water note: Rinse with fresh water and dry following exposure to water, particularly salt or chlorinated.

Technical Information

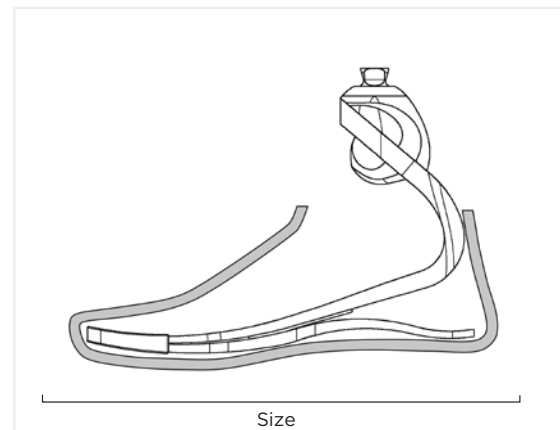

Feet - PROTEOR - Freedom Agilix

Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
166kg	L M	515g	22-31	14.3 (size 25)	3 years Foot shell - 6 months

Build Height and Weight

Size (cm)	Build Height (cm)	Foot inc. Shell
22-25	14.3	553g
26-28	16.1	609g
29-31	17.2	772g

Spare Parts Ordering Information

Part Number	Foot Shell Type
NFTC-3M-0 ▲●4-R	Standard, with CAP
NFTC-2M-1 ▲●4-S	Sandal Toe, without CAP

Spare Parts Ordering Guide

- ▲ = Size in cm (22-31 = Regular, 22-28 = Sandal Toe)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = "L" or "R" for Left or Right

Spectra® Sock	Size (cm)	Sock Size (cm)
NSO-NPS-10030-00	22-25	30
NSO-NPS-10035-00	26-28	35
NSO-NPS-10040-00	29-31	40

X denotes colour. Insert 1 for white, or 2 for black. **Note:** Default colour is black

Bumper Kit	Size (cm)
NKIT-00-11363-00	22-25
NKIT-00-11364-00	26-28
NKIT-00-11365-00	29-31