

Freedom Kinnex 2.0

The Kinnex 2.0 utilises cutting-edge microprocessor technology to provide low to moderate impact ambulators (K3) superior ground compliance for everyday activities. The new hydraulic system has been optimised with improved seal technology and a robust cylinder design. With increased durability, the Kinnex 2.0 delivers reliable energy storage and return throughout the day.

Clinicians can lock the ankle in a neutral position for alignment and calibration, customise the setup of the device, and monitor outcomes using the GaitLab app.


The carbon fibre split heel and keel design provides optimal support and range of motion to reduce internal socket pressure, helping to accommodate obstacles and uneven terrain.



Features and Benefits

- Sensor-controlled ankle movement reads and processes data at a rate of 200 times per second, and Rapid Response™ valve technology responds at a rate of 5 milliseconds
- New hydraulic oil reservoir and seal technology for increased oil capacity and longevity
- IP67 rated - up to 30 minutes in 1 metre of fresh water
- Dynamic lock provides a firm foundation, and can be engaged/disengaged using the app or via the button panel on the front
- Apparel Assist Mode locks the foot at an angle, making it easier for the user to change clothes or footwear
- When charge is below 20%, users will be alerted by a soft vibration and low profile light on the button panel
- Up to 2" of heel height adjustment allows users to program up to 100 shoe styles
- 20° plantar flexion and 10° dorsiflexion range of motion means the foot is ideally positioned with the ground
- User App allows the user to make limited fine resistance adjustments

Cosmetics

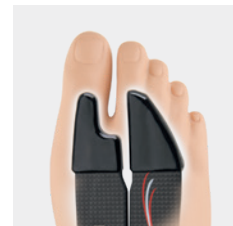

Sandal toe: sizes 24-28 only



Includes foot shell, available in three skin tones



Water resistant - safe for contact with water and submersion*

***Water note:** The Kinnex should not be immersed in salt water or chlorinated water - following exposure to salt water or chlorinated water, immediately rinse with fresh water and dry.

Ordering Information

Trial Period - Please speak to Customer Services or a Steeper Product Manager to organise a trial. **Note:** Once the definitive foot has been purchased there is no trial period.

Impact Level Descriptions

Low: Daily activities include mostly level ground walking, moving around in the home and the community.

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

Ordering Example - NF14526SL

| Part Number | Category | Size | Toe Configuration | Side |
|-------------|----------|------|-------------------|------|
| NF14 | 5 | 26 | S | L |

Includes foot shell (light skin tone as standard, medium and dark also available if specified) and Spectra® Sock

Category Selection

| Weight (kg) | Impact Level | |
|-------------|--------------|-----|
| | LOW | MOD |
| 44-52 | 1 | 1 |
| 53-59 | 1 | 2 |
| 60-68 | 2 | 3 |
| 69-77 | 3 | 4 |
| 78-88 | 4 | 5 |
| 89-100 | 5 | 6 |
| 101-116 | 6 | 7 |
| 117-125 | 7 | - |

Toe Configuration

| Code | Description |
|-------------|--------------|
| Leave Blank | Regular |
| S | Sandal Toe** |

**Sandal toe only available in sizes 24-28cm

| Side | Description |
|------|-------------|
| L | Left |
| R | Right |


Additional Warranty Options

To order a Kinnex 2.0 foot with a five year warranty add the suffix '-5Y' to the end of the part number.

5 year warranty ordering example:

Kinnex 2.0 Foot & 5 Year Warranty, Cat. 5 Size 26 Regular Toe Left **NF14526L-5Y**

Technical Information

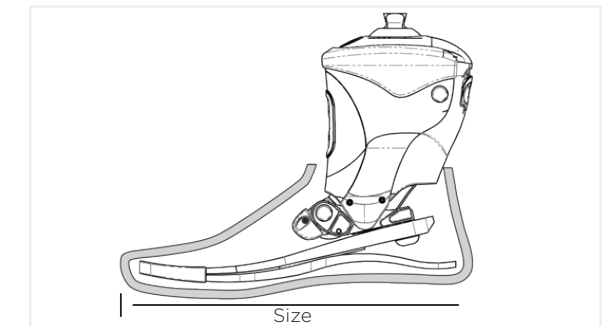

Build Height and Weight

| Size (cm) | Build Height (cm) | Foot inc. Shell |
|-----------|-------------------|-----------------|
| 24-25 | 18.4 | 1432g |
| 26-28 | 18.6 | 1494g |
| 29-30 | 18.6 | 1553g |

Product Specification

| Weight Limit | Impact Level | Weight of Unit | Sizes (cm) | Build Height (cm) | Battery Life | Warranty |
|--------------|--------------|-----------------|------------|-------------------|--------------|-------------------------------|
| 125kg | L M H | 1432g (size 25) | 24-30 | 18.4 (size 25) | 24 hours | 3 years Foot shell - 6 months |

Spare Parts Ordering Information



| Part Number | Foot Shell Type |
|---------------------|------------------------|
| NFTC-3K-1 ▲ ● 4-R ■ | Regular without CAP |
| NFTC-2K-1 ▲ ● 4-S ■ | Sandal Toe without CAP |

Foot Shell Ordering Guide

- ▲ = Size in cm (24-30 = Regular, 24-28 = Sandal Toe)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = Side (R = Right, L = Left)

| Spectra® Sock | Size (cm) | Sock Size | Note: Only available in white colour |
|------------------|-----------|-----------|--------------------------------------|
| NSO-F14-00024-00 | 24 | 24cm | |
| NSO-F14-00025-00 | 25 | 25cm | |
| NSO-F14-00026-00 | 26 | 26cm | |
| NSO-F14-00027-00 | 27 | 27cm | |
| NSO-F14-00028-00 | 28 | 28cm | |
| NSO-F14-00029-00 | 29 | 29cm | |
| NSO-F14-00030-00 | 30 | 30cm | |

| Bumper Kit | Size |
|------------------|----------|
| NKIT-00-11403-00 | 24-25 cm |
| NKIT-00-11404-00 | 26-28 cm |
| NKIT-00-11405-00 | 29-30 cm |