

## LP Symes™

The LP Symes accommodates the build height requirements of patients with Symes amputations. The carbon fibre design delivers energy storage and return, smooth and natural roll-over, and durability in a remarkably lightweight design.

The LP Symes foot incorporates a direct lamination cup, which keeps the distal mass of the prosthesis as light as possible, reducing amputee energy expenditure during the gait cycle. An enhanced material tie-in groove located on the direct lamination cup has also been enlarged to improve the ease of socket connection and lamination.

### Features and Benefits

- Carbon fibre for flexibility and comfort
- Designed specifically for the Symes amputee
- Lightweight
- Direct lamination cup
- Stiffness customised to patient weight



Water resistant - safe for contact with water and submersion\*

### Cosmetics



Foot shell



### Ordering Information

Trial Period - There is no trial period on this foot.

#### Impact Level Descriptions

**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

**High:** Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

**Note:** Although the 'high' impact level states 'running' as an appropriate activity, the LP Symes foot is not ideal for running.

#### Ordering Example - NLP2526L


Part Number	Category	Size	Left / Right
NLP2	5	26	L

Includes foot shell and Spectra® Sock

#### Category Selection

Weight (kg)	Impact Level		
	LOW	MOD	HIGH
44-52	1	1	2
53-59	1	2	3
60-68	2	3	4
69-77	3	4	5
78-88	4	5	6
89-100	5	6	7
101-116	6	7	8
117-130	7	8	9
131-147	8	9	
148-166	9		

### Technical Information

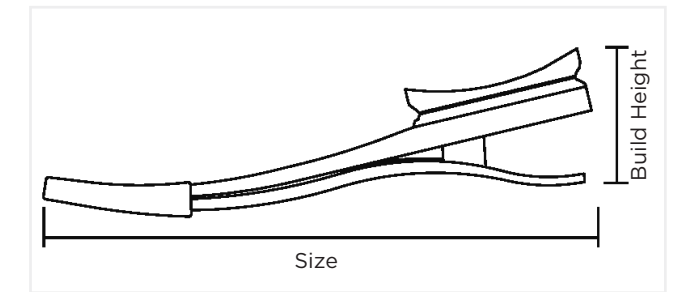

#### Build Height and Weight

Size (cm)	Build Height (cm)	Foot inc. Shell
22-25	4.3	420g
26-28	4.8	545g
29-31	5.3	670g

### Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
147kg		420g (Size 25)	22-31	4.3 (Size 25)	3 years Foot shell - 6 months

\*Water note: Rinse with fresh water and dry following exposure to water, particularly salt or chlorinated.

### Spare Parts Ordering Information

Part Number	Foot Shell Type
NFTC-3M-1▲●4-R■	Standard, without CAP
<b>Lamination plate</b>	
NKIT-00-16200-00	

#### Spare Parts Ordering Guide

- ▲ = Size in cm (22-31)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = "L" or "R" for Left or Right

Bumper Kit	Size (cm)
NKIT-00-11100-00	22-25
NKIT-00-11300-00	26-28
NKIT-00-11500-00	29-31

Spectra® Sock	Size (cm)	Sock Size (cm)
NSO-NPS-X0025-00	22-25	25
NSO-NPS-X0030-00	26-28	30
NSO-NPS-X0035-00	29-31	35

X denotes colour. Insert 1 for white, or 2 for black. **Note:** Default colour is black