

## Maverick Comfort AT

For K3 users who are looking for the ultimate fibreglass walking foot and/or have clearance challenges, the Maverick Comfort AT offers a low profile combined with a split keel, to ensure that users can stride with ease over all terrains.

The Comfort AT offers high energy storage/return and ± 16° inversion/eversion for ample ground compliance.


### Features and Benefits

- Low build height for those with clearance challenges
- Split keel for increased ground compliance and versatility
- Increased durability without compromising the desired energy storage/return
- Smooth roll-over during stance phase of gait
- Lightweight fibreglass
- Innovative heel/keel design



Water resistant - safe for contact with water and submersion\*

### Cosmetics

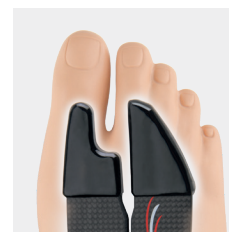

Sandal toe: sizes 23-28 only



### Ordering Information

Trial Period - This foot has a 30 day trial period - only on sizes 23-30 and categories 2-8.

#### Impact Level Descriptions

**Low:** Daily activities include mostly level ground walking, moving around in the home and the community.

**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

**High:** Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

#### Ordering Example - NF22526L

Part Number	Category	Size	Left / Right
NF22	5	26	L

#### Sandal Toe Ordering Example - NF22526SL

Part Number	Category	Size	Sandal Toe	Left / Right
NF22	5	26	S	L


Includes foot shell and Spectra® Sock

#### Category Selection

Weight (kg)	Impact Level		
	LOW	MOD	HIGH
44-52	1	1	2
53-59	1	2	3
60-68	2	3	4
69-77	3	4	5
78-88	4	5	6
89-100	5	6	7
101-116	6	7	8
117-130	7	8	9
131-147	8	9	
148-166	9		

Stiffening Bumpers	Size
NKIT-00-111433-00	23-25 cm
NKIT-00-111434-00	26-28 cm
NKIT-00-111435-00	29-30 cm

### Technical Information

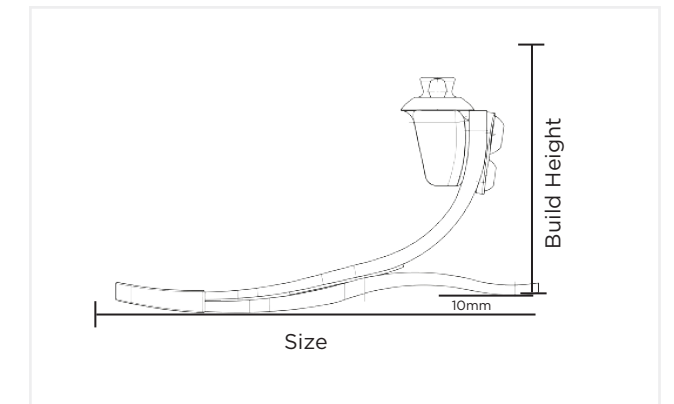

#### Build Height and Weight

Size (cm)	Build Height (cm)	Foot inc. Shell
23-25	11.6	626g
26-28	11.7	763g
29-30	11.9	841g

### Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
166kg	L M H	626g (size 25)	23-30	11.6 (size 25)	3 years Foot shell - 6 months

\*Water note: Rinse with fresh water and dry following exposure to water, particularly salt or chlorinated.

### Spare Parts Ordering Information

Part Number	Foot Shell Type
NFTC-3A-1 ▲●4-R ■	Without CAP
NFTC-2A-1 ▲●4-S ■	With sandal toe, without CAP

#### Foot Shell Ordering Guide

- ▲ = Size in cm (22-31 = Regular, 23-28 = Sandal Toe)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = "L" or "R" for Left or Right

Spectra® Sock	Size (cm)	Sock Size (cm)
NSO-NPS-X0025-00	23-24	25
NSO-NPS-X0030-00	25-27	30
NSO-NPS-X0035-00	28-30	35

X denotes colour. Insert 1 for white, or 2 for black. **Note:** Default colour is black

Bumper Kit	Size (cm)
NKIT-00-11433-00	23-25
NKIT-00-11434-00	26-28
NKIT-00-11435-00	29-30