

## Maverick Xtreme™ AT

The Maverick Xtreme AT fibreglass foot comes with a split keel design that provides excellent inversion and eversion for enhanced ground compliance and patient stability on uneven terrains.

This foot was designed for active K3 and K4 level ambulators who require a product with increased durability without compromising the desired energy storage/return and smooth roll-over during stance phase of gait. The aerospace grade fibreglass combined with the innovative heel/keel design allows for a natural, flexible feel during activities of daily living and various strenuous exercises.


### Features and Benefits

- Multi-axial performance: Split-keel design delivers ± 16° inversion/eversion for all-terrain control and enhanced stability
- The lightest fibreglass foot available, providing all the durability without the compromise of added weight
- Increased durability without compromising the desired energy storage/return
- Smooth roll-over during stance phase of gait
- Lightweight fibreglass
- Innovative heel/keel design
- Regular and sandal toe option
- Designed for high impact, extreme use



Water resistant - safe for contact with water and submersion\*

### Cosmetics

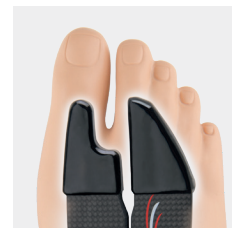

Sandal toe: sizes 23-28 only



### Ordering Information

Trial Period - This foot has a 30 day trial period - only on sizes 23-30 and categories 2-8.

#### Impact Level Descriptions

**Low:** Daily activities include mostly level ground walking, moving around in the home and the community.

**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

**High:** Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

#### Ordering Example - NF21526L

Part Number	Category	Size	Left / Right
NF21	5	26	L

#### Sandal Toe Ordering Example - NF21526SL

Part Number	Category	Size	Sandal Toe	Left / Right
NF21	5	26	S	L


Includes foot shell and Spectra® Sock

#### Category Selection

Weight (kg)	Impact Level		
	LOW	MOD	HIGH
44-52	1	1	2
53-59	1	2	3
60-68	2	3	4
69-77	3	4	5
78-88	4	5	6
89-100	5	6	7
101-116	6	7	8
117-130	7	8	9
131-147	8	9	
148-166	9		

Stiffening Bumpers	Size
NKIT-00-111423-00	23-25 cm
NKIT-00-111424-00	26-28 cm
NKIT-00-111425-00	29-30 cm

### Technical Information

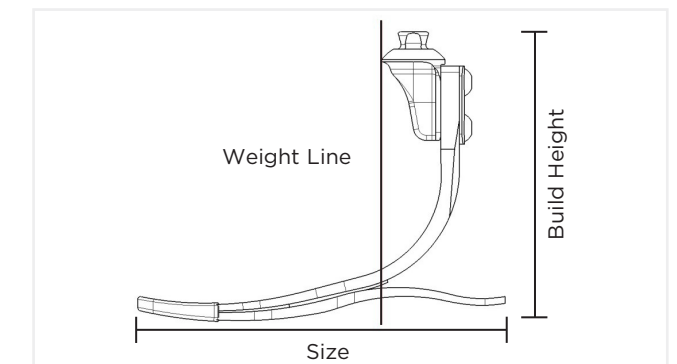

#### Build Height and Weight

Size (cm)	Build Height (cm)	Foot inc. Shell
23-25	16.2	646g
26-28	17	798g
29-30	18.3	891g

### Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
166kg	L M H	646g (size 26)	23-30	16.2 (size 25)	3 years Foot shell - 6 months

\*Water note: Rinse with fresh water and dry following exposure to water, particularly salt or chlorinated.

### Spare Parts Ordering Information

Part Number	Foot Shell Type
NFTC-3A-1▲●4-R■	Without CAP
NFTC-2A-1▲●4-S■	With sandal toe, without CAP

#### Foot Shell Ordering Guide

- ▲ = Size in cm (22-31 = Regular, 23-28 = Sandal Toe)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = "L" or "R" for Left or Right

Spectra® Sock	Size	Sock Size (cm)
NSO-NPS-X0030-00	23-25 cm	30
NSO-NPS-X0035-00	26-28 cm	35
NSO-NPS-X0040-00	29-30 cm	40

X denotes colour. Insert 1 for white, or 2 for black. **Note:** Default colour is black

Bumper Kit	Size (cm)
NKIT-00-11423-00	23-25
NKIT-00-11424-00	26-28
NKIT-00-11425-00	29-30