

## Maverick Xtreme™

The Maverick Xtreme foot was designed for active K3 and K4 level ambulators who require a product with increased durability without compromising the desired energy storage/return and smooth roll-over during stance phase of gait. The increased durability is a result of the fibreglass material used to manufacture the foot and a result of the innovative heel/keel design.

The material combined with the design allows for a natural, flexible feel during activities of daily living and various strenuous exercises. This new, lighter weight fibreglass option is water resistant and is available in a sandal toe option for greater lifestyle versatility.


### Features and Benefits

- The lightest fibreglass foot available, providing all the durability without the compromise of added weight
- Increased durability without compromising the desired energy storage/return
- Smooth roll-over during stance phase of gait
- Lightweight fibreglass
- Innovative heel/keel design
- Regular and sandal toe option
- Designed for high impact, extreme use



Water resistant - safe for contact with water and submersion\*

### Cosmetics

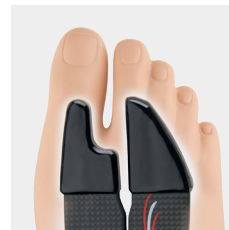

Sandal toe: sizes 23-28 only



### Ordering Information

Trial Period - This foot has a 30 day trial period - only on sizes 23-30 and categories 2-8.

#### Impact Level Descriptions

**Low:** Daily activities include mostly level ground walking, moving around in the home and the community.

**Moderate:** Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

**High:** Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

#### Ordering Example - NF11526L

Part Number	Category	Size	Left / Right
NF11	5	26	L

#### Sandal Toe Ordering Example - NF11526SL

Part Number	Category	Size	Sandal Toe	Left / Right
NF11	5	26	S	L


Includes foot shell and Spectra® Sock

#### Category Selection

Weight (kg)	Impact Level		
	LOW	MOD	HIGH
44-52	1	1	2
53-59	1	2	3
60-68	2	3	4
69-77	3	4	5
78-88	4	5	6
89-100	5	6	7
101-116	6	7	8
117-130	7	8	9
131-147	8	9	
148-166	9		

Stiffening Bumpers	Size
NKIT-00-111413-00	23-25 cm
NKIT-00-111414-00	26-28 cm
NKIT-00-111415-00	29-30 cm

### Technical Information

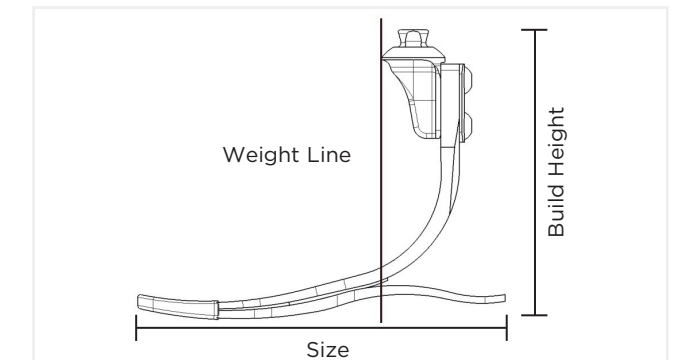

Build Height and Weight		
Size (cm)	Build Height (cm)	Foot inc. Shell
23-25	16.2	650g
26-28	17	821g
29-30	18.3	893g

### Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty
166kg	L M H	650g (size 25)	23-30	16.2 (size 25)	3 years Foot shell - 6 months

\*Water note: Rinse with fresh water and dry following exposure to water, particularly salt or chlorinated.

### Spare Parts Ordering Information

Part Number	Foot Shell Type
NFTC-3A-1▲●4-R■	Without CAP
NFTC-2A-1▲●4-S■	With sandal toe, without CAP

#### Foot Shell Ordering Guide

- ▲ = Size in cm (22-31 = Regular, 23-28 = Sandal Toe)
- = Skin Tone (L = Light, M = Medium, D = Dark)
- = "L" or "R" for Left or Right

Spectra® Sock	Size	Sock Size (cm)
NSO-NPS-X0030-00	23-25 cm	30
NSO-NPS-X0035-00	26-28 cm	35
NSO-NPS-X0040-00	29-30 cm	40

X denotes colour. Insert 1 for white, or 2 for black. **Note:** Default colour is black

Bumper Kit	Size (cm)
NKIT-00-11413-00	23-25
NKIT-00-11414-00	26-28
NKIT-00-11415-00	29-30