


Soleus®

The Soleus foot integrates multiple composite springs through a tough, flexible plantar ligament. This unique design produces a superior quality range of motion while also returning high levels of energy. Providing a very smooth transition between heel strike and toe-off, the Soleus offers stability and unsurpassed comfort.

CES CPI F06: Use this code to read the Soleus clinical evaluation summary on the Steeper website.

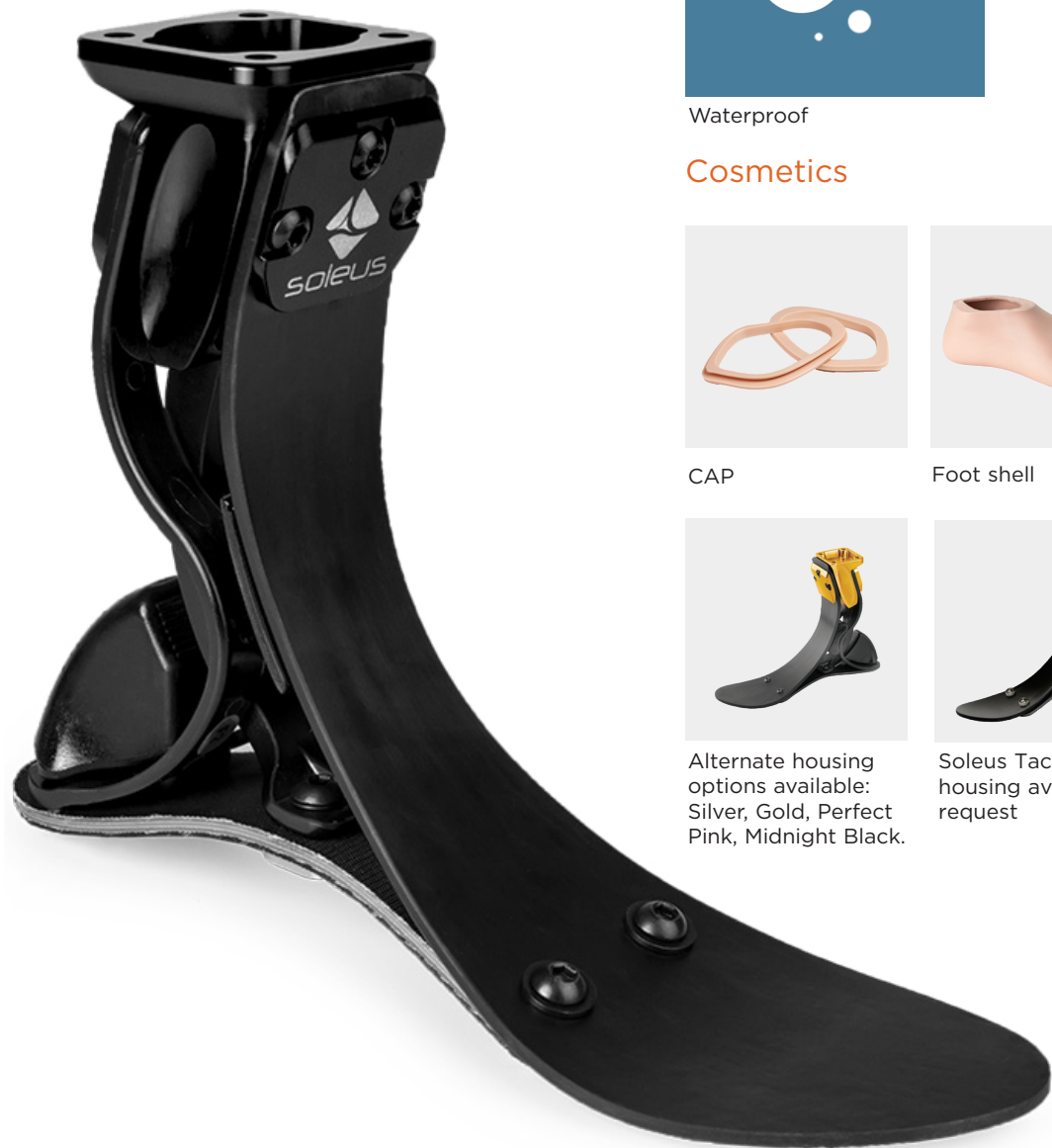

Features and Benefits

- Waterproof (see spare parts for waterproof adapter)
- Multiple spring elements in Integrated Spring Technology (IST) link together to function as one, providing a seamless progression from heel to toe
- Polycentric ankle design provides superior stance phase control
- Vertical shock absorption accommodates high impact activities
- Highly dynamic springs provide multi-axial terrain compliance



Waterproof


Cosmetics




CAP



Foot shell

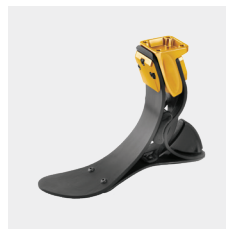

Alternate housing options available: Silver, Gold, Perfect Pink, Midnight Black.

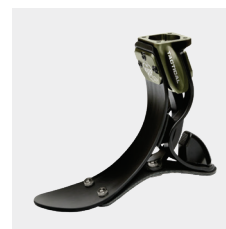

Soleus Tactical housing available on request

Ordering Information

Trial Period - This foot has a 30 day trial period.


Impact Level Descriptions

Moderate: Daily activities include up to unlimited walking, climbing stairs and occasional moderate lifting (this does not include running).

High: Daily activities include fast walking, jogging, running, lifting heavy objects and/or recreational sports.

Soleus foot includes sock, shell and fasteners to mount the 4-hole adapter of your choice. See Order Forms on page 38 for additional ordering information.

Technical Information

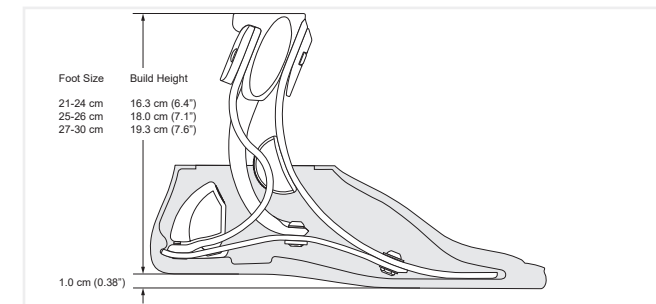

Product Specification

Weight Limit	Impact Level	Weight of Unit	Sizes (cm)	Build Height (cm)	Warranty	Foot Shell
100kg (Sizes 21-24)						
113kg (Sizes 25-26)	MH	667g (Size 26)	21-30	16.3 - 19.3	3 years Foot shell - 6 months	C T B J Caucasian, Tan, Brown or Jet Black
125kg (Sizes 27-30)						Narrow or wide width

Shelltread - Ordering Information

Shelltread	Size
NVT ■●▲	21-30cm
NVT ■●▲ W	25-30cm

Shelltread Ordering Guide
 ■ = "L" or "R" for Left or Right
 ● = Colour (N = Neutral, J = Jet Black)
 ▲ = Foot size in cm (21-30)
 W = Wide

Feet - College Park - Soleus

Spare Parts Ordering Information

Titanium Waterproof Adapter
NTI500B



Foot Shell	Size
NV S ■●▲	21-30
NV S ■●▲ W	24-30

W = Wide



CPI Sock	Size
NSO SS 2126	21-26
NSO SS 2730	27-30

Spare Parts Ordering Guide
 ■ = "L" or "R" for Left or Right
 ● = Colour (C = Caucasian, T = Tan, B = Brown, J = Jet Black)
 ▲ = Foot size in cm (21-30)
 W = Wide