

	Name	Impact Level	K Level	Features and Benefits	Foot Size	Build Height	Weight based on size 25cm inc foot shell	Weight Limit	Warranty	Part Number	Trial Period	Sandal Toe	Foot Shell
K1-K2	Seattle SACH	L M H	K1	<ul style="list-style-type: none"> <li>Easy roll-over</li> <li>Lightweight</li> <li>Natural looking Lightfoot cosmesis</li> <li>Available in both light and dark foot cosmeses</li> <li>Sandal toe</li> </ul>	22-30cm	9cm (25-26cm)	549g	102kg	2 Years Foot Shell - 1 Year	NSCH110- (9.5mm Heel, Light, Lightfoot cosmesis) NSCH113- (9.5mm Heel, Dark, Lightfoot cosmesis)	None	Yes	Light/Dark
	Seattle Natural	L M H	K1-K2	<ul style="list-style-type: none"> <li>Lightweight</li> <li>One standard Delrin keel per foot size</li> <li>Propels the user through cadence</li> <li>Natural looking Lightfoot cosmesis</li> <li>Available in both light and dark foot cosmesis</li> <li>Sandal toe only</li> </ul>	22-29cm	8.5cm	448g	102kg	2 Years Foot Shell - 1 Year	NSNF160- (9.5mm Heel, Light, Lightfoot cosmesis) NSNF163- (9.5mm Heel, Dark, Lightfoot cosmesis)	None	Yes	Light/Dark
	Seattle Kinetic Light	L M H	K2	<ul style="list-style-type: none"> <li>Multi-axial ankle within the foot shell</li> <li>Five exchangeable bumper options through the foot shell using the kinetic tooling kit</li> <li>Excellent foot inversion/eversion and plantar flexion</li> <li>Full-length flexible keel</li> <li>Integral sandal toe cosmesis in two colour options</li> <li>Integral male pyramid</li> </ul>	22-30cm	7.3cm	481g	160kg	2 Years Foot Shell - 1 Year	NSKF100- (9.5mm Heel, Light, Lightfoot cosmesis) NSKF103- (9.5mm Heel, Dark, Lightfoot cosmesis)	None	Yes	Light/Dark
K2(3)	Zumo	L M H	K2(3)	<ul style="list-style-type: none"> <li>Durable composite design for low to moderate users</li> <li>Pre-assembled with a male pyramid for quick and easy installation</li> <li>Split composite keel for increased stability</li> <li>Dynamic toe and cushioned heel response</li> <li>Removable, lifelike, sandal toe foot shell available in light or dark</li> <li>Safe for occasional contact with water, but not intentional, frequent contact or submersion</li> </ul>	25-30cm	9cm	750g	205kg	2 Years Foot Shell - 6 months	NSZM450- (9.5mm Heel, Light) NSZM453- (9.5mm Heel, Dark)	None	Yes	Light/Dark
	Seattle Lightfoot with Pyramid	L M H	K2(3)	<ul style="list-style-type: none"> <li>Low profile selectable keel</li> <li>Ideal for the longer residuum</li> <li>Appropriate for low, medium and high activity levels</li> <li>Natural appearance</li> <li>Available in both light and dark foot shells</li> <li>Sandal toe</li> </ul>	22-30cm	9cm	480g	136kg	2 Years Foot Shell - 1 Year	NSLF195- (9.5mm Heel, Light) NSLF198- (9.5mm Heel, Dark)	None	Yes	Light/Dark
	Seattle Energy	L M H	K2(3)	<ul style="list-style-type: none"> <li>Dynamic toe and cushioned heel response</li> <li>Full length, split composite keel for increased stability</li> <li>Removable, lifelike, sandal toe foot shell</li> <li>Pre-assembled with male pyramid for easy installation</li> <li>Pyramid and spacer can be removed for Symes patient applications</li> <li>Available in both light and dark foot shells</li> <li>Safe for occasional contact with water, but not intentional, frequent contact or submersion</li> </ul>	22-30cm	10cm	330g (Unit only)	136kg	2 Years Foot Shell - 1 Year	NSEF190- (9.5mm Heel, Light) NSEF193- (9.5mm Heel, Dark)	None	Yes	Light/Dark
	Seattle Kinetic Edge	L M H	K2(3)	<ul style="list-style-type: none"> <li>All the benefits of the Kinetic Light with the addition of a removable, sandal toe foot shell and a more lively keel</li> <li>Low profile composite keel with multi-axial ankle</li> <li>Ideal for patients on the edge of advancing to a K3 product</li> <li>Male pyramid attachment</li> </ul>	22-30cm	7cm (27cm)	590g (27cm)	160kg	2 Years Foot Shell - 6 months	NSKF300- (Light) NSKF303- (Dark)	None	Yes	Light/Dark
	Solas	L M H	K3-K4	<ul style="list-style-type: none"> <li>Smooth rollover and controlled plantar flexion</li> <li>Dynamic full length composite keel, split for compliance and stability on uneven surfaces</li> <li>Dynamic response</li> <li>Optional toe wedges are provided to increase dorsiflexion resistance and mid-stance stability</li> <li>Flexible shank provides axial rotation and multiaxial articulation</li> </ul>	22-30cm	18.3cm	670g (26cm)	160kg	3 Years Foot Shell - 6 months	NLSL490- (Light) NLSL493- (Dark)	60 Days	Yes	Light/Dark
Paediatric	Seattle Child's Play Energy Foot	L M H	K3-4	<ul style="list-style-type: none"> <li>Low profile selectable keel design</li> <li>Smooth roll-over and dynamic response</li> <li>Natural appearance</li> <li>Lightweight</li> <li>Choice of foot shell colours</li> <li>Sandal toe (sizes 16-21)</li> </ul>	13-21cm	4.5cm (Size 17)	150g (17cm)	65kg	2 Years Foot Shell - 1 Year	NSEF132- (Light) NSEF144- (Dark)	None	Yes	Light/Dark
Lindhe Xtend	Xtend Foot	L M H	K3	<ul style="list-style-type: none"> <li>Adapts to uneven surfaces through its lightweight, flexible structure</li> <li>+/- 17° of inversion and eversion</li> <li>Enables a more natural gait on slopes and inclines</li> <li>Strong and effective energy return in every step</li> <li>Waterproof* in fresh, chlorinated and salt water</li> </ul>	23-29cm	See Below	695g (26cm)	150kg	3 Years Foot Shell - 6 months	See Below	30 Days	No	Caucasian/Brown


Ordering Information (Xtend Foot)

How to Order:  
For left or right, add L or R after the article code. Example: NH100-426L  
For Caucasian or Brown foot-shell, add FC or FB as a suffix. Example: NH100-425L-FB

Weight Limit	Foot Size (cm)							
	23	24	25	26	27	28	29	
60kg	NH100-323	NH100-324	-	-	-	-	-	
80kg	NH100-423	NH100-424	NH100-425	NH100-426	NH100-427	NH100-428	NH100-429	
100kg	NH100-523	NH100-524	NH100-525	NH100-526	NH100-527	NH100-528	NH100-529	
125kg	-	-	NH100-625	NH100-626	NH100-627	NH100-628	NH100-629	
150kg	-	-	-	NH100-726	NH100-727	NH100-728	NH100-729	

Build Heights

Weight Limit	Foot Size (cm)																					
	23			24			25			26			27			28			29			
	A	B	C	A	B	C	A	B	C	A	B	C	A	B	C	A	B	C	A	B	C	
60kg	81	95	111	82	95	111	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-	-
80kg	81	95	111	82	96	112	80	104	120	90	106	122	88	106	122	89	109	125	92	109	125	
100kg	81	96	112	82	97	113	80	105	121	90	106	122	88	107	123	89	110	126	92	111	127	
125kg	-	-	-	-	-	-	80	106	122	90	107	123	88	107	123	89	112	128	92	112	128	
150kg	-	-	-	-	-	-	-	-	-	90	108	124	88	108	124	89	113	129	92	113	129	

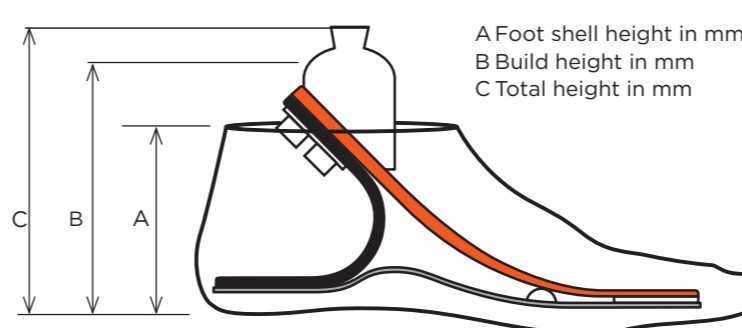

Category Selection Chart (Seattle Lightfoot with Pyramid and Seattle Energy)

To optimise the selection and ensure patient safety, follow the steps below to determine the appropriate keel.

- Locate the column that corresponds with the patient's impact level.
- Within the selected column, locate the patient's weight.
- Place the relevant H number/category number into the part number for ordering.

Note: If the patient has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Seattle Lightfoot with Pyramid		
Cat.	Foot Size (cm)	Low to Medium Weight
H9	24-39	114-136kg
H8	23-29	91-113kg
H7	22-29	69-90kg
H6	22-28	46-68kg
H5	22-27	< 45kg

Seattle Energy		
Cat.	Foot Size (cm)	Low to Medium Weight
4	24-30	101-136kg
3	23-29	81-100kg
2	22-28	61-80kg
1	22-27	< 60kg

Category Selection Instructions (Seattle Child's Play Energy Foot)

To determine which keel is appropriate for your patient, check the below table for keel availability by foot size. Please note that size 13-15 feet do not have sandal toe foot shells.

Important Note: Choose the next higher keel level if your patient has a long BK limb, regularly carries heavy loads, or is at the upper limit for weight or impact level.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please browse the catalogue for other foot options or contact our customer services department.

Category Selection Seattle Child's Play Energy Foot											
Weight Limit	Keel	Foot Size									
		13	14	15	16	17	18	19	20	21	
<15kg	1	✓	✓	✓	-	-	-	-	-	-	-
16-25kg	2	✓	✓	✓	✓	✓	✓	✓	-	-	-
26-35kg	3	✓	✓	✓	✓	✓	✓	✓	✓	✓	✓
36-45kg	4	-	-	-	✓	✓	✓	✓	✓	✓	✓
46-55kg	5	-	-	-	-	-	✓	✓	✓	✓	✓
56-65kg	6	-	-	-	-	-	-	✓	✓	✓	✓

Category Selection Chart (Solus)

To optimise the selection and ensure amputee's safety, follow the steps below to determine the appropriate keel.

- Locate the column that corresponds with the amputee's impact level.
- Within the selected column, locate the amputee's weight.
- Select the relevant foot length.
- Place the relevant number into the part number for ordering.

Note: If the patient has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Category Selection Solas				
Cat.	Foot Length (cm)	Low	Medium	High
4	26-30	137-160kg	124-145kg	106-125kg
3	24-30	101-136kg	89-123kg	76-105kg
2	23-29	81-100kg	71-88kg	61-75kg
1	22-28	<80kg	<70kg	<60kg

Bumper Selection Chart (Seattle Kinetic Light and Kinetic Edge)

Bumpers are coded according to Shore Durometer and based on foot size and patient weight. To optimise selection, follow the steps below to determine the appropriate category.

- Locate the column associated with the patient's foot length.
- Locate the row corresponding to the patient's weight.
- The area where the column and row intersect lists the appropriate category.

Note: If the patient has a long BK or carries heavy loads choose the next category higher.

Patient Weight	Seattle Kinetic Light and Seattle Kinetic Edge									
	Foot Size (cm)									
	22	23	24	25	26	27	28	29	30	
137-160 kg	-	-	-	-	5	5	5	5	5	
101-136 kg	-	-	4	4	4	4	4	4	4	
81-100 kg	-	2	3	3	3	3	3	3	-	
<80 kg	1	1	2	2	2	2	2	-	-	

- Order Guide
- = Foot size in cm (22-30)
  - = "L" or "R" for Left or Right
  - ▲ = Category/Bumper Selection

Note: Category range differs for each foot, please refer to the Steeper Lower Limb Catalogue for category tables.

= Waterproof  
Following exposure to salt and chlorinated water, rinse immediately with fresh water and dry.