

This lightweight polycentric knee is a robust unit with all the features you'd expect from a child's knee, but with one very useful additional feature. It has the facility to be set up as a semi auto or manual knee lock. This makes it particularly suitable for a child with bilateral limb loss, or with other functional issues, such as muscle weakness in the sound limb.

The lock can easily be removed should the child's ability not require it, or when an appropriate level of rehabilitation has been achieved. It can be just as easily replaced should there be a deterioration in function, such as following surgery.

Features and Benefits

- Polycentric design
- Lightweight
- Friction control of swing phase
- Extension assist
- Adjustable stance stability
- Optional manual lock
- Knee disarticulation adapter available (replaces the removable pyramid)



Ordering Information

Loaner units are available for this knee, please contact Customer Services for further details.

Part Number	Description	Weight Limit
NSSK610A	Child's Play Knee with Manual Lock	55kg

Additional Parts Ordering Information

Image	Description	Weight Limit
	Lamination Bracket NSSK610A-DIS	55kg

Image	Description	Weight Limit
	Child's Play Energy Foot NSEF132, 144	65kg

Additional Parts Ordering Information

Image	Description	Weight Limit
	Paediatric Rotatable Receiver - M36 Thread NJR-50FS-A	60kg


Image	Description	Weight Limit
	3 Prong Laminating Base NSS-800BASE	136kg

Product Specification

Weight Limit	Weight	Build Height	Warranty	Distance to Knee Centre
55kg	310g	12cm	2 years	3cm

Technical Information

Front View

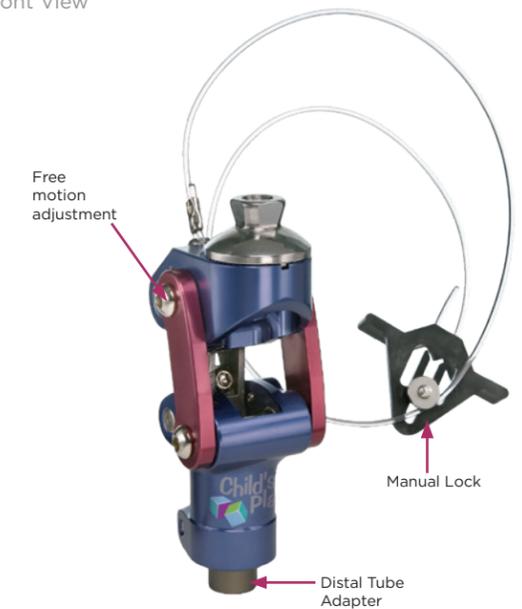

Rear View

