

Kinegen Stream

The Kinegen Stream is a polycentric unit, with hydraulic swing phase control and a pneumatic extension assist, designed for active users weighting up to 150kg. It is geometrically stable in the stance phase and transitions smoothly into the well controlled swing phase, with little effort required and with excellent control of the final 10 degrees of extension, to prevent terminal impact

The Kinegen Stream has a 'cycling' mode enabled by the user by means of a lateral push button with a bicycle symbol.

Features and Benefits

- For active patients weighing up to 150kg
- 'Cycling' mode enabled by the user by means of a lateral push button with a bicycle symbol
- High performance hydraulic unit with precise adjustment allows an optimal, individual customisation of the joint
- Individually adjustable end position damping to reduce terminal impact
- Convenient front access to easily adjustable hydraulic swing phase damper



Convenient front access to damper

'Cycling' mode button




Ordering Information

Trial Period - There is no trial period on this knee, however please contact Customer Services or a Steeper Product Manager to enquire.

Mobility Class & Weight Limit

For active patients weighing up to 150kg.

Part Number	Description
NMA3A2500	Kinegen Stream knee with pyramid adapter (150kg)
NMA3A2500/4	Kinegen Stream knee with M36 (150kg) for use with 4 arm adapter NMA10A25/S
NMA3A2500/3	Kinegen Stream knee with M36 (125kg) for use with 3 arm adapter NMA10A21/S
NMA22A9	Lamination Dummy (suitable for either adapter)

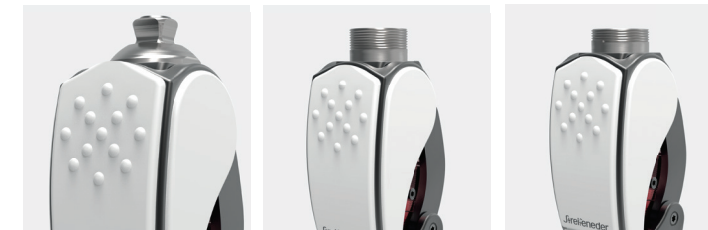

NMA3A2500
Pyramid Adapter

NMA3A2500/4
M36 for 4-Arm
Socket

NMA3A2500/3
M36 for 3-Arm
Socket

Technical Information

	NMA3A2500	NMA3A2500/4	NMA3A2500/3
Proximal Connection	Pyramid adapter	M36 for 4-arm socket	M36 for 3 arm socket
Activity Level	K3	K3	K3
Maximum Patient Weight	150kg	150kg	125kg (the 3 prong socket adapter has a maximum weight limit of 125kg)
Distal Connection	Pyramid adapter	Pyramid adapter	Pyramid adapter
Knee Flexion Angle	136°	136°	136°
Weight	1110g	1150g	1120g
Effective Assembly Height	196mm	209mm	203mm
Total Height Joint Head	36mm	36mm	30mm
Total Height	224mm	222mm	217mm
Material	Aluminium	Aluminium	Aluminium
Warranty	3 years	3 years	3 years