

The Streifeneder silicone liners offer solutions to meet a wide range of user requirements. Each liner uses robust, first-class materials, alongside Streifeneder's new knit technology, to provide a great balance of comfort and function.

The AKControl.Sil is an Above-Knee 2mm silicone liner made with moderately firm silicone.

Features and Benefits


- Ideal silicone liner for distal suspension on residual limbs with reasonable soft tissue coverage
- Suitable for low, moderate and high activity levels
- Lined with antiseptic additives called 'SteriOne', to help keep liners free of bacteria and micro-organisms
- Vertical panels resist elongation
- For above knee patients, alternatively may be useful for a larger transtibial residual limb
- Easy to don and doff due resilient polyamide/elastane cover
- Optional anti-rotation silicone strips prevent twisting between liner and socket (NMA60L10 only)
- Integrated coloured vertical stripe for easier donning orientation
- Distal silicone-rotation control
- Available with distal connection



Antibacterial Technology (SteriOne)



Suitable for Water Use

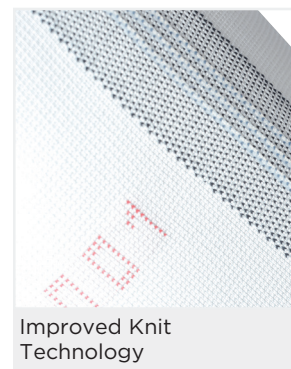

AKControl.Sil



Ordering Information

AKControl.Sil		
Ordering Example NMA60L10/32 (with silicone strips, size 32)		
Size	Silicone Strips	Part Number
25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45, 50, 55	With	NMA60L10/Size
25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45, 50, 55	Without	NMA60L50/Size



Sizing Note: Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

Suitability Chart

Amputation Type	AKControl.Sil
Pear-shaped below-knee (Transtibial)	●
Normal, cylindrical below-knee (Transtibial)	●● for large BK-stumps
Bony, conical below-knee (Transtibial)	●
Cylindrical above-knee (Transfemoral)	●●●
Conical above-knee (Transfemoral)	●

Key: Optimal: ●●● Good: ●● Suitable: ●

Technical Information

Activity Level	K1 to K4
Wall Thickness	2mm
Total Length	Approx 370mm
Max. Patient Weight	150kg
Material	Medium/firm silicone, medical grade, approx. 50 shore 00
Amputation Level	transfemoral/large transtibial
Distal Connector	Pin Adapter or without
Warranty	6 Month
Colour	Silver Grey

Supplied with a drying stand.