

The Streifeneder silicone liners offer solutions to meet a wide range of user requirements. Each liner uses robust, first-class materials, alongside Streifeneder's new knit technology, to provide a great balance of comfort and function.

The ComfortSil prosthetic liner is a 3mm/6mm silicone liner made with soft silicone.



## Features and Benefits

- Ideal silicone liner for those with pressure sensitive residual limbs
- Suitable for low to moderate activity levels
- Designed for those looking for additional comfort and compliance
- Provides additional protection for residual limbs with reduced soft tissue or partially prominent bony anatomy (.soft only)
- Lined with antiseptic additives called 'SteriOne', to help keep liners free of bacteria and micro-organisms
- Integrated 10cm knitted matrix reduces distal elongation under tensile load
- Easy to don and doff due to resilient polyamide/elastane cover
- Integrated coloured vertical stripe for easier donning orientation
- Distal silicone-rotation control
- Available with or without distal connections



Antibacterial Technology (SteriOne)




Suitable for Water Use




Improved Knit Technology

### ComfortSil.basic

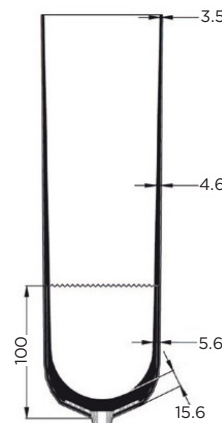

### ComfortSil.soft

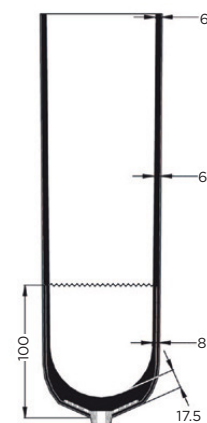

## ComfortSil.basic

### Ordering Information

ComfortSil.basic				
Ordering Example NMA40L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	10cm	with	NMA40L10/Size
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	without	without	NMA40L20/Size



**Sizing Note:** Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

### Technical Information

Activity Level	K1 to K3
Wall Thickness	3mm
Total Length	Approx 370mm
Material	Especially soft silicone, medical grade, approx. 26 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter or without
Warranty	6 Month
Colour	Silver Grey

Supplied with a drying stand.

## ComfortSil.soft

### Ordering Information

ComfortSil.soft				
Ordering Example NMA41L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	6mm	10cm	with	NMA41L10/Size
18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	6mm	without	without	NMA41L20/Size



**Sizing Note:** Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

### Technical Information

Activity Level	K1 to K3
Wall Thickness	6mm
Total Length	Approx 370mm
Material	Especially soft silicone, medical grade, approx. 26 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter or without
Warranty	6 Month
Colour	Silver Grey

Supplied with a drying stand.

## Suitability Chart

Amputation Type	ComfortSil.basic/soft
Pear-shaped below-knee (Transtibial)	● 3mm basic
Normal, cylindrical below-knee (Transtibial)	●●● 3mm basic
Bony, conical below-knee (Transtibial)	●●● 6mm soft

Key: Optimal: ●●● Good: ●● Suitable: ●