

Contex Gel Liner without Matrix

Utilising Streifeneder's unique polymer formula, the ContexGel liner is available for both transfemoral and transtibial application. Supplied in three different wall thicknesses of 3mm, 6mm and 9mm, the liners are also thermo-plastically conformable at 60°C, allowing them to be applied to a wide variety of residual limb shapes and conditions.



Features and Benefits


- Increased comfort as a result of the antibacterial, skin friendly ContexGel
- Protection and good shock absorption over the whole residual limb
- Excellent circumferential elasticity compensates for changes in the residual limb volume
- Thermo-formable at 60°C (140°F) to accommodate the shape of the residual limb, when required
- Textile cover for ease of donning and doffing

Ordering Information

Note: For ContexGel Liner Without Distal Connector change L100 to L200

Part Number	Description	Distal (cm)	Proximal (cm)
NMA20L100/S	ContexGel Liner Small 3mm Distal Connection	16-23	25-36
NMA20L100/M	ContexGel Liner Medium 3mm Distal Connection	18-30	24-41
NMA20L100/L	ContexGel Liner Large 3mm Distal Connection	24-38	38-62
NMA20L100/XL	ContexGel Liner X Large 3mm Distal Connection	28-41	41-65
NMA20L100/XXL	ContexGel Liner XX Large 3mm Distal Connection	32-54	43-73
NMA21L100/S	ContexGel Liner Small 6mm Distal Connection	16-23	25-36
NMA21L100/M	ContexGel Liner Medium 6mm Distal Connection	18-30	24-41
NMA21L100/L	ContexGel Liner Large 6mm Distal Connection	24-38	38-62
NMA21L100/XL	ContexGel Liner X Large 6mm Distal Connection	28-41	41-65
NMA21L100/XXL	ContexGel Liner XX Large 6mm Distal Connection	32-54	43-73
NMA22L100/S	ContexGel Liner Small 9mm Distal Connection	16-23	25-36
NMA22L100/M	ContexGel Liner Medium 9mm Distal Connection	18-30	24-41
NMA22L100/L	ContexGel Liner Large 9mm Distal Connection	24-38	38-62
NMA22L100/XL	ContexGel Liner X Large 9mm Distal Connection	28-41	41-65
NMA22L100/XXL	ContexGel Liner XX Large 9mm Distal Connection	32-54	43-73

Size Information

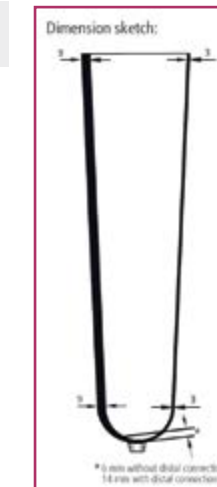

Measure the distal circumference at 4cm from the distal end of the residual limb. Choose the liner size that accommodates that dimension, whilst also taking into consideration the proximal circumference. Having defined the size required, please refer to the table adjacent for part number information.

Technical Information

Material	Polymer Gel - dermatologically tested, inhibits bacteria, hypo allergenic
Matrix	None
Length	Approx 42cm
Amputation Level	Transtibial and Transfemoral*
Distal Connector	Pin Adapter or Cushion
Impact Level	Medium/High
Warranty	6 months

*Note: For transfemoral use it is advised that only the 3mm option should be chosen.

Contex Gel Liner with Matrix

Utilising Streifeneder's unique polymer formula, the ContexGel liner is available for both transfemoral and transtibial application. Supplied in three different wall thicknesses of 3mm, 6mm and 9mm, the liners are also thermo-plastically conformable at 60°C, allowing them to be applied to a wide variety of residual limb shapes and conditions. The addition of a 10cm distal matrix makes this version particularly applicable to those patients requiring improved control of soft residual limb tissue, or where reduction in the distension of a locking liner is required to help prevent socket pistoning.

Features and Benefits

- Increased comfort as a result of the antibacterial, skin friendly ContexGel
- Stitched distal matrix prevents undue longitudinal extension
- Excellent circumferential elasticity compensates for changes in the residual limb volume
- Thermo-formable at 60°C (140°F) to accommodate the shape of the residual limb, when required
- Textile cover for ease of donning and doffing



Ordering Information

Note: For ContexGel Liner Without Distal Connector change L100 to L200

Part Number	Description	Distal (cm)	Proximal (cm)
NMA25L100/S	ContexGel Liner Small 3mm Distal Connection & Matrix	16-23	25-36
NMA25L100/M	ContexGel Liner Medium 3mm Distal Connection & Matrix	18-30	24-41
NMA25L100/L	ContexGel Liner Large 3mm Distal Connection & Matrix	24-38	38-62
NMA25L100/XL	ContexGel Liner X Large 3mm Distal Connection & Matrix	28-41	41-65
NMA25L100/XXL	ContexGel Liner XX Large 3mm Distal Connection & Matrix	32-54	43-73
NMA26L100/S	ContexGel Liner Small 6mm Distal Connection & Matrix	16-23	25-36
NMA26L100/M	ContexGel Liner Medium 6mm Distal Connection & Matrix	18-30	24-41
NMA26L100/L	ContexGel Liner Large 6mm Distal Connection & Matrix	24-38	38-62
NMA26L100/XL	ContexGel Liner X Large 6mm Distal Connection & Matrix	28-41	41-65
NMA26L100/XXL	ContexGel Liner XX Large 6mm Distal Connection & Matrix	32-54	43-73
NMA27L100/S	ContexGel Liner Small 9mm Distal Connection & Matrix	16-23	25-36
NMA27L100/M	ContexGel Liner Medium 9mm Distal Connection & Matrix	18-30	24-41
NMA27L100/L	ContexGel Liner Large 9mm Distal Connection & Matrix	24-38	38-62
NMA27L100/XL	ContexGel Liner X Large 9mm Distal Connection & Matrix	28-41	41-65
NMA27L100/XXL	ContexGel Liner XX Large 9mm Distal Connection & Matrix	32-54	43-73

Size Information



Measure the distal circumference at 4cm from the distal end of the residual limb. Choose the liner size that accommodates that dimension, whilst also taking into consideration the proximal circumference. Having defined the size required, please refer to the table adjacent for part number information.

Technical Information

Material	Polymer Gel - dermatologically tested, inhibits bacteria, hypo allergenic
Matrix	10cm (approx length)
Length	Approx 42cm
Amputation Level	Transtibial and Transfemoral*
Distal Connector	Pin Adapter or Cushion
Impact Level	Medium/High
Warranty	6 months

*Note: For transfemoral use it is advised that only the 3mm option should be chosen.