

Genesee Gel Sleeve

College Park's Genesee Gel Sleeve is suitable as the primary suspension method, or for use with a suction or vacuum suspension. The sleeve is made with hypo-allergenic material and infused with mineral oil to help protect against shear and friction, as well as moisturising skin and scar tissue. With a 3mm dimension, the sleeve provides comfort and knee flexion.

Features and Benefits

- Reinforced interior for puncture-resistance
- Durable fabric-to-gel bond
- Available in 5 sizes
- 2 neutral tones



Reinforced interior



Ordering Information

Ordering Information
Part Number
NGSS- ■ ●
■ - Size (1, 2, 3, 4, 5)
● - Colour (G = Grey, T = Tan)


Technical Information

Proximal Circumference (cm)	Distal Circumference (cm)	Size
22-33	20-30	1
27-38	25-35	2
35-50	30-40	3
40-55	33-48	4
43-60	38-50	5

Product Specification

Material	Thermoplastic Elastomer (TPE) Gel
Length	38cm
Colours	Grey, Tan
Warranty	90 days

Size Selection

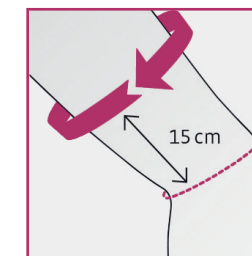

Determining the size
Measure the thigh circumference 15cm above the patella and choose the appropriate size of suspension sleeve.