

1

## Genu Sleeve

The Genu is a 3mm thick soft silicone suspension sleeve, for use on transtibial prostheses. It is gently tapered and, though not pre-flexed, the elasticated fabric outer cover has a more openly knitted section over the patella to allow easy flexion of the knee.

2

This high level of flexibility is achieved thanks to Össur's expertise in round-knit technology. As a consequence the sleeve is virtually wrinkle free even when flexed, and the freedom of movement available to the user is very high.

3

Though this sleeve was designed to be used with a suction system, providing an excellent level of suspension, it can also be used as a stand-alone suspension sleeve and is particularly suitable for patients who are less active, or require a sleeve that allows easier knee flexion.

4

5

The Genu Sleeve is supplied with a sleeve protector as standard, though these can also be purchased separately. Made from a durable fabric, they help prevent damage from the edge of the socket and thereby improve the longevity of the sleeve. Silicone pimples on the bottom edge of the sleeve hold it securely in place on the socket.

6

7

8

9

10

### Features and Benefits

- Highly flexible due to Össur round-knit technology
- Virtually wrinkle-free fit
- High freedom of movement for the user
- Protective sleeve included




### Ordering Information

GENU Sleeve with Protective Sleeve			
Part Number	Size	Circumference (cm)	Length (cm)
NMK00202	Small	30-40	Approx. 35
NMK00203	Medium	38-50	Approx. 35
NMK00204	Large	48-64	Approx. 35

### Product Specification

Thickness	3mm
Warranty	3 months

### Size Selection

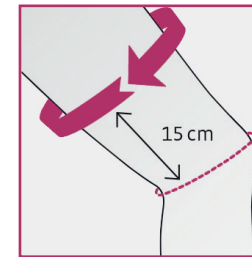

**Determining the size**  
Measure the thigh circumference 15cm above the patella and choose the appropriate size of suspension sleeve.

## Liners - Össur - Genu Sleeve

1

### Additional Parts Ordering Information



Protective Inner Sleeve  
Produced by Össur

Made from a durable fabric, these sleeve protectors help prevent damage from the edge of the socket, thereby improving the longevity of the suspension sleeve. Silicone pimples on the bottom edge of the sleeve hold it securely in place on the socket.

2

3

4

5

6

7

8

9

10

Protective Inner Sleeve	
Part Number	Size
NM000202	Small
NM000203	Medium
NM000204	Large