

1

Product Features


2



silicone gel

- High water vapour, permeability and very good climate exchange characteristics reducing perspiration and odours.
- Very high dimensional stability of the silicone structure for an exceptional fit
- Constant compression gradient, even with long-term use
- Microscopically smooth surface to prevent bacterial colonisation and to allow thorough cleaning.

3

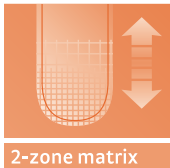

2-zone matrix

- Very high radial stretch for a step-free, gentle and comfortable transition in the matrix area
- Progressive stretch profile for very high soft-tissue grip


4



2-zone-I matrix

- Step-free, gentle transition from the matrix area to the elastic area
- Progressive stretch profile for optimum soft-tissue grip
- Proprioception improved by internal 2-zone matrix

5

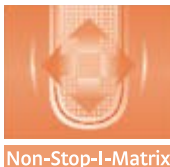

Non-Stop-I-Matrix

- Very good adaptation to various residual limb shapes due to optimum radial elasticity (transverse stretch)
- Gentle to soft-tissue structures by preventing longitudinal stretch

6

7



Easy Glide PLUS

- Proprioception improved by internal 2-zone-matrix
- Unique gliding surface coating for simple donning and doffing
- Simple, residue-free cleaning for the best possible hygiene
- Excellent stretch characteristics through elimination of textile cover

8

9



No Spray

- More independence thanks to Easy Glide PLUS coating
- No alcohol spray needed for turning the liner inside out or putting it on

10




Textile cover

- Durable, special textile cover for easy comfortable donning
- High quality, sheer, elastic Össur knit for an optimum fit
- Glides easily over the skin without a donning aid or spray



Flexible cup

- For simplified, shape-moulding and secure donning of the liner
- Very good damping properties for absorbing and thus reducing the resulting heel-strike impact in the socket
- No impairment of the fit or the wearing comfort because the distal end of the residual limb is ideally moulded and contained

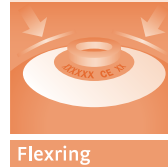

Flexring

- Improved fit which promotes the benefits of the flexible cup
- Neat distal edge, no fraying



Size marking

- Guaranteed follow-up supply
- Size and serial number coded for traceability and documentation



Umbrellan

- Umbrellan is a patented knitted fabric, combining very high elasticity with electromagnetic shielding
- Increased well-being and mobility by alleviation and/or removal of phantom pain, sensations and idiopathic residual limb pain



TF-Matrix

- Adaptation to different residual thigh shapes thanks to very good radial elasticity (transverse stretch)
- Gentle to the soft-tissue structures by inhibiting longitudinal stretch
- More accurate control of the prosthesis with the residual limb due to improved proprioception