



SiLiTOUCH Liners
PROTEOR

 **steeper**
by eqwal.



SILiTOUCH CLASSIC


The SiLiTOUCH Classic Liner integrates an inner layer of silicone and textile coating. This combination promotes elasticity and compliance with a variety of residual limb shapes.

The liner without the distal connection is recommended for residual limbs with bony prominences prone to pressure peaks.


With the distal attachment, the liner provides effective support and suspension for the user, improving comfort as a softer socket interface.

Features and Benefits

- Can be hand or machine washed at 30°C
- Suitable for high activity levels (K2-K4)
- Choose between 3mm or 6mm liner thickness
- 10cm Matrix integrated with distal connection
- Available with or without distal connections
- SiLiTOUCH Classic Liner profiles are equivalent to the Aloe for ease in changing prescription

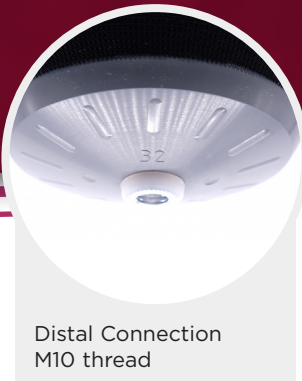

Distal Connection
M10 thread

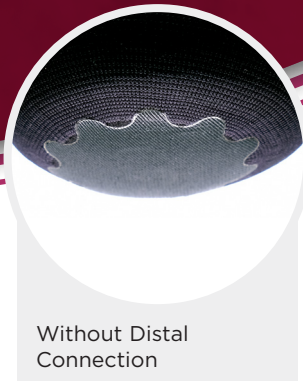


Without Distal
Connection

Ordering Information

| SiLiTOUCH Classic Liner with Distal Connection | | |
|--|-----------|--------------|
| Size | Thickness | Part Number |
| 16 | 3mm | NIS712-3-160 |
| 18 | | NIS712-3-180 |
| 20 | | NIS712-3-200 |
| 22 | | NIS712-3-220 |
| 23.5 | | NIS712-3-235 |
| 25 | | NIS712-3-250 |
| 26.5 | | NIS712-3-265 |
| 28 | | NIS712-3-280 |
| 30 | | NIS712-3-300 |
| 32 | | NIS712-3-320 |
| 34 | 6mm | NIS712-3-340 |
| 36 | | NIS712-3-360 |
| 40 | | NIS712-3-400 |
| 45 | | NIS712-3-450 |
| 20 | | NIS712-6-200 |
| 22 | | NIS712-6-220 |
| 23.5 | | NIS712-6-235 |
| 25 | | NIS712-6-250 |
| 26.5 | | NIS712-6-265 |
| 28 | | NIS712-6-280 |

| SiLiTOUCH Classic Liner without Distal Connection | | |
|---|-----------|--------------|
| Size | Thickness | Part Number |
| 16 | 3mm | NIS722-3-160 |
| 18 | | NIS722-3-180 |
| 20 | | NIS722-3-200 |
| 22 | | NIS722-3-220 |
| 23.5 | | NIS722-3-235 |
| 25 | | NIS722-3-250 |
| 26.5 | | NIS722-3-265 |
| 28 | | NIS722-3-280 |
| 30 | | NIS722-3-300 |
| 32 | | NIS722-3-320 |
| 34 | 6mm | NIS722-3-340 |
| 36 | | NIS722-3-360 |
| 40 | | NIS722-3-400 |
| 45 | | NIS722-3-450 |
| 16 | | NIS722-6-160 |
| 18 | | NIS722-6-180 |
| 20 | | NIS722-6-200 |
| 22 | | NIS722-6-220 |
| 23.5 | | NIS722-6-235 |
| 25 | | NIS722-6-250 |

Size Selection

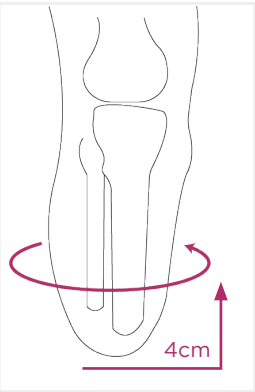

Determining the size
Measure the circumference of the residual limb 4cm above the distal end of the residual limb.

Technical Information

| | SiLiTOUCH Classic Liner | |
|-------------------|-------------------------|--------------------------|
| | With Distal Connector | Without Distal Connector |
| Material | Silicone | |
| Shore Hardness | 35 | |
| Thickness | 3mm or 6mm | |
| Matrix | 10cm | - |
| Amputation Level | Transtibial | |
| Activity Level | K2-K4 | |
| Warranty | 6 Months | |
| Length of Liner | 440mm | 450mm |
| Distal Connection | M10 Thread | - |

SILiTOUCH ALOE


The SiLiTOUCH Aloe liner creates a comfortable fit for users. With the silicone and Aloe Vera combination, this makes this liner most suitable for those with sensitive skin.

The liner without the distal connection is recommended for residual limbs with bony prominences prone to pressure peaks.


With the distal attachment, the liner provides a comfortable, effective support and suspension for the user.

Features and Benefits

- Can be hand or machine washed at 30°C
- Suitable for high activity levels (K2-K4)
- Choose between 3mm or 6mm liner thickness
- 10cm Matrix integrated with distal connection
- Available with or without distal connections
- SiLiTOUCH Aloe Liner profiles are equivalent to the Classic for ease in changing prescription
- Aloe Vera-infused silicone supports skin at the interface

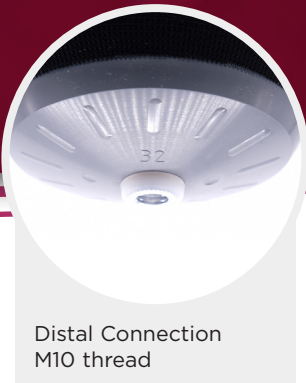

Distal Connection
M10 thread

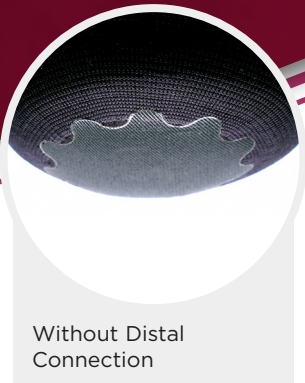

Without Distal
Connection



Added Aloe Vera to prevent skin irritation




Ordering Information

| SiLiTOUCH Aloe Liner with Distal Connection | | |
|---|-----------|--------------|
| Size | Thickness | Part Number |
| 16 | 3mm | NIS711-3-160 |
| 18 | | NIS711-3-180 |
| 20 | | NIS711-3-200 |
| 22 | | NIS711-3-220 |
| 23.5 | | NIS711-3-235 |
| 25 | | NIS711-3-250 |
| 26.5 | | NIS711-3-265 |
| 28 | | NIS711-3-280 |
| 30 | | NIS711-3-300 |
| 32 | | NIS711-3-320 |
| 34 | 6mm | NIS711-3-340 |
| 36 | | NIS711-3-360 |
| 40 | | NIS711-3-400 |
| 45 | | NIS711-3-450 |
| 16 | | NIS711-6-160 |
| 18 | | NIS711-6-180 |
| 20 | | NIS711-6-200 |
| 22 | | NIS711-6-220 |
| 23.5 | | NIS711-6-235 |
| 25 | | NIS711-6-250 |
| 26.5 | | NIS711-6-265 |
| 28 | | NIS711-6-280 |
| 30 | | NIS711-6-300 |
| 32 | | NIS711-6-320 |
| 34 | | NIS711-6-340 |
| 36 | | NIS711-6-360 |
| 40 | | NIS711-6-400 |
| 45 | | NIS711-6-450 |

| SiLiTOUCH Aloe Liner without Distal Connection | | |
|--|-----------|--------------|
| Size | Thickness | Part Number |
| 16 | 3mm | NIS721-3-160 |
| 18 | | NIS721-3-180 |
| 20 | | NIS721-3-200 |
| 22 | | NIS721-3-220 |
| 23.5 | | NIS721-3-235 |
| 25 | | NIS721-3-250 |
| 26.5 | | NIS721-3-265 |
| 28 | | NIS721-3-280 |
| 30 | | NIS721-3-300 |
| 32 | | NIS721-3-320 |
| 34 | 6mm | NIS721-3-340 |
| 36 | | NIS721-3-360 |
| 40 | | NIS721-3-400 |
| 45 | | NIS721-3-450 |
| 16 | | NIS721-6-160 |
| 18 | | NIS721-6-180 |
| 20 | | NIS721-6-200 |
| 22 | | NIS721-6-220 |
| 23.5 | | NIS721-6-235 |
| 25 | | NIS721-6-250 |
| 26.5 | | NIS721-6-265 |
| 28 | | NIS721-6-280 |
| 30 | | NIS721-6-300 |
| 32 | | NIS721-6-320 |
| 34 | | NIS721-6-340 |
| 36 | | NIS721-6-360 |

Size Selection

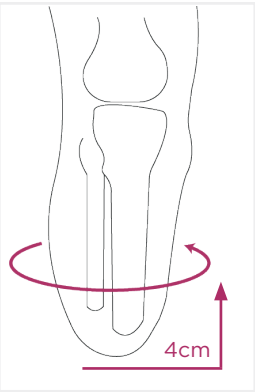

Determining the size
Measure the circumference of the residual limb 4cm above the distal end of the residual limb.

Technical Information

| | SiLiTOUCH Aloe Liner | |
|-------------------|-----------------------|--------------------------|
| | With Distal Connector | Without Distal Connector |
| Material | Silicone | |
| Shore Hardness | 20 | |
| Thickness | 3mm or 6mm | |
| Matrix | 10cm | - |
| Amputation Level | Transtibial | |
| Activity Level | K2-K4 | |
| Warranty | 6 Months | |
| Length of Liner | 440mm | 450mm |
| Distal Connection | M10 Thread | - |

Prosthetist Notes

Steeper Group
Unit 3 Stourton Link, Intermezzo Drive
Leeds, LS10 1DF

Email: customerservices@steepergroup.com
Call: +44 (0) 870 240 4133

www.steepergroup.com

STPPR340 V1, April 2025
©2025 Steeper Group

