

The Streifeneder silicone liners offer solutions to meet a wide range of user requirements. Each liner uses robust, first-class materials, alongside Streifeneder's new knit technology, to provide a great balance of comfort and function.

The SkinCareSil prosthetic liner is a 3mm/6mm silicone liner made with soft silicone.



Features and Benefits

- Ideal silicone liner for those with pressure sensitive residual limbs
- Suitable for low activity levels
- Designed with Aloe Vera to promote skin integrity
- Provides additional protection for residual limbs with reduced soft tissue or partially prominent bony anatomy (.soft only)
- Lined with antiseptic additives called 'SteriOne', to help keep liners free of bacteria and micro-organisms
- Integrated 10cm knitted matrix reduces distal elongation under tensile load
- Easy to don and doff due resilient polyamide/elastane cover
- Integrated coloured vertical stripe for easier donning orientation
- Distal silicone rotation control
- Available with or without distal connection




Made using Aloe Vera Extract



Antibacterial Technology (SteriOne)




Suitable for Water Use

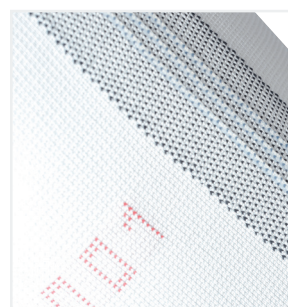


Improved Knit Technology

SkinCareSil.basic

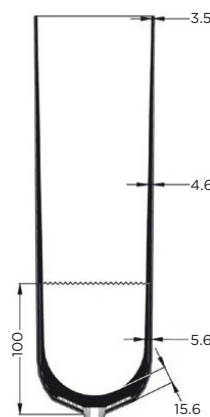

SkinCareSil.soft

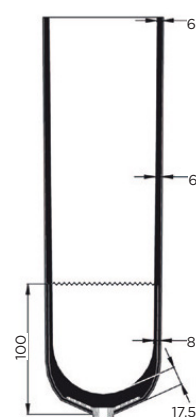

SkinCareSil.basic Ordering Information

SkinCareSil.basic				
Ordering Example NMA46L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	10cm	with	NMA46L10/Size
16, 18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	3mm	without	without	NMA46L20/Size



Sizing Note: Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

Technical Information

Activity Level	K1 to K2
Wall Thickness	3mm
Total Length	Approx 370mm
Material	Especially soft silicone, medical grade, approx. 28 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter
Warranty	6 Month
Colour	Silver Grey

Supplied with a drying stand.

SkinCareSil.soft Ordering Information

SkinCareSil.soft				
Ordering Example NMA48L10/18 (with distal connection, size 18)				
Size	Thickness	Outer Matrix	Distal Connection	Part Number
18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	6mm	10cm	with	NMA48L10/Size
18, 20, 21, 22, 23.5, 25, 26.5, 28, 30, 32, 34, 36, 38, 40, 42, 45	6mm	without	without	NMA48L20/Size



Sizing Note: Measure the distal circumference approx. 4cm above the distal end and select the liner in the same size or one size smaller.

Technical Information

Activity Level	K1 to K2
Wall Thickness	6mm
Total Length	Approx 370mm
Material	Especially soft silicone, medical grade, approx. 28 shore 00
Amputation Level	Transtibial
Distal Connector	Pin Adapter or without
Warranty	6 Month
Colour	Silver Grey

Supplied with a drying stand.

Suitability Chart

Amputation Type	SkinCareSil.basic/soft
Pear-shaped below-knee (Transtibial)	●●● 3mm basic
Normal, cylindrical below-knee (Transtibial)	●●●● 3mm basic
Bony, conical below-knee (Transtibial)	●●●● 6mm soft

Key: Optimal: ●●●● Good: ●●● Suitable: ●