

Uniprox Liner Overview

SoftskinAir®

Breathable Silicone Liner

Softskin

Silicone Liner

With distal connector

Without distal connector

With distal connector

Without distal connector

Firm silicone on the outside, Soft Silicone on the inside

NS-2L-AIR - NS-2L -



NS30-AIR NS30CL-AIR NS30 NS30CL



Soft Silicone

NS40-AIR NSPORT-AIR NSC40-AIR NS40CL-AIR NS40 NSC40 NS40CL



Medium Silicone

NS50-AIR NS50CL-AIR NS50 NS50CL



Firm Silicone

NSTF-AIR* - NSTF* -



Transfemoral Liners


with consistent matrix for transfemoral users*



The SoftSkin standard silicone liner range by Uniprox, offers a selection of shore hardnesses with an optional pin connection, or with a smooth distal contour. The liners include a matrix to reduce pistoning whilst the design of the outer knitted fabric allows the socket to be donned more easily.

Features and Benefits

- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Machine-washable (30°C)
- Available in three shore hardnesses: from the softer NS30 for more sensitive residual limbs, through to the stiffer NS50 liners suitable for higher activity
- The smooth contour (cushion) liners are particularly suitable for bony residual limbs, through-knee amputations and above knee residua with deep scars

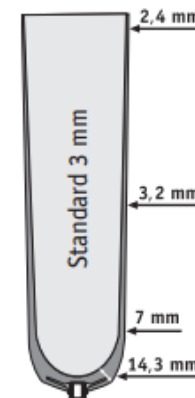

Ordering Information

All liners are 3mm thick, 370mm long.

Note: The circumference is measured 4cm above the distal end of the residual limb

| Size | With Distal Connector | | | Without Distal Connector | | | Smooth Contour |
|------|-----------------------|--------------|--------------|--------------------------|---------------|---------------|---------------------|
| | SoftSkin S30 | SoftSkin S40 | SoftSkin S50 | SoftSkin CL30 | SoftSkin CL40 | SoftSkin CL50 | SoftSkin Cushion 40 |
| 16 | NS30-16 | NS40-16 | NS50-16 | NS30CL-16 | NS40CL-16 | NS50CL-16 | - |
| 18 | NS30-18 | NS40-18 | NS50-18 | NS30CL-18 | NS40CL-18 | NS50CL-18 | - |
| 20 | NS30-20 | NS40-20 | NS50-20 | NS30CL-20 | NS40CL-20 | NS50CL-20 | - |
| 21 | NS30-21 | NS40-21 | NS50-21 | NS30CL-21 | NS40CL-21 | NS50CL-21 | - |
| 22 | NS30-22 | NS40-22 | NS50-22 | NS30CL-22 | NS40CL-22 | NS50CL-22 | - |
| 23.5 | NS30-23 | NS40-23 | NS50-23 | NS30CL-23 | NS40CL-23 | NS50CL-23 | NSC40-23 |
| 25 | NS30-25 | NS40-25 | NS50-25 | NS30CL-25 | NS40CL-25 | NS50CL-25 | NSC40-25 |
| 26.5 | NS30-26 | NS40-26 | NS50-26 | NS30CL-26 | NS40CL-26 | NS50CL-26 | NSC40-26 |
| 28 | NS30-28 | NS40-28 | NS50-28 | NS30CL-28 | NS40CL-28 | NS50CL-28 | NSC40-28 |
| 30 | NS30-30 | NS40-30 | NS50-30 | NS30CL-30 | NS40CL-30 | NS50CL-30 | NSC40-30 |
| 32 | NS30-32 | NS40-32 | NS50-32 | NS30CL-32 | NS40CL-32 | NS50CL-32 | NSC40-32 |
| 34 | NS30-34 | NS40-34 | NS50-34 | NS30CL-34 | NS40CL-34 | NS50CL-34 | NSC40-34 |
| 36 | NS30-36 | NS40-36 | NS50-36 | NS30CL-36 | NS40CL-36 | NS50CL-36 | NSC40-36 |
| 38 | NS30-38 | NS40-38 | NS50-38 | NS30CL-38 | NS40CL-38 | NS50CL-38 | NSC40-38 |
| 40 | NS30-40 | NS40-40 | NS50-40 | NS30CL-40 | NS40CL-40 | NS50CL-40 | - |
| 42 | NS30-42 | NS40-42 | NS50-42 | NS30CL-42 | NS40CL-42 | NS50CL-42 | - |
| 45 | NS30-45 | NS40-45 | NS50-45 | NS30CL-45 | NS40CL-45 | NS50CL-45 | - |

Technical Information

| | SoftSkin NS30/NS40/NS50 | SoftSkin NS30CL/NS40CL/NS50CL | SoftSkin Cushion 40 | | |
|------------------|--------------------------|-------------------------------|---|---------------------|-----------|
| Material | Silicone | Silicone | Silicone | | |
| Shore Hardness | 30/40/50 | 30/40/50 | 40 | | |
| Thickness | 3mm | 3mm | 3mm | | |
| Length | 370mm | 370mm | 370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<) | | |
| Amputation Level | Transtibial | Transtibial | Transtibial/ Transfemoral | | |
| Distal Connector | With Distal Connector | Without Distal Connector | Smooth Contour | | |
| Mobility Grades | NS30 & NS40 1,2 and 3 | NS50 2,3 and 4 | NS30CL & NS40CL 1,2 and 3 | NS50CL 2,3 and 4 | 1,2 and 3 |
| Warranty | 6 months | 6 months | 6 months | | |

4 Silicone Liners

4 Silicone Liners

The SoftSkinAir® is a unique, breathable silicone liner featuring a distinctive structure and special micropores to reduce the build up of sweat and draw moisture away from the residual limb; providing a unique solution to patients who previously have been unable to wear a liner due to high sweat levels. It has proven effectiveness in wearing tests and moisture level analyses.

Available in multiple shore hardnesses, with an optional pin connection, or with a smooth distal contour.

The SoftSkinAir® Sport liner offers an additional distal matrix to further improve the reduction of longitudinal movements.


Features and Benefits

- Wicks sweat via micropores
- Unique functional and progressive outside texture
- Available with or without distal connector
- Machine-washable (30°C)
- Breathable and comfortable solution
- Improved grip contact
- Unique outer textile allows for easier socket donning
- Decreased skin irritation
- Reduced temperature inside the liner
- Available in three shore hardnesses: from the softer NS30 for more sensitive residual limbs, through to the durable NS50 liners suitable for higher activities

Air Sport

Additional Features and Benefits

- Excellent wicking of sweat through micropores
- Shore hardness: 40 Shore A (soft)
- Additional distal matrix for more stability
- Sizes: 20 - 30
- Dual layer with silky inner surface
- Reduction of elongation

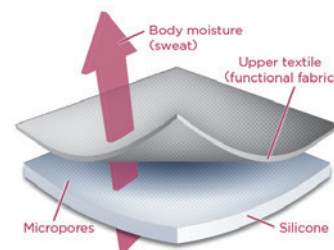

Ordering Information

All liners are 3mm thick, 370mm long.

Note: The circumference is measured 4cm above the distal end of the residual limb

| Size | With Distal Connector | | | Without Distal Connector | | | Smooth Contour |
|------|-----------------------|------------------|------------------|--------------------------|-------------------|-------------------|-------------------------|
| | SoftSkin Air S30 | SoftSkin Air S40 | SoftSkin Air S50 | SoftSkin Air CL30 | SoftSkin Air CL40 | SoftSkin Air CL50 | SoftSkin Air Cushion 40 |
| 16 | NS30-16-AIR | NS40-16-AIR | NS50-16-AIR | NS30CL-16-AIR | NS40CL-16-AIR | NS50CL-16-AIR | |
| 18 | NS30-18-AIR | NS40-18-AIR | NS50-18-AIR | NS30CL-18-AIR | NS40CL-18-AIR | NS50CL-18-AIR | |
| 20 | NS30-20-AIR | NS40-20-AIR | NS50-20-AIR | NS30CL-20-AIR | NS40CL-20-AIR | NS50CL-20-AIR | |
| 21 | NS30-21-AIR | NS40-21-AIR | NS50-21-AIR | NS30CL-21-AIR | NS40CL-21-AIR | NS50CL-21-AIR | |
| 22 | NS30-22-AIR | NS40-22-AIR | NS50-22-AIR | NS30CL-22-AIR | NS40CL-22-AIR | NS50CL-22-AIR | |
| 23.5 | NS30-23-AIR | NS40-23-AIR | NS50-23-AIR | NS30CL-23-AIR | NS40CL-23-AIR | NS50CL-23-AIR | NSC40-23-AIR |
| 25 | NS30-25-AIR | NS40-25-AIR | NS50-25-AIR | NS30CL-25-AIR | NS40CL-25-AIR | NS50CL-25-AIR | NSC40-25-AIR |
| 26.5 | NS30-26-AIR | NS40-26-AIR | NS50-26-AIR | NS30CL-26-AIR | NS40CL-26-AIR | NS50CL-26-AIR | NSC40-26-AIR |
| 28 | NS30-28-AIR | NS40-28-AIR | NS50-28-AIR | NS30CL-28-AIR | NS40CL-28-AIR | NS50CL-28-AIR | NSC40-28-AIR |
| 30 | NS30-30-AIR | NS40-30-AIR | NS50-30-AIR | NS30CL-30-AIR | NS40CL-30-AIR | NS50CL-30-AIR | NSC40-30-AIR |
| 32 | NS30-32-AIR | NS40-32-AIR | NS50-32-AIR | NS30CL-32-AIR | NS40CL-32-AIR | NS50CL-32-AIR | NSC40-32-AIR |
| 34 | NS30-34-AIR | NS40-34-AIR | NS50-34-AIR | NS30CL-34-AIR | NS40CL-34-AIR | NS50CL-34-AIR | NSC40-34-AIR |
| 36 | NS30-36-AIR | NS40-36-AIR | NS50-36-AIR | NS30CL-36-AIR | NS40CL-36-AIR | NS50CL-36-AIR | NSC40-36-AIR |
| 38 | NS30-38-AIR | NS40-38-AIR | NS50-38-AIR | NS30CL-38-AIR | NS40CL-38-AIR | NS50CL-38-AIR | NSC40-38-AIR |
| 40 | NS30-40-AIR | NS40-40-AIR | NS50-40-AIR | NS30CL-40-AIR | NS40CL-40-AIR | NS50CL-40-AIR | |
| 42 | NS30-42-AIR | NS40-42-AIR | NS50-42-AIR | NS30CL-42-AIR | NS40CL-42-AIR | NS50CL-42-AIR | |
| 45 | NS30-45-AIR | NS40-45-AIR | NS50-45-AIR | NS30CL-45-AIR | NS40CL-45-AIR | NS50CL-45-AIR | |

| SoftSkinAir Sport® (With Distal Connector) | | | | | | | |
|--|---------------|-----------------|---------------|-----------------|---------------|---------------|--|
| Size | | | | | | | |
| 20 | 22 | 23.5 | 25 | 26.5 | 28 | 30 | |
| NSPORT-20-AIR | NSPORT-22-AIR | NSPORT-23.5-AIR | NSPORT-25-AIR | NSPORT-26.5-AIR | NSPORT-28-AIR | NSPORT-30-AIR | |

Technical Information

| | SoftSkin Air NS30/NS40/NS50 | SoftSkin Air Sport | SoftSkin Air NS30CL/NS40CL/ NS50CL | SoftSkin Air Cushion 40 |
|------------------|-----------------------------|--------------------------|------------------------------------|---|
| Material | Silicone with micropores | Silicone with micropores | Silicone with micropores | Silicone with micropores |
| Shore Hardness | 30/40/50 | 40 | 30/40/50 | 40 |
| Thickness | 3mm | 3mm | 3mm | 3mm |
| Length | 370mm | 370mm | 370mm | 370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<) |
| Amputation Level | Transtibial | Transtibial | Transtibial | Transtibial/ Transfemoral |
| Distal Connector | With Distal Connector | With Distal Connector | Without Distal Connector | Smooth Contour |
| Mobility Grades | NS30 & NS40 1,2 and 3 | NS50 2,3 and 4 | NS30CL & NS40CL 1,2 and 3 | NS50CL 2,3 and 4 |
| Warranty | 6 months | 6 months | 6 months | 6 months |

SoftSkin 2L / SoftSkinAir® 2L

The SoftSkin 2L and SoftSkinAir® 2L are dual layer liners featuring a soft silicone inner for gentle skin contact, with a firmer silicone outer to provide stability.

The SoftSkinAir® 2L offers the additional benefit of sweat wicking micropores to reduce moisture, resulting in drier skin and lower temperature within the liner for residual limb comfort and reduced skin irritation.



Features and Benefits

- Soft on the inside, stable on the outside (Shore hardness 30/50)
- Soft residual limb padding
- Innovative functional knitted fabric to reduce the longitudinal movement without an internal matrix
- Sizes: 16 – 45
- Unique functional outer fabric
- Machine-washable (30°C)
- Improved grip contact
- Decreased skin irritation

SoftSkinAir® 2L

Additional Features and Benefits

- Breathable and comfortable solution
- Wicks sweat via micropores

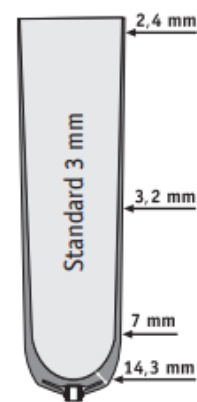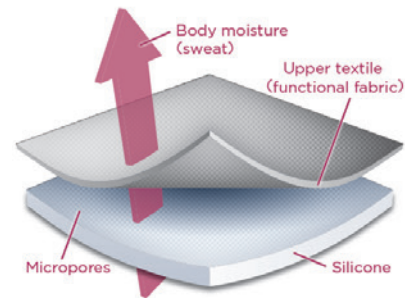

Liners - Uniprox® - SoftSkin 2L / SoftSkinAir® 2L

Ordering Information

All liners are 3mm thick, 370mm long.

Note: The circumference is measured 4cm above the distal end of the residual limb

| Size | With Distal Connector | |
|------|-----------------------|-----------------|
| | SoftSkin 2L | SoftSkinAir® 2L |
| 16 | NS-2L-16 | NS-2L-16-AIR |
| 18 | NS-2L-18 | NS-2L-18-AIR |
| 20 | NS-2L-20 | NS-2L-20-AIR |
| 21 | NS-2L-21 | NS-2L-21-AIR |
| 22 | NS-2L-22 | NS-2L-22-AIR |
| 23.5 | NS-2L-23 | NS-2L-23-AIR |
| 25 | NS-2L-25 | NS-2L-25-AIR |
| 26.5 | NS-2L-26 | NS-2L-26-AIR |
| 28 | NS-2L-28 | NS-2L-28-AIR |
| 30 | NS-2L-30 | NS-2L-30-AIR |
| 32 | NS-2L-32 | NS-2L-32-AIR |
| 34 | NS-2L-34 | NS-2L-34-AIR |
| 36 | NS-2L-36 | NS-2L-36-AIR |
| 38 | NS-2L-38 | NS-2L-38-AIR |
| 40 | NS-2L-40 | NS-2L-40-AIR |
| 42 | NS-2L-42 | NS-2L-42-AIR |
| 45 | NS-2L-45 | NS-2L-45-AIR |

Technical Information

| | SoftSkin 2L | SoftSkinAir® 2L |
|------------------|-----------------------|--------------------------|
| Material | Silicone | Silicone with micropores |
| Shore Hardness | 30 inner/50 outer | 30 inner/50 outer |
| Thickness | 3mm | 3mm |
| Length | 370mm | 370mm |
| Amputation Level | Transtibial | Transtibial |
| Distal Connector | With Distal Connector | With Distal Connector |
| Mobility Grades | 1,2 and 3 | 1,2 and 3 |
| Warranty | 6 months | 6 months |

The SoftSkin TF and SoftSkinAir® TF silicone liners offer the benefits of the Softskin range to a wider range of transfemoral users thanks to their increased size range.

With a shore hardness of 40, these liners allow long lasting comfort to the user with high transverse elasticity to adapt to residuum shape.

The TF Air offers the additional benefit of sweat-wicking micropores to reduce moisture, resulting in drier skin and lower temperature within the liner for residual limb comfort and skin irritation reduction.



Features and Benefits

- Soft residual limb padding
- Perfect sweat evacuation through micropores
- Innovative functional knitted fabric to reduce the longitudinal movement without an internal matrix
- Sizes: 28 - 55
- Unique functional outer fabric
- Machine-washable (30°C)
- Improved grip contact
- Decreased skin irritation

SoftSkinAir® TF

Additional Features and Benefits

- Breathable and comfortable solution
- Wicks sweat via micropores

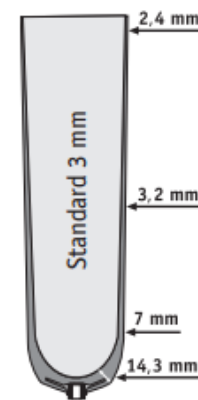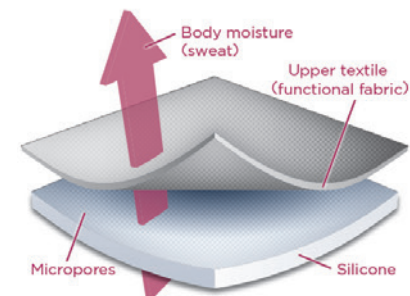

Ordering Information

All liners are 3mm thick, 350mm long.

Note: The circumference is measured 4cm above the distal end of the residual limb

| Size | SoftSkin TF | SoftSkinAir® TF |
|------|-------------|-----------------|
| 28 | NSTF-28 | NSTF-28-AIR |
| 30 | NSTF-30 | NSTF-30-AIR |
| 32 | NSTF-32 | NSTF-32-AIR |
| 34 | NSTF-34 | NSTF-34-AIR |
| 36 | NSTF-36 | NSTF-36-AIR |
| 38 | NSTF-38 | NSTF-38-AIR |
| 40 | NSTF-40 | NSTF-40-AIR |
| 42 | NSTF-42 | NSTF-42-AIR |
| 45 | NSTF-45 | NSTF-45-AIR |
| 50 | NSTF-50 | NSTF-50-AIR |
| 55 | NSTF-55 | NSTF-55-AIR |

Technical Information

| | SoftSkin TF | SoftSkinAir® TF |
|------------------|-----------------------|--------------------------|
| Material | Silicone | Silicone with micropores |
| Shore Hardness | 40 | 40 |
| Thickness | 3mm | 3mm |
| Length | 350mm | 350mm |
| Amputation Level | Transfemoral | Transfemoral |
| Distal Connector | With Distal Connector | With Distal Connector |
| Mobility Grades | 1,2,3 and 4 | 1,2,3 and 4 |
| Warranty | 6 months | 6 months |



**SoftSkinAir®
NS30/NS40/NS50**
Breathable Silicone-Liners

**SoftSkinAir®
NS30CL/NS40CL/NS50CL**
Breathable Silicone-Liners

**SoftSkinAir®
2L**
Breathable Silicone-Liners

**SoftSkinAir®
SPORT**
Breathable Silicone-Liner

**SoftSkinAir®
Cushion 40**
Breathable Silicone-Liners

**SoftSkinAir®
TF**
Comfortable & Breathable
Transfemoral Liners

| Material | Silicone (with Micropores) | | | | Silicone (with Micropores) | | Silicone (with Micropores) | | Silicone (With Micropores) | | Silicone (with Micropores) | | Silicone (with Micropores) | | | |
|-------------------|----------------------------|--|--------------------------|--|----------------------------|--|----------------------------|--|---|--|----------------------------|--|----------------------------|--|-------------|--|
| Shore Hardness | 30/40/50 | | 30/40/50 | | 30/50 | | 40 | | 40 | | 40 | | 40 | | | |
| Thickness | 3mm | | | | 3mm | | 3mm | | 3mm | | 3mm | | 3mm | | | |
| Length | 370mm | | 370mm | | 370mm | | 370mm | | 370mm (23.5 - 28mm) 400mm (28 - 36mm) 430mm (36mm<) | | 350mm | | 350mm | | | |
| Amputation Level | Transtibial | | Transtibial/Transfemoral | | | | Transtibial | | Transtibial | | Transtibial/Transfemoral | | Transfemoral | | | |
| Distal Connection | With Distal Connector | | Without Distal Connector | | | | With Distal Connector | | With Distal Connector | | Without Distal Connector | | With Distal Connector | | | |
| Mobility Grades | NS30 & NS40 1,2,3 | | NS50 2,3 and 4 | | NS30CL & NS40CL 1,2,3 | | NS50CL 2,3 and 4 | | 1,2 and 3 | | 2,3 and 4 | | 1,2 and 3 | | 1,2,3 and 4 | |
| Warranty | 6 months | | 6 months | | | | 6 months | | 6 months | | 6 months | | 6 months | | 6 months | |