

TASKA CX™

Faster. Stronger. More connected.



 **steeper**  **TASKA™**

by eqwal.

*This new hand has
been a real jump
in **quality and**
trust for me.*

Korrin Barrett

TASKA CX™

TASKA's most advanced prosthetic hand

The TASKA CX makes daily tasks easier, giving the user an enhanced sense of connection via fine control, contoured aesthetics, and personalisation features. Robust and waterproof, this prosthesis can perform a range of grips including the default general grip, flexitool, pincer, key and hook.

Not only is the TASKA CX lighter when compared to other multi-articulating hands, it also has a 20kg carry load capacity making it a strong hand with faster, more accurate control. The hand has built in wrist flexion to help reduce strain and fatigue from compensatory movements in the shoulder or torso.

Features and Benefits

- Available in black and sand with a choice of five different coverplate options
- Small (7¼ inch) and Medium (7¾ inch) sizing available
- MyTASKA App allows users to customise grip patterns
- Waterproof IP67 rating
- Quick Disconnect or Low Profile Wrist options



23 Grips

Complete most tasks with 3 core grips



20kg

Overall carry load



0 - 110mm*

Maximum opening range
*dependent upon hand size



IP67

Waterproof rating

IP67 Waterproof Rating

With an IP67 waterproof rating, users can wash dishes, water the garden or simply walk in the rain, safe in the knowledge they're not putting the hand at risk. It also makes cleaning the prosthesis easier as users can carry out the task with both hands together, resulting in significantly improved hygiene.

Also available are wrist options featuring TASKA's HydroSeal™ Waterproofing Technology, allowing submersion past the wrist. Please note, the hand is only waterproof to the wrist if used with a seal ring. Any proximal components are not protected.

Waterproof*

***Waterproof note:** Always wash the TASKA hand after using it in situations that could make it dirty (including salt water). Clean the hand before using it to eat, prepare food, or treat injuries using only water, soft soap, or light disinfectant.



Available in a Range of Finishes

Users can customise their TASKA CX with a choice of two base colours and five interchangeable coverplates. Coverplates are supplied as single plates in a single colour and side orientation, along with two removal tools.



Easily Configure in the MyTASKA App



Remote Connection

The Remote Connection feature in MyTASKA allows clinicians to configure user settings efficiently, conveniently, and 100% remotely. Thus, eliminating the need for travel, giving users more time to focus on what truly matters.



Clinician Settings

Input sensor readings are displayed on an accurate and easy to read graph to provide clearer insight of an individual's control. This supports clinicians with configuration adjustments, ensuring each TASKA hand can be tailored to the needs of each individual.




Customisable Grips

The app features a library of 18 preconfigured grips with the opportunity to add a further five custom options. From precision grips to power grips, users can simply select a grip of choice and customise it to their needs.




Built-in Games

Users can hone signal accuracy and improve control of their TASKA hand through app based games. There are three games, each designed to help the user refine their input signal clarity and practice proportional control.

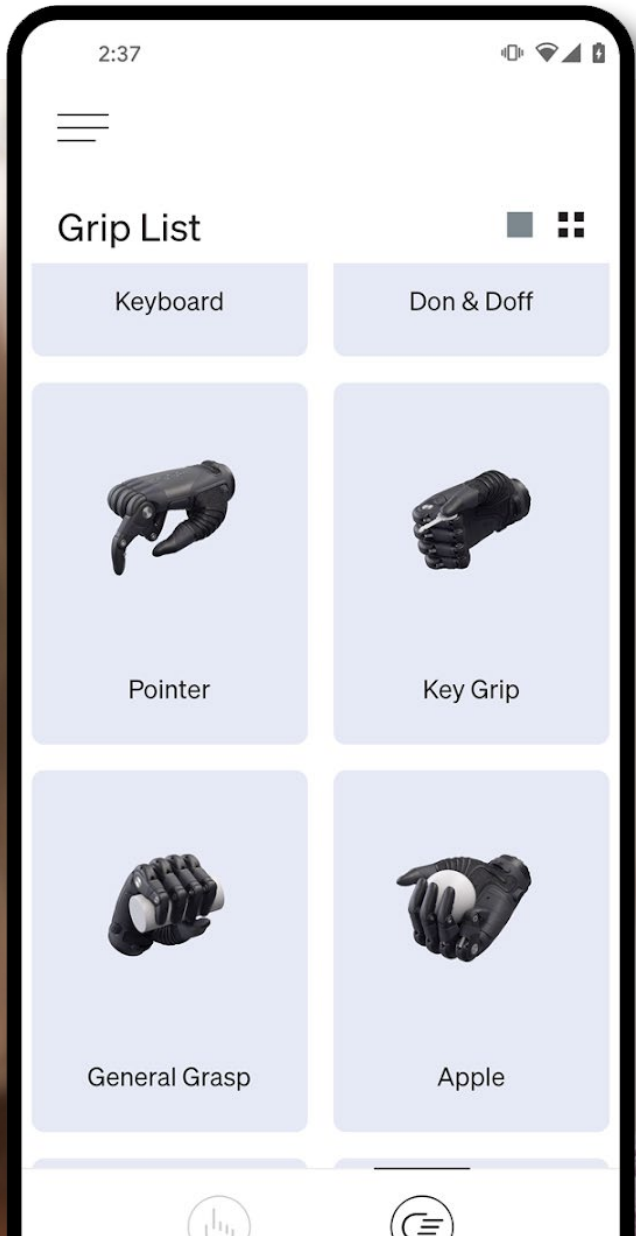

Download the MyTASKA app from your App Store

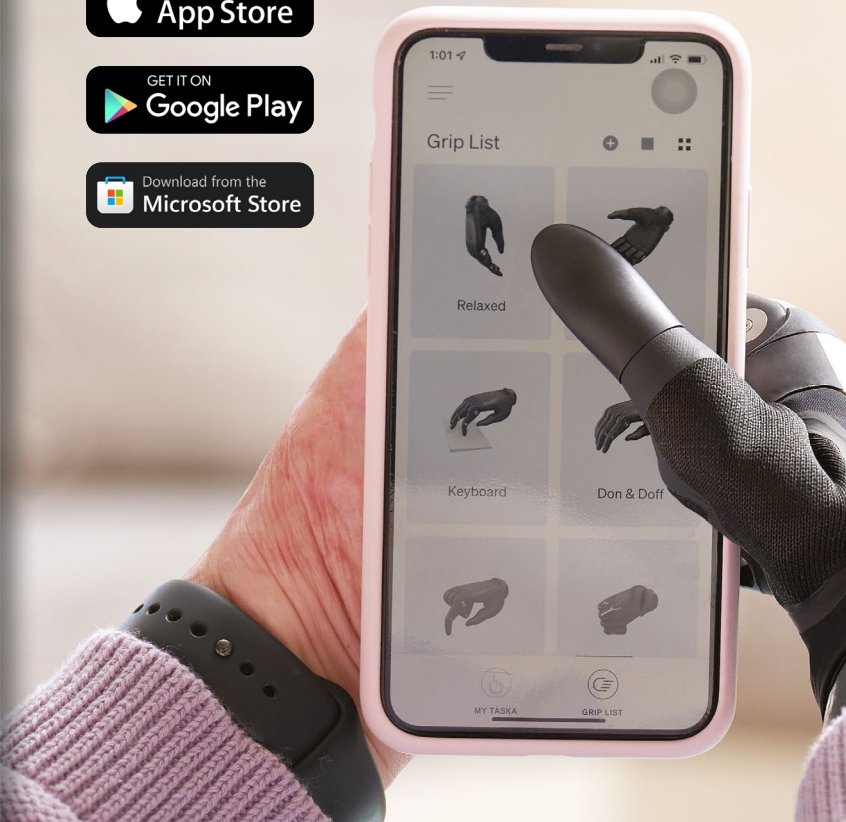

“It’s the **ultimate resource** for anyone wanting to fit a TASKA Hand.”


Mark Godenho, Head of Marketing

Become a TASKA Certified Clinician

Whether you want to become a TASKA Certified Clinician to fit TASKA hands, or to simply refresh your pre-existing knowledge on TASKA and their suite of products, the TASKA Learning Academy is the ideal online learning platform.


By registering, you’ll gain exclusive access to in-depth training materials, troubleshooting guides, and practical insights from experts in prosthetic technology. Equip yourself with the essential knowledge and confidence to fit the CX.

Visit www.taskaprosthetics.com/support/taska-learning-academy



Scan to start learning today

Technical Information

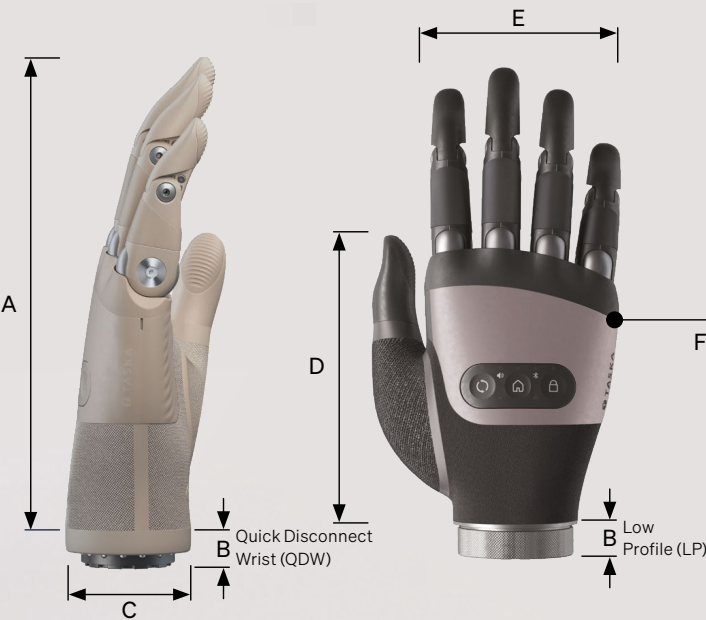

Dimensions			Small	Medium
A	Hand length	(QDW) (LP)	172mm 168mm	180mm 176mm
B	Hand base to wrist base	(QDW) (LP)	16mm 12.5mm	16mm 12.5mm
C	Diameter at wrist		48mm	54mm
D	Thumb tip length	(QDW) (LP)	112mm 108mm	114mm 110mm
E	Palm width		76mm	83mm
F	Knuckle circumference		184mm (7¼")	197mm (7¾")

Maximum carry load		20kg	
Grip forces	Thumb	24N	
	Finger	13N	
Digit speed	Thumb	0.52s or 120°/s	
	Finger	0.7s or 120°/s	
	Thumb rotator	0.55s or 160°/s	
Maximum opening	Small	100mm (3.94")	
	Medium	110mm (4.33")	
Wrist rotation range	(QDW) (LP)	Unlimited 300°	
Wrist extension lock position		15°	
Wrist Flexion lock position		25°	
Waterproof rating		IP67	
Operating temperature range		-10°C to 45°C (14°F to 113°F)	
Operating voltage		7.4V (nominal)	
Unit weight	Small	(QDW) (LP)	471g 430g
	Medium	(QDW) (LP)	533g 481g
Warranty		2 years (extendable to 5 years)	



Espire Elbow pictured with a TASKA HandGen2.

Compatible with a variety of myoelectric components

The TASKA CX is compatible with most industry leading control, elbow and power systems and passive or powered wrists. When paired with the **Steeper Espire Pro** and **COAPT** system, users can harness the full capabilities of myoelectric power.

For a full list of compatible components, visit www.taskaprosthetics.com/support/compatibility

Ordering Information



TASKA CX Hand

Included with each TASKA CX is a power system, charger, spare coverplate and replacement grip tips. When purchased with a low profile wrist, a lamination collar kit is also included.

Trial Period - Please speak to Customer Services or your Steeper Product Manager to organise a trial.

Ordering example: NTASKA-D1CB31-U1

Part ID		Size	Side	Wrist	Colour	Part ID	
NTASKA	-					31-U1	
7¼	D	Left	2	QDW	F	Black	B
7¾	C	Right	1	LP	C	Sand	C

Spare Parts Ordering Information



Coverplates

Supplied as single plates in a single colour and side orientation.

Ordering example: N193-20-109-GR

Size		Part ID		Side	Size		Colour
	-	20	-			-	
7¼	N193	Left	10	7¼	9	Black	BK
7¾	N203	Right	00	7¾	3	Grey	GR
						Dusk	DS
						White	WH
						Sand	SN



TASKA Glove

Size (In)	Side		Colour		Product Code
	Left	Right	Black	Sand	
7¼	•		•	•	N193-20-131-BK
				•	N193-20-132-SN
		•	•		N193-20-031-BK
				•	N193-20-032-SN
7¾	•		•		N203-20-106-BK
				•	N203-20-106-SN
		•	•		N203-20-006-BK
				•	N203-20-006-SN



Grip Tip Kit

Size (Inch)	Colour	Part Number
7¼	Black	N193-20-007-BK
	Sand	N193-20-007-SN
7¾	Black	N203-20-001-BK
	Sand	N203-20-001-SN
7¾ non-capacitive	Black	N203-20-007-BK
	Sand	N203-20-007-SN



Proximal Grip Kit

Size (Inch)	Colour	Part Number
7¼	Black	N193-20-008-BK
	Sand	N193-20-008-SN
7¾	Black	N203-20-002-BK
	Sand	N203-20-002-SN





Spare Parts Ordering Information - Wrists

Quick Disconnect Wrist - Lamination Collar Kits

1	Plug Type	Product Code
	6-band (45mm, black)	NTASKA-12BB01
	6-band (45mm, sand)	NTASKA-12BC01
	6-band (50mm, black)	NTASKA-13AB01
	6-band (50mm, sand)	NTASKA-13AC01

Quick Disconnect Wrist - Seal Ring Kits, Black

2a	Seal Ring Type	Product Code
	TASKA Seal Ring Kit (45mm, black)	N193-20-025
	TASKA Seal Ring Kit (50mm, black)	N195-20-003
	MC Seal Ring Kit (small, black)	N193-20-010
	MC Seal Ring Kit (medium/large, black)	N195-20-004

Quick Disconnect Wrist - Seal Ring Kits, Sand

2b	Seal Ring Type	Product Code
	TASKA Seal Ring Kit (45mm, sand)	N193-20-028
	TASKA Seal Ring Kit (50mm, sand)	N195-20-003-SN
	MC Seal Ring Kit (small, sand)	N193-20-029
	MC Seal Ring Kit (medium/large, sand)	N195-20-004-SN

Low Profile Wrist

	Description	Product Code
3	Low Profile Lamination Kit, CX, 45mm	N193-20-011
4	Low Profile Lamination Kit, CX, 50mm	NTASKA-13CD01
5	Low Profile Maintenance Kit, CX	N193-20-026
6	Switch Block	N152-20-014

Spare Parts Ordering Information - Accessories

User Accessories

	Description	Product Code
7	Mains Charger	NTASKA-MC-01
8	Car Charger	NTASKA-CC-01
9	Hand Travel Case	N192-20-002

Clinician Accessories

	Description	Product Code
10	TASKA Battery	NTASKA-BIG1-1
11	Power Switch	NTASKA-PWR-01
12	Battery Y-cable	N152-00-064
13	Lamination Dummy Blocks	NTASKA-MDB-01



Steeper Group
Unit 3 Stourton Link
Intermezzo Drive
Leeds, LS10 1DF
Email: customerservices@steepergroup.com
Tel: +44 (0) 870 240 4133

www.steepergroup.com

STPPR335

Official UK distributor of TASKA Prosthetics