

W Wrist Technology


### Nylon Child's Wrist

The lightweight Nylon Wrist is available in two sizes and features an adjustable collar which provides constant friction.

Part Number	Description
NINFANTWRIST	Infant's Wrist
NTRCHILDWRIST	Child's Wrist

Technical Information

Model	Infant's Wrist	Child's Wrist
Style	Round	Round
Diameter	3.2cm	3.8cm

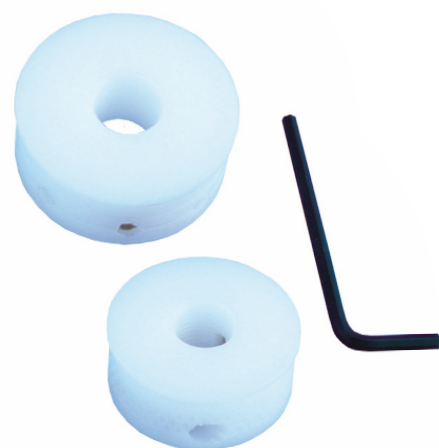

Features

- 1/2" x 20TPI thread
- Friction can easily be adjusted to users requirements with the supplied Allen key

W Wrist Technology

### Omega Friction /Disconnect Wrist

The Omega Friction/Disconnect Wrist series provides adjustable, rotational friction control with twelve locking positions. It uses reliable Omega disconnect and insert capture technology, to ensure a precise and solid fit. Disconnect features include lock-in, lock-out and lock-open features.

This device is manufactured using surgical stainless steel and aluminium, and features a standard laminating ring for fabrication.


Part Number	Model
NOMEGA-R	Omega Round Large
NOMEGA-L	Omega Oval Large
NOMEGA-M	Omega Oval Medium
NOM003	Omega Wrist Insert with 1/2" x 20TPI internal thread
NOM005	Omega Hook Adapter
NPTA005	Omega 'O' Ring

■ = Colour (Aluminium = Leave blank, B = Black)



Features

- Hook adapter for use on hooks and other narrow-based terminal devices to provide the broader platform necessary for the insert frictional adjustment function

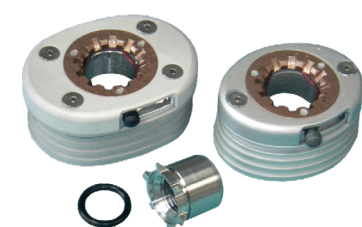

Omega Split Hook Adapter

Technical Information

Model	Omega Round Large	Omega Oval Large	Omega Oval Medium
Major Diameter	5.18cm	6cm	5cm
Minor Diameter	N/A	4.5cm	3.8cm
Total Height with Laminating Ring	2.62cm	2.5cm	2.5cm
Height (Face Plate Only)	0.8cm	0.8cm	0.8cm
Weight	103.1g	111.2g	89.4g

AC Waterproof Arm Cosmeses

### Arm Skinz


Arm Skinz are a waterproof arm cosmeses, designed to provide an improved cosmetic covering for body-powered and passive systems. Manufactured using a stain-resistant, tan coloured silicone material, they have optimal elasticity, can be easily cleaned, and are similar to donning a regular cosmetic glove.

Please note: these are not designed for manual labour or activities which may cause abrasion.

Part Number	Description
NARMSKINZ-AD	Arm Skinz Adult
NARMSKINZ-PD	Arm Skinz Paediatric

Technical Information

Model	Paediatric	Adult
Circumference	10.2 - 22cm	10.2 - 35.6cm
Length	43cm	
Thickness	Approx. 0.15cm	

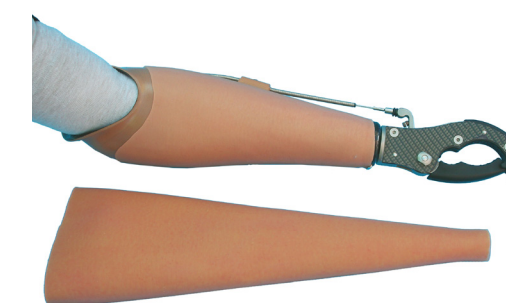

AC Pro Cuff

### Pro Cuff Prosthesis

Complete and ready to use, the Pro Cuff allows almost all TRS terminal devices to be positioned in an anatomically correct location within minutes, adding little to no additional length to the arm. It utilises a built-in BOA® cable lacing system for tension and closure, to allow for easy adjustability and quick donning/doffing.

Ideal for activities such as golfing, kayaking, swimming and ball/bat sports.



Features

- Stainless steel terminal device disconnect bracket with multiple mounting options
- Thermomouldable cuff
- 1/2" x 20TPI internal thread on adapter insert
- Order extra Adapter Inserts for quick and easy interchange of terminal devices (see page 25 for details)

Part Number
NPRO-CUFF-■

■ = Side (L = Left, R = Right)

Technical Information

Length	12.7cm
Weight (with stainless bracket)	187g
Weight (without bracket)	102g
Distal End Internal Circumference (adjustable)	-14.6 - -20.3cm
Proximal End Internal Circumference (adjustable)	-20.3 - -25.4cm

Applications

- Easily modified and adjustable, the universal size will fit most users
- Suitable for those with a temporary or permanent hand injury, paralysis, partial hand absence or wrist disarticulation
- Used to restore gripping dexterity when using the Pro Cuff with an appropriate terminal device