




We have adapted and developed our Elegance silicone gloves, and are pleased to introduce the new Elegance Plus cosmetics to our cosmetic glove range. These enhanced gloves are designed to fit the Steeper Myo Kinisi myoelectric hand, Steeper cable operated hands, and functional hands from other manufacturers.

Crafted using multiple layers of advanced silicone material and constructed with an integrated, reinforced inner mesh layer, these gloves are highly tear-resistant. This silicone material allows for additional stretch in comparison to a PVC glove, and so reduces the energy expended by the myoelectric hand motor each time it opens - making this a highly robust cosmetics particularly suited to higher activity patients.

Key Features and Benefits

- New integrated, reinforced inner mesh layer provides increased resilience against tearing and is highly robust
- Increased longevity, making this range highly cost effective
- Enhanced cosmetic finish with additional detailing on the palms, knuckles, nails and joints
- An easy-glide coating enables a patient to dress with ease
- Silicone surface can be easily wiped down when required making the glove suitable for use in a variety of environments
- Easier to don in comparison to the first generation Elegance silicone gloves
- Innovative reinforced mesh construction offers the toughness associated with PVC, along with the ease of cleaning provided by silicone
- Compatible with Steeper and other manufacturers' myoelectric and mechanical hands
- Available in the advanced micro-pigmentation TrueFinish™ for added realism and depth

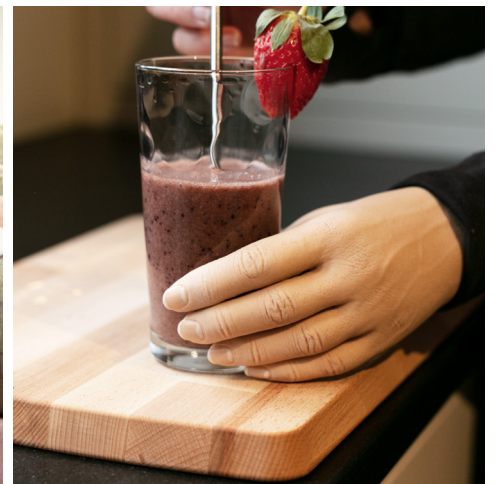

Elegance Plus Ordering Information

Elegance Plus Cosmetic Glove - Silicone TrueFinish™						
Hand Size	Approx Length	Application	Left	Right	Colour Range	Elegance Plus Suffix
7 ¼	645mm (25 ¾")	Adolescent/Adult Male	TSG306	TSG307	E0 to E18	-EP
7 ¼	645mm (25 ¾")	Adolescent/Adult Female	TSG304	TSG305	E0 to E18	-EP
7 ¾	620mm (24 ½")	Female	TSG312	TSG313	E0 to E18	-EP
7 ¾	645mm (25 ¾")	Male	TSG302	TSG303	E0 to E18	-EP
8 ¼	645mm (25 ¾")	Male	TSG300	TSG301	E0 to E18	-EP

Ordering Examples - E* denotes colour


TSG300/E13-EP - Orders an 8 ¼ Left Male Silicone Elegance Plus Glove in TrueFinish™ colour E13

TSG305/E7-EP - Orders and 7 ¼ Right Female Silicone Elegance Plus Glove in TrueFinish™ colour E7

Elegance Plus gloves can be requested with the following non-standard finishes: uncut, unpainted nails. If a non-standard finish or bilateral gloves are required, please specify when placing an order.

19 TrueFinish™ Colour Shades

For Elegance Plus gloves, the following 19 shades are available in the advanced TrueFinish™ micro-pigmentation, for enhanced realism and depth:

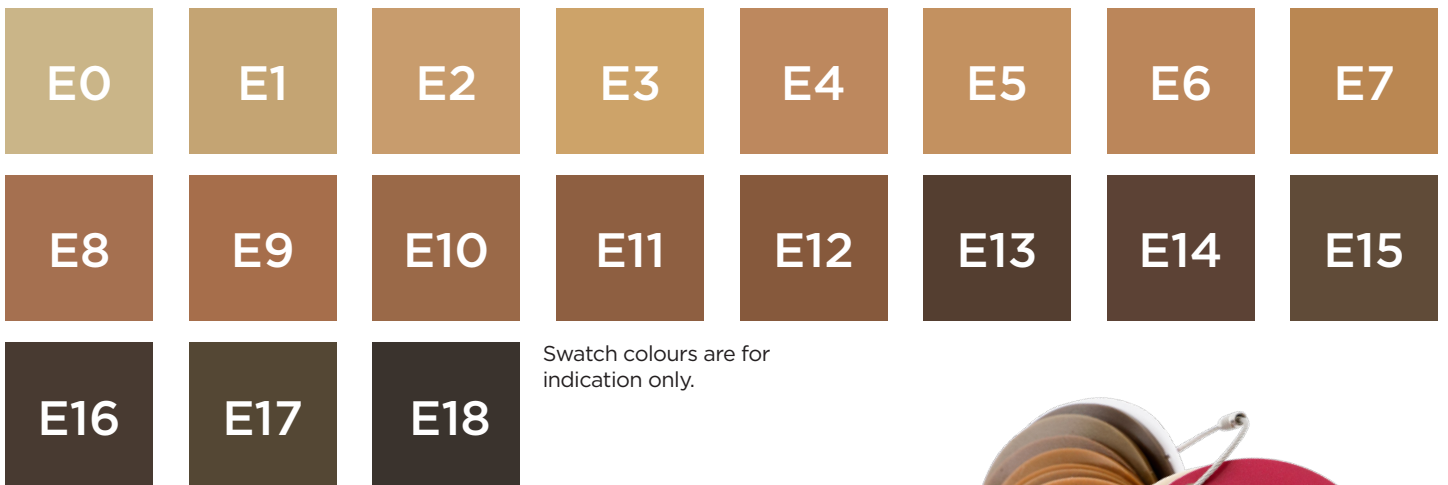

To help you select the most appropriate colour of glove, swatches are available and supplied without charge. These are manufactured from production silicone for the most accurate colour matching and will not fade or discolour over time.

TrueFinish™ Swatch Part Number: E24655

To request a swatch sample or to place an order for Elegance Plus Gloves, please contact customerservices@steepergroup.com.